

LISTEN.
THINK.
SOLVE.®

PharmaSuite®



RECIPE AND WORKFLOW MANAGEMENT

RELEASE 8.4

FUNCTIONAL REQUIREMENT SPECIFICATION

PUBLICATION PSFRSRD-RM008E-EN-E-DECEMBER-2017

Supersedes publication PSFDRD-RM008D-EN-E



Allen-Bradley • Rockwell Software

**Rockwell
Automation**

Contact Rockwell See contact information provided in your maintenance contract.

Copyright Notice © 2017 Rockwell Automation Technologies, Inc. All rights reserved.
This document and any accompanying Rockwell Software products are copyrighted by Rockwell Automation Technologies, Inc. Any reproduction and/or distribution without prior written consent from Rockwell Automation Technologies, Inc. is strictly prohibited. Please refer to the license agreement for details.

Trademark Notices FactoryTalk, PharmaSuite, Rockwell Automation, Rockwell Software, and the Rockwell Software logo are registered trademarks of Rockwell Automation, Inc.

The following logos and products are trademarks of Rockwell Automation, Inc.:

FactoryTalk Shop Operations Server, FactoryTalk ProductionCentre, FactoryTalk Administration Console, FactoryTalk Automation Platform, and FactoryTalk Security.
Operational Data Store, ODS, Plant Operations, Process Designer, Shop Operations, Rockwell Software CPGSuite, and Rockwell Software AutoSuite.

Other Trademarks ActiveX, Microsoft, Microsoft Access, SQL Server, Visual Basic, Visual C++, Visual SourceSafe, Windows, Windows 7 Professional, Windows Server 2008, Windows Server 2012, and Windows Server 2016 are either registered trademarks or trademarks of Microsoft Corporation in the United States and/or other countries.

Adobe, Acrobat, and Reader are either registered trademarks or trademarks of Adobe Systems Incorporated in the United States and/or other countries.

ControlNet is a registered trademark of ControlNet International.

DeviceNet is a trademark of the Open DeviceNet Vendor Association, Inc. (ODVA).

Ethernet is a registered trademark of Digital Equipment Corporation, Intel, and Xerox Corporation.

OLE for Process Control (OPC) is a registered trademark of the OPC Foundation.

Oracle, SQL*Net, and SQL*Plus are registered trademarks of Oracle Corporation.

All other trademarks are the property of their respective holders and are hereby acknowledged.

Warranty This product is warranted in accordance with the product license. The product's performance may be affected by system configuration, the application being performed, operator control, maintenance, and other related factors. Rockwell Automation is not responsible for these intervening factors. The instructions in this document do not cover all the details or variations in the equipment, procedure, or process described, nor do they provide directions for meeting every possible contingency during installation, operation, or maintenance. This product's implementation may vary among users.

This document is current as of the time of release of the product; however, the accompanying software may have changed since the release. Rockwell Automation, Inc. reserves the right to change any information contained in this document or the software at any time without prior notice. It is your responsibility to obtain the most current information available from Rockwell when installing or using this product.

-
-
- Rockwell Software PharmaSuite® 8.4 - Functional Requirement Specification Recipe and Workflow Management
-
-

Chapter 1	Introduction	1
	Typographical Conventions	1
Chapter 2	Recipe and Workflow Management (SR3146+)	3
	Recipe Designer - Batch (SR3071.6+).....	4
	Recipe Designer - Device (SR3071.9+).....	5
	Workflow Designer (SR3071.7+)	5
Chapter 3	Management of Master Recipes and Master Workflows (SR3146.4+)	7
	Master Recipe Model - S88 (SR3146.9.3).....	9
	Master Recipe - Header Attributes (SR3146.9.3.1).....	9
	Master Recipe - Review Mode (SR3146.9.3.2).....	11
	Master Workflow Model - S88 (SR3146.9.12)	11
	Master Workflow - Header Attributes (SR3146.9.12.1)	11
	Master Workflow - Review Mode (SR3146.9.12.2)	12
	Procedural Structure (SR3146.9.4)	13
	Process Inputs (SR3146.9.5)	15
	Process Outputs (SR3146.9.6+)	24
	Process Parameters (SR3146.9.7+).....	25
	Locking of Parameters (SR3146.9.13+).....	27
	Locking of Transitions (SR3146.9.14+)	27
	Workbench Object Operations (SR3146.9.10+)	28
	Workbench Object Operations - Master Recipe (SR3146.9.10.1)	28
	Workbench Object Operations - Master Workflow (SR3146.9.10.9)	29
	Workbench Object Operations - Element (SR3146.9.10.6)	31
	Workbench Object Operations - Link and Transition	32

Workbench Object Operations - Smart Replace (SR3146.9.10.10)	32
Workbench Object Operations - Replace Material (Product) (SR3146.9.10.11)	36
Material Flow Control (MFC) (SR3146.9.8).....	37
Graphical View (SR3146.9.8.1)	37
Table View (SR3146.9.8.1.1)	38
Tooltips (SR3146.9.8.2).....	39
Automatic Merge - Default MFC Definition (SR3146.9.8.3)	39
Manual Merge	40
Split (SR3146.9.8.8)	42
Reset (SR3146.9.8.9).....	42
Automated Update of MFC Representation (SR3146.9.8.10)	42
MFC Definition Linked to Procedure Level (SR3146.9.8.12)	42
Information Flow (SR3146.9.9+)	42
Output Mechanism (SR3146.9.9.1).....	43
Referencing (SR3146.9.9.2)	43
Validator (SR3146.9.9.3)	43
Standard Statuses for Version Control (SR3146.6.4)	43
Default Version Graph for Master Recipes (SR3146.6.4.1)	44
Default Version Graph for Master Workflows (SR3146.6.4.4).....	46
Validator Consistency Checks (SR3146.9.4.3+)	48
Checks that Apply to SFC Modeling	48
Checks that Apply to Transition Conditions	49
Checks that Apply to Expressions within Parameter Attributes.....	50
Checks that Apply to MFC Modeling of a Procedure within Recipe and Workflow Designer	50
Checks that Apply to Master Recipes in General	53
Checks that Apply to Master Recipes Based on an ERP BOM.....	54
Checks that Apply to Material Parameters of a Phase.....	56
Checks that Apply to Specific Dispense Building Blocks.....	60
Checks that Apply to Inline Weighing	62
Checks that Apply to Output Weighing.....	64

Checks that Apply to Non-production Master Recipes in General	65
Checks that Apply to Master Workflows in General.....	65
Checks that Apply to Signature Privileges of a Phase	65
Checks that Apply to Equipment Requirements of a Phase.....	66
Checks that Apply to Work Center Assignments	68
Checks that Apply to Capability Assignments in General	68
Checks that Apply to Escalation-enabled Capability Assignments.....	70
Checks that Apply to Trigger-enabled Capability Assignments.....	70
Checks that Apply to Server-run Capability Assignments	71
Chapter 4 Management of Building Blocks (SR3146.10+).....	73
Header Attributes (SR3146.10.4)	73
Workbench Object Operations (SR3146.9.10.7)	74
Workbench Object Operations - New Building Block Dialog (SR3146.9.10.7.2)	75
Workbench Object Operations - View Locked Building Block (SR3146.9.10.7.1)	75
Status Handling of Building Blocks (SR3146.10.1)	75
Revisioning (SR3146.10.2).....	76
Usage List (SR3146.10.3)	76
Content (SR3146.10.3.1)	77
Objects Information (SR3146.10.3.2)	77
Filter (SR3146.10.3.3)	77
Chapter 5 Comparing Master Recipes, Master Workflows, and Custom Building Blocks (SR3146.15+)	79
Comparison Algorithm (SR3146.15.1+)	79
Comparing Header Attributes (SR3146.15.1.1+).....	79
Comparing Element Attributes (SR3146.15.1.2)	80
Comparing Source Attributes (SR3146.15.1.3).....	80
Comparing Parameter Attributes (SR3146.15.1.4+).....	80
Comparing Procedural Elements (SR3146.15.1.5+)	81
Defining a Comparison Baseline (SR3146.15.2)	82

	Performing a Comparison (SR3146.15.3)	82
Chapter 6	Management of Change Requests (SR3146.12+)	83
	Header Attributes (SR3146.12.1)	83
	Workbench Object Operations (SR3146.12.6).....	84
	Open Dialog for Protected Objects (SR3146.12.6.2).....	84
	Workbench Object Operations - View Locked Change Request (SR3146.12.6.1)	85
	Status Management of Change Requests (SR3146.12.2+)	85
	Master Recipe-specific Change Requests (SR3146.12.2.2)	85
	Master Workflow-specific Change Requests (SR3146.12.2.3).....	87
	Scope Validity Check (SR3146.12.2.1)	89
	Scope Definition (SR3146.12.3+)	89
	Define Old/New Building Blocks (SR3146.12.3.1)	89
	Select from Usage List (SR3146.12.3.2).....	90
	Action List (SR3146.12.4)	90
	Execute Change Request (SR3146.12.5)	91
	Execution Status (SR3146.12.5.2).....	93
	Cancel Change Request (SR3146.12.5.1)	93
Chapter 7	Master Recipe Report - Batch (SR3146.11+).....	95
	Print MR Report (SR3146.11.1)	95
	MR Report - SFC Pagination (SR3146.11.1.1).....	95
	MR Report Header (SR3146.11.2)	95
	MR Report Footer (SR3146.11.3)	95
	MR Report Cover Page (SR3146.11.5)	96
	MR Report Sections (SR3146.11.4+)	96
Chapter 8	Master Recipe Report - Device (SR3146.14+)	103
	Print MR Report (SR3146.14.1)	103
	MR Report - SFC Pagination (SR3146.14.1.1).....	103
	MR Report Header (SR3146.14.2)	103
	MR Report Footer (SR3146.14.3)	103

MR Report Cover Page (SR3146.14.5)	104
MR Report Sections (SR3146.14.4+)	104
Chapter 9 Master Workflow Report (SR3146.13+)	111
Print MWF Report (SR3146.13.1)	111
MWF Report - SFC Pagination (SR3146.13.1.1)	111
MWF Report Header (SR3146.13.2)	111
MWF Report Footer (SR3146.13.3)	111
MWF Report Cover Page (SR3146.13.5)	112
MWF Report Sections (SR3146.13.4+)	112
Chapter 10 Designer Workbench (SR3146.9+)	117
Workbench Usability (SR3146.9.1+)	117
Building Block Image	117
Process View (SR3146.9.1.1)	119
Navigation (SR3146.9.1.2)	119
Context menus (SR3146.9.1.3)	120
Docking (SR3146.9.1.4)	120
Tooltips (SR3146.9.1.5)	121
Workbench UI Elements (SR3146.9.2+)	121
Work Area (SR3146.9.2.1)	121
Open Dialog (SR3146.9.2.20+)	122
Library Universe (SR3146.9.2.4+)	122
Setlist (SR3146.9.2.9)	124
Parameter Panel (SR3146.9.2.7)	126
Property Windows	127
Messages (SR3146.9.2.16+)	128
Comparison (SR3146.9.2.19+)	130
Preview Panel (SR3146.9.2.11)	131
Expression Editor (SR3146.9.9.4)	131
Other Editors (SR3146.9.2.17+)	138
Statistics (SR3146.9.2.18)	140

Version Control (SR3146.4.5)	141
Status Change (SR3146.4.5.1)	142
Version History (SR3146.4.5.3)	142
Configurable Version Control (SR3146.6)	142
Flexible Version Graph Assignment (SR3146.6.1).....	143
Flexible Definition of Approval Details (SR3146.6.7)	143
Library Concept (SR3146.9.11+)	143
Basic Building Block Import (SR3146.9.11.1)	143
Custom Building Block (SR3146.9.11.3)	143
Chapter 11 Export/Import of Master Recipes, Master Workflows, and Building Blocks (SR1075.1+).....	145
General Prerequisites and Usage of B2MML (SR1075.1.1).....	145
Prerequisites For Import (SR1075.1.1.1).....	145
Checksum (SR1075.1.1.2)	145
Handling of Status and Version (SR1075.1.1.3)	146
Handling of Status Transition History (SR1075.1.1.4)	146
Chapter 12 Reference Documents	147
Chapter 13 Document Information	149
Approval	149
Version Information	149
Revision History.....	149
Index	159

Figure 1: Objects relevant to recipe management	8
Figure 2: Objects relevant to workflow management.....	9
Figure 3: Master recipe - status graph	45
Figure 4: Master workflow - status graph	47
Figure 5: Building block - status graph	76
Figure 6: Master recipe-specific change request - status graph	87
Figure 7: Master workflow-specific change request - status graph	89
Figure 8: Building block image	117
Figure 9: Operation building block image with capabilities.....	118
Figure 10: Unit procedure building block image with capabilities	118

-
-
- Rockwell Software PharmaSuite® 8.4 - Functional Requirement Specification Recipe and Workflow Management
-
-

Introduction

This document details the requirements of the functions implemented related to recipe and workflow management with PharmaSuite Recipe and Workflow Designer.

Each requirement is composed of a name and a unique identifier (e.g. Procedural structure (SR3146.9.4)). If a requirement's meaning is for requirement grouping only, the identifier is appended by a plus sign (e.g. Messages (SR3146.9.2.16+)).

In some cases, additional context information is available, indicated in the document by a frame and a gray background color. This context information is related to the respective requirement, but not part of the formal requirement description.

The revision history (page [149](#)) lists the changes made to the document with PharmaSuite 8.3 as the comparison baseline. Changes related to a requirement are marked as "Editorial", "Update", "New", or "Deleted", changes to the additional context information are marked as "Context information-related".

Typographical Conventions

This documentation uses typographical conventions to enhance the readability of the information it presents. The following kinds of formatting indicate specific information:

Bold typeface	Designates user interface texts, such as
	■ window and dialog titles
	■ menu functions
	■ panel, tab, and button names
	■ box labels
	■ object properties and their values (e.g. status).

-
-
- Rockwell Software PharmaSuite® 8.4 - Functional Requirement Specification Recipe and Workflow Management
-
-

Recipe and Workflow Management (SR3146+)

➤ For recent changes, see revision history (page 149).

In PharmaSuite, recipe and workflow management is done by means of the **Designer Workbench (SR3146.9+)** tool (page 117).

S95 defines a master recipe as follows: The master recipe is that level of recipe that is targeted to a process cell or a subset of the process cell equipment.

In its role as container for both the what (to manufacture) and the how (to manufacture), the master recipe is the place where the flow of the materials through the production process is specified. Thus a master recipe typically contains all pharmaceutically-relevant data required to support GxP-compliant production. Master recipes are under version and status control.

Apart from master recipes PharmaSuite supports the definition and execution of workflows that also follow the S88 structure of a master recipe (master workflow).

Recipe and Workflow Designer of PharmaSuite is a graphical workbench for building and maintaining master recipes, master workflows, and their component building blocks. With its functions for status and version control, it covers the entire life cycle of a master recipe, master workflow, or building block. The workbench provides material flow control on the basis of material parameters, information flows for values from process parameters, privilege parameters for access rights, capability parameters, and equipment parameters.

For each structure level of a master recipe, master workflow, or building block, the system presents an SFC graph (sequential function chart).

PharmaSuite Recipe and Workflow Designer is available in several modes: as Recipe Designer - Batch, Recipe Designer - Device, and Workflow Designer.

Most of the requirements that are related to Recipe Designer - Batch also apply to Recipe Designer - Device and Workflow Designer. If a requirement applies only to a subset (Recipe Designer - Batch, Recipe Designer - Device, Workflow Designer), this is explicitly noted within the requirement.

CONFIDENTIAL OBJECTS

The concept of confidential objects protects the intellectual property of recipes, workflows, building blocks, orders, ERP BOMs, and related data from unauthorized access. The system provides a specific access rights type to define access privileges that allow users to maintain protected recipes, workflows, and orders.

A master recipe or master workflow can be protected by assigning an explicit access privilege or inherit the access privilege from the used ERP BOM or the selected material (product). Orders and workflows based on a protected master recipe or master workflow inherit the access privilege automatically. A user without a protected-specific access privilege is not able to

- view a protected master recipe or master workflow,
- start a protected order or workflow for execution

- log in at a station on which a protected order or workflow runs,
- see a protected order or workflow in the Production Response Client
- print the batch report, weighing report, or workflow report of a protected order or workflow,
- see any transaction history data created by a protected order or workflow.

In case several confidential objects are linked with each other (e.g. a confidential building block is used in a recipe or a confidential workflow is assigned to an order), only the access privilege defined for the master object is required to access the linked objects with the master object as starting point. In order to link two confidential objects, they do not need to be protected by the same access privilege. The user, however, needs to have both access privileges.

Example scenario: A workflow and an order have different access privileges. A user with access to both can assign the workflow to the order. Afterwards another user who only has the access privilege of the order is able to print the batch report and can review exceptions of the workflow if the review is started with the order as starting point.

STRUCTURE OF THE DOCUMENT

This description is structured in the following sections:

- **Management of master recipes and master workflows (SR3146.4+) requirements** (page [7](#))
- **Comparing master recipes, master workflows, and custom building blocks (SR3146.15+) requirements** (page [79](#))
- **Management of building blocks (SR3146.10+) requirements** (page [73](#))
- **Management of change requests (SR3146.12+) requirements** (page [83](#))
- **Master recipe report - Batch (SR3146.11+) requirements** (page [95](#))
- **Master recipe report - Device (SR3146.14+) requirements** (page [103](#))
- **Master workflow report (SR3146.13+) requirements** (page [111](#))
- **Designer Workbench (SR3146.9+) requirements** (page [117](#))
- **Export/Import of Master Recipes, Master Workflows, and Building Blocks (SR1075.1+) requirements** (page [145](#))

Recipe Designer - Batch (SR3071.6+)

Access control (SR3071.6.2)

The access control of the Recipe Designer - Batch workbench enables the definition of user-specific access rights.

Upon user login, the Recipe Designer - Batch workbench checks if a user is authorized to edit recipes. If not, the system provides restricted access that does not allow to save manipulated data. Only status transitions of master recipes, building blocks, and change requests can be triggered by authorized users.

Recipe Designer - Device (SR3071.9+)

Access control (SR3071.9.2)

The access control of the Recipe Designer - Device workbench enables the definition of user-specific access rights.

Upon user login, the Recipe Designer - Device workbench checks if a user is authorized to edit recipes. If not, the system provides restricted access that does not allow to save manipulated data. Only status transitions of master recipes, building blocks, and change requests can be triggered by authorized users.

Workflow Designer (SR3071.7+)

Access control (SR3071.7.2)

The access control of the Workflow Designer workbench enables the definition of user-specific access rights.

Upon user login, the Workflow Designer workbench checks if a user is authorized to edit workflows. If not, the system provides restricted access that does not allow to save manipulated data. Only status transitions of master workflows, building blocks, and change requests can be triggered by authorized users.

-
-
- Rockwell Software PharmaSuite® 8.4 - Functional Requirement Specification Recipe and Workflow Management
-
-

Management of Master Recipes and Master Workflows (SR3146.4+)

- As requirements have been restructured compared to the original DRS document, this section also contains requirements of other requirement groups than SR3146.4+.

This section provides the requirements relevant to the management of master recipes and master workflows:

- **Master recipe model - S88 (SR3146.9.3) requirements** (page 9)
- **Master workflow model - S88 (SR3146.9.12) requirements** (page 11)
- **Procedural structure (SR3146.9.4) requirements** (page 13)
- **Workbench object operations (SR3146.9.10+) requirements** (page 28)
- **Material flow control (MFC) (SR3146.9.8) requirements** (page 37)
- **Information flow (SR3146.9.9+) requirements** (page 42)
- **Standard statuses for version control (SR3146.6.4) requirements** (page 43)
- **Validator consistency checks (SR3146.9.4.3+) requirements** (page 48)

A master recipe has relationships to the following objects:

- **Objects to build a procedural structure - building blocks** (page 13):
Procedure, unit procedure, operation, phase
- **Capabilities to control the behavior during execution**
 Detachable - Operation (SR3146.9.5.9.1) capability (page 20)
 Event-triggered (SR3146.9.5.9.2) capability (page 21)
 Escalation-enabled (SR3146.9.5.9.3) capability (page 21)
 Trigger-enabled (SR3146.9.5.9.4) capability (page 22)
 Server-run (SR3146.9.5.9.5) capability (page 22)
 Auto-startable (SR3146.9.5.9.6) capability (page 22)
 Pause-enabled (SR3146.9.5.10.1) capability (page 23)
 Detachable - Unit procedure (SR3146.9.5.10.2) capability (page 23)
 Abort-and-reactivate-enabled (SR3146.9.5.10.3) capability (page 23)
- **Resource master data**
 Material (process input (SR3146.9.5.1) data (page 16), process output (SR3146.9.6.1) data (page 24))
 Equipment requirement (class) (SR3146.9.5.5) data (page 19)
 Work center (SR3146.9.5.4) data (page 18)

Station (SR3146.9.5.8) data (page 18)

- Material flow control (SR3146.9.8) data (page 37)
- Signature privilege (SR3146.9.5.3) data (page 17)

The diagram illustrates the connections between the objects.

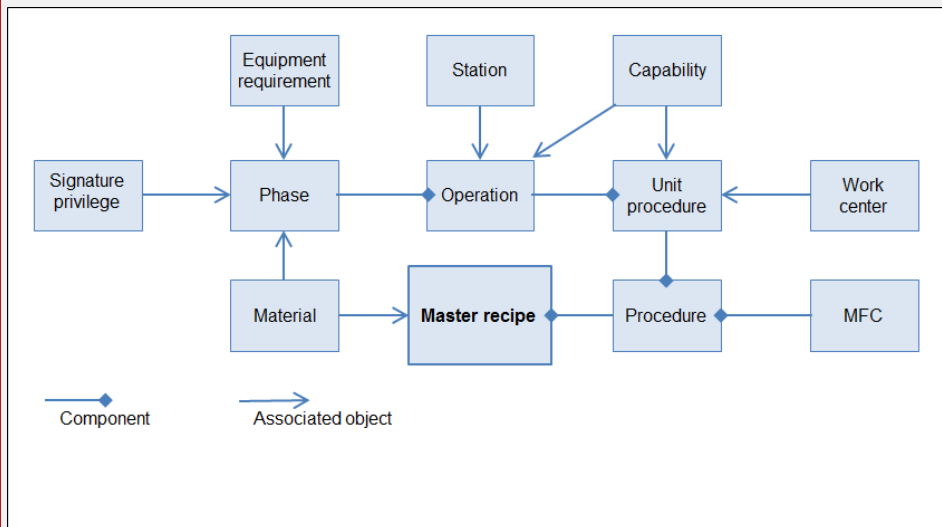


Figure 1: Objects relevant to recipe management

A master workflow has relationships to the following objects:

- Objects to build a procedural structure - building blocks (page 13):
Procedure, unit procedure, operation, phase
- Resource master data
Material (process input (SR3146.9.5.1) data (page 16), process output (SR3146.9.6.1) data (page 24))
Equipment requirement (class) (SR3146.9.5.5) data (page 19)
Work center (SR3146.9.5.4) data (page 18)
Station (SR3146.9.5.8) data (page 18)
- Capabilities to control the behavior during execution
Detachable - Operation (SR3146.9.5.9.1) capability (page 20)
Event-triggered (SR3146.9.5.9.2) capability (page 21)
Escalation-enabled (SR3146.9.5.9.3) capability (page 21)
Trigger-enabled (SR3146.9.5.9.4) capability (page 22)
Server-run (SR3146.9.5.9.5) capability (page 22)
Auto-startable (SR3146.9.5.9.6) capability (page 22)
Pause-enabled (SR3146.9.5.10.1) capability (page 23)
Detachable - Unit procedure (SR3146.9.5.10.2) capability (page 23)
- Material flow control (SR3146.9.8) data (page 37)
- Signature privilege (SR3146.9.5.3) data (page 17)

The diagram illustrates the connections between the objects.

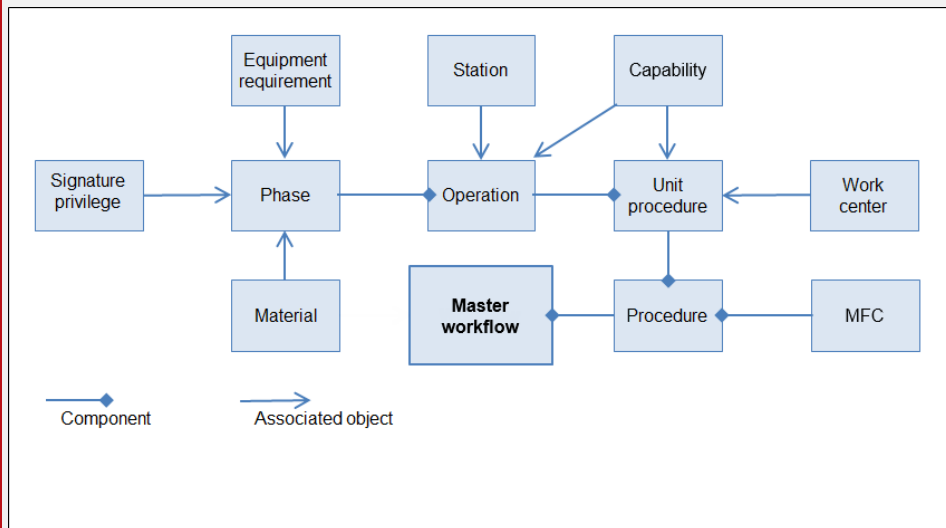


Figure 2: Objects relevant to workflow management

Master Recipe Model - S88 (SR3146.9.3)

The applied master recipe model is S88-compliant and includes

- **header attribute** (SR3146.9.3.1) data (page 9),
- a **procedural structure** (SR3146.9.4) hierarchy (page 13), and
- the following formula-related data:
 - **process input** (SR3146.9.5) data (page 15),
 - **process output** (SR3146.9.6+) data (page 24),
 - **process parameter** (SR3146.9.7+) data (page 25).

The model also applies to the respective top-level custom BB in support of **BB management** (SR3146.10+) capabilities (page 73).

Master Recipe - Header Attributes (SR3146.9.3.1)

- For recent changes, see revision history (page 149).

The master recipe object supports the following header attributes:

- Basic data:
 - Identifier
 - Version
 - Description

- Comparison baseline
Version is displayed as <major revision>.x.
- Method
- Method description
- Reason for creation
- Usage type (Production, Cost center)
Does not apply to Recipe Designer - Device.
- Access privilege
- Status data:
 - Status
 - Valid from/until
- Additional data:
 - Planned quantity
For Recipe Designer - Device: 1 ea (or any unit that is convertible into ea, read-only).
 - Minimum quantity
 - Maximum quantity
 - Registration number
 - Review mode (Automatic, Manual) (page [11](#))
 - Comment
- Material (product):
 - Material identifier
 - [Material] Short description
 - ERP BOM alternative
 - ERP BOM base quantity
 - ERP BOM released
- Packaging Levels (editable with the **Packaging Level Data (SR3146.9.2.17.3)** editor (page [140](#))):
 - L6 (Meaning, Contained number, Inventory level, Hide during execution)
 - ...
 - L9 (Meaning, Contained number, Inventory level, Hide during execution)

The attributes are maintained in the **header component-specific (SR3146.9.2.15)** property window (page 128).

Master Recipe - Review Mode (SR3146.9.3.2)

The system allows to define how an order that is based on a master recipe can be reviewed:

- Automatic
The order will be reviewed automatically.

The status of the order is automatically changed to **Reviewed**, once the order has been set to **Finished** and no exceptions with a risk other than **None** have been recorded for the order or any of its appended workflows.
See also **Automatic Review of an Order (SR1084.4.1.1)** in "Functional Requirement Specification Runtime Data Management" [A3] (page 147).

- Manual
The order should be reviewed manually.

Master Workflow Model - S88 (SR3146.9.12)

The applied master workflow model is S88-compliant and includes

- **header attribute (SR3146.9.12.1)** data (page 11),
- a **procedural structure (SR3146.9.4)** hierarchy (page 13), and
- the following formula-related data:
 - **process input (SR3146.9.5)** data (page 15),
 - **process output (SR3146.9.6+)** data (page 24),
 - **process parameter (SR3146.9.7+)** data (page 25).

The model also applies to the respective top-level custom BB in support of **BB management (SR3146.10+)** capabilities (page 73).

Master Workflow - Header Attributes (SR3146.9.12.1)

- For recent changes, see revision history (page 149).

The master workflow object supports the following header attributes:

- Basic data:
 - Identifier
 - Version

- Processing name
Preset with the master workflow identifier. Once the Processing name differs from the workflow identifier, the system does not change it automatically when the master workflow is copied or renamed.
- Short description
- Description
- Comparison baseline
Version is displayed as <major revision>.x.
- Reason for creation
- Access privilege
- Status data:
 - Status
 - Valid from/until
- Additional data:
 - Review mode (Automatic, Manual) (page 12)
 - Type
 - One-click startable
 - Execution prefix
 - Production-relevant
 - Appendable during processing
 - Comment

The attributes are maintained in the **header component-specific (SR3146.9.2.15)** property window (page 128).

The **Type** attribute defines to which group the one-click startable workflow belongs in the Cockpit of the Production Execution Client. See also **Grouping (SR3071.4.1.1.4)** in "Functional Requirement Specification Execution Framework" [A2] (page 147).

The **Appendable during processing** attribute can only be set if the **Production-relevant** attribute is set.

Master Workflow - Review Mode (SR3146.9.12.2)

The system allows to define how a workflow can be reviewed:

- Automatic
The workflow be reviewed automatically.

The status of the workflow is automatically changed to **Reviewed**, once the workflow has been set to **Finished** and no exceptions with a risk level other

than **None** have been recorded for the workflow.
See also **Automatic Review of an Order (SR1084.4.1.1)** in "Functional Requirement Specification Runtime Data Management" [A3] (page 147).

- **Manual**

The workflow should be reviewed manually.

Procedural Structure (SR3146.9.4)

The supported procedural structure includes

- **Building blocks** (page 13)
(procedure, unit procedure, operation, phase)
- **Links and branches** (page 14)
(connector, loop, selection branch, simultaneous branch)
- **Transitions** (page 14)

For master workflows, only one procedure and one unit procedure can be defined.

Procedural structure - Building blocks

The building blocks are:

- **Procedure**
The procedure is the highest level in the (recipe) hierarchy and defines the strategy for carrying out a major processing action such as making a batch. It is defined in terms of an ordered set of unit procedures.
- **Unit procedure**
According to S88, a unit procedure consists of an ordered set of operations that causes a contiguous production sequence to take place within a unit. Only one operation is presumed to be active in a unit at any time. An operation is carried to completion in a single unit. However, multiple unit procedures of one procedure may run concurrently, each in different units.
However, PharmaSuite allows modeling and execution of parallel operations in compliance with general SFC rules.
- **Operation**
An operation is an ordered set of phases that defines a major processing sequence that takes the material being processed from one state to another, usually involving a chemical or physical change. It is often desirable to locate operation boundaries at points in the procedure where normal processing can safely be suspended.
See also "Procedural Structure - Operation Type (SR3146.9.4.4)" (page 15).
- **Phase**
The smallest element of procedural control that can accomplish a process-oriented task is a phase.

Procedural structure - Links and Branches

The Designer workbench uses various types of links to connect building blocks:

- **Connector**
A link is a simple connector between two components.
- **Loop**
A loop is a specific type of link that points back to a previous step in the graph. The loop endpoint must be on the same sequence as the loop start point, so you cannot draw a loop out of a branch.
- **Selection branch**
A selection branch opens and closes two or more exclusive alternative sequences of steps in the graph. During execution, only one of the sequences must be executed before an operator can proceed with the step that follows after the sequences have been joined again. The result of the transition conditions that precede a branch determines the subsequent step.
- **Simultaneous branch**
A simultaneous branch opens and closes two or more parallel sequences of steps in the graph. The leftmost sequence in the graph is the first to be displayed during execution according to the default execution order of SFCs (left to right). During execution, all parallel sequences must be executed before an operator can proceed with the step that follows after the sequences have been joined again.

Procedural structure - Transitions

A transition consists of a condition that defines which step is the next to be executed when there is more than one potential successor step. Thus a transition predetermines which step to choose in a selection branch during execution or whether a loop needs to be executed. Transitions with conditions and or a description are equipped with transition identifiers. Any transition can be locked. For details about transition conditions, see "Procedural Structure - Conditions (SR3146.9.4.2)" (page 15), "Locking of Transitions (SR3146.9.14+)" (page 27), and "Expressions for Transitions (SR3146.9.9.4.2)" (page 132).

Procedural structure - Parallel execution (SR3146.9.4.1)

Building blocks can be executed in parallel.

Procedural structure - Conditions (SR3146.9.4.2)

Transitions can include conditions based on calculations.

- Transition conditions are based on the results of a previous step.
- Each transition is preset with a default condition, which defines that the step(s) preceding a transition must have been finished before the step(s) after the transition can be started.
Transition conditions can be edited in an **Expression editor (SR3146.9.9.4)** UI element (page 131).
- Non-default conditions are marked by an indicator.

Procedural structure - Hidden phase (SR3146.9.4.5)

During recipe or workflow authoring, the system provides a Hidden Phase, which is a structural phase that is required to build adjacent branches or loops in valid SFC syntax. Thus it is only visible in Recipe and Workflow Designer and does not appear as executable phase during order or workflow processing.

Two different Hidden Phases are available, one for the execution on a client and one for the execution on a server.

Procedural structure - Operation type (SR3146.9.4.4)

- Does not apply to Recipe Designer - Device.
- Does not apply to Workflow Designer.

Operations are equipped with a **Dispense** attribute that marks an operation as being intended for dispensing if set. The recipe author has to define the operation type.

- The default setting is **non-Dispense**.
- The **Dispense** attribute is automatically populated to the superordinate unit procedure. For unit procedures, the **Dispense** attribute is read-only.

The following rules apply:

- A Dispense operation cannot be combined with other non-Dispense operations.
- A Dispense unit procedure must contain exactly one Dispense operation.

Process Inputs (SR3146.9.5)

According to S88, process inputs are supported as part of the formula.

The attributes of process input parameters are maintained in the **Parameter Panel (SR3146.9.2.7)** UI element (page [126](#)).

PROCESS INPUTS - MATERIAL (SR3146.9.5.1)

Materials are supported as process input parameters.

The minimum and maximum numbers of process inputs are taken over from the phase building block.

PROCESS INPUTS - MATERIAL ATTRIBUTES (SR3146.9.5.1.1)

➤ For recent changes, see revision history (page [149](#)).

The material input parameters support the following attributes:

- Basic data
 - Lock
 - Material identifier (read-only)
 - [Material] Short description (read-only)
 - Planned quantity mode (None, As defined, As produced)
 - Planned quantity (with UoM)
 - Usage type (e.g. Input, Output) (read-only)
 - Usage name
 - MFC-relevant (Yes/No), default = No
 - Position
 - Hierarchy level (read-only)
 - Path (read-only)
- Weighing:
 - Default weighing method
 - Allowed weighing methods
 - Tolerances (lower and upper, absolute and relative)
 - Weighing material type
 - Planned potency
 - Fixed quantity (Yes/No)
 - Target weight
(not relevant to material inputs and ignored during execution)

- Target weight tolerances (lower and upper, absolute and relative)
(not relevant to material inputs and ignored during execution)
- Packaging Levels (editable with the **Packaging Level Data (SR3146.9.2.17.3)** editor (page 140)):
 - (not relevant to material inputs and ignored during execution)
 - L0 (Meaning, Contained number, Inventory level, Hide during execution)
(the Contained number attribute is not relevant to the L0 level)
 - ...
 - L5 (Meaning, Contained number, Inventory level, Hide during execution)

- The **Basic** attributes are displayed in both the horizontal and vertical panes of the **Parameter Panel (SR3146.9.2.7)** UI element (page 126).
- Materials can be flagged as MFC-relevant and thus be available for MFC management (page 37).
- Weighing-specific attributes can be defined.
Weighing methods and tolerances are populated automatically from the material. Weighing material type and planned potency are preset with default values and will be populated into the respective list of materials item.
- Positions of MFC-relevant materials (process input) can be defined.
- For master recipes based on an ERP BOM, the ERP BOM items can be assigned as material parameters via **Setlist (SR3146.9.2.9)** UI element (page 124).

PROCESS INPUTS - SIGNATURE PRIVILEGE (SR3146.9.5.3)

Signature privileges are supported as process input parameters.

- They apply to phase-specific exceptions, explicit capability-related exceptions, and the phase completion.
- The correlation between signature privilege and exception-specific risk level must be defined manually.
- Per PharmaSuite installation, a default signature privilege can be configured which will be used if no signature privilege has been assigned to a phase for a specific risk level.

PROCESS INPUTS - SIGNATURE PRIVILEGE ATTRIBUTES (SR3146.9.5.3.1)

The signature privilege input parameters support the following attributes:

- Read-only:
 - Identifier

- Description (only in the **Universe (SR3146.9.2.4+)** UI element (page [122](#)))
- Signature type (Single/Double, only in the **Universe (SR3146.9.2.4+)** UI element (page [122](#)))
- Default reason (1st)
- Default reason (2nd)
- Category
- Editable:
 - Usage (High risk exception, High risk exception (comment mandatory), Medium risk exception, Medium risk exception (comment mandatory), Low risk exception, Low risk exception (comment mandatory), No risk exception, Phase completion, Phase completion (sequential))
 - Usage name
 - Reason (1st)
 - Maximum length is 80 characters.
 - Reason (2nd)
 - Maximum length is 80 characters.

PROCESS INPUTS - WORK CENTER (SR3146.9.5.4)

Work centers are supported as process input parameters on unit procedure level.

The operations and phases of the unit procedure have to be processed at one of the defined work centers.

PROCESS INPUTS - WORK CENTER ATTRIBUTES (SR3146.9.5.4.1)

The work center input parameters support the following attributes:

- Identifier
- Description
- Storage area

PROCESS INPUTS - STATION (SR3146.9.5.8)

Stations are supported as process input parameters on operation level.

The phases of the operation have to be processed at one of the defined stations.

PROCESS INPUTS - STATION ATTRIBUTES (SR3146.9.5.8.1)

The station input parameters support the following attributes:

- Identifier
- Description
- Work center

PROCESS INPUTS - EQUIPMENT REQUIREMENT (CLASS) (SR3146.9.5.5)

Equipment requirements of the **Class** type are supported as process input parameters.

PROCESS INPUTS - REQUIREMENT ATTRIBUTES (CLASS) (SR3146.9.5.5.1)

The equipment class input parameters support the following attributes:

- Read-only:
 - Identifier
 - Short description
 - Description
 - Level
 - Type (e.g. Class, Item)
- Editable:
 - BOE position
 - Usage name

PROCESS INPUTS - EQUIPMENT REQUIREMENT (PROPERTY TYPE) (SR3146.9.5.6)

Property types are supported as part of an already existing equipment requirement (e.g. equipment class).

PROCESS INPUTS - REQUIREMENT ATTRIBUTES (PROPERTY TYPE) (SR3146.9.5.6.1)

The (equipment) property type input parameters support the following attributes:

- Editable:
 - Identifier (unique)
 - Description
 - Rule (only editable for property types of the **Specification** usage type and of the **Runtime** usage type, not editable for property types of the **Automation** or **Historian** usage type, representation always shows the related property type identifier)

- Read-only:
 - Type (Dependent, Independent)

PROCESS INPUTS - EQUIPMENT REQUIREMENT (RULE) (SR3146.9.5.7)

Rules are supported as part of an already existing equipment requirement. Flexible rules apply to a single equipment entity and the parent entity of an equipment entity group whereas conditional rules apply to a single equipment entity and to all entities of an equipment entity group.

PROCESS INPUTS - REQUIREMENT ATTRIBUTES (RULE) (SR3146.9.5.7.1)

The rules of an equipment class input parameter supports the following attributes:

- Editable (with the **Expression editor (SR3146.9.9.4)** UI element (page [131](#))):
 - Identifier (unique)
 - Description
 - Rule
- Read-only:
 - Type (Flexible, Conditional rule (group-enabled))

PROCESS INPUTS - TRANSITION (SR3146.9.5.11)

Transitions are supported as process input parameters.

PROCESS INPUTS - TRANSITION ATTRIBUTES (SR3146.9.5.11.1)

The transition input parameters support the following attributes:

- Identifier
- Description
- Condition

Transition conditions can be edited in an **Expression editor (SR3146.9.9.4)** UI element (page [131](#)).

PROCESS INPUTS - CAPABILITIES ON OPERATION LEVEL (SR3146.9.5.9)

Capabilities are supported as process input parameters for operations.

PROCESS INPUTS - DETACHABLE CAPABILITY - OPERATION (SR3146.9.5.9.1)

The system provides the Detachable capability, which supports the following attributes:

- Read-only:
 - Identifier

- Action
- Editable:
 - Exception enabled
 - Risk assessment
 - Exception text

This capability allows an operator to detach a running operations from its work center and station (see **Detaching an Operation (SR1089.6.1+)** in "Functional Requirement Specification Execution Framework" [A2] (page [147](#))).

This capability cannot be used in Dispense operations.

PROCESS INPUTS - EVENT-TRIGGERED CAPABILITY (SR3146.9.5.9.2)

The system provides the Event-triggered capability, which supports the following attributes:

- Read-only:
 - Identifier
 - Action (New run, Cancel run, Cancel all runs, Remove template)
- Editable (for each action):
 - Exception enabled
 - Risk assessment
 - Exception text

This capability cannot be used in Dispense operations.

PROCESS INPUTS - ESCALATION-ENABLED CAPABILITY (SR3146.9.5.9.3)

The system provides the Escalation-enabled capability, which supports the following attributes:

- Read-only:
 - Identifier
 - System action (Run due reminder, Run overdue, Run expired)
- Editable (for each system action):
 - Duration before action (since trigger)
 - Automatic exception enabled
 - Risk assessment

■ Exception text

This capability is only supported for event-triggered operations.
This capability cannot be used in Dispense operations.

PROCESS INPUTS - TRIGGER-ENABLED CAPABILITY (SR3146.9.5.9.4)

The system provides the Trigger-enabled capability, which supports the following attributes:

- Read-only:
 - Identifier
- Editable (for each new line item):
 - Trigger (unique)
 - Description
 - Trigger phase

The system allows to add multiple line items for multiple trigger references.

This capability is only supported for event-triggered operations.
This capability cannot be used in Dispense operations.

PROCESS INPUTS - SERVER-RUN CAPABILITY (SR3146.9.5.9.5)

The system provides the Server-run capability, which supports the following attributes:

- Read-only:
 - Identifier

This capability cannot be combined with other capabilities.

This capability cannot be used in Dispense operations.

PROCESS INPUTS - AUTO-STARTABLE CAPABILITY (SR3146.9.5.9.6)

The system provides the Auto-startable capability, which supports the following attributes:

- Read-only:
 - Identifier

This capability cannot be used in Dispense operations.

PROCESS INPUTS - CAPABILITIES ON UNIT PROCEDURE LEVEL (SR3146.9.5.10)

Capabilities are supported as process input parameters for unit procedures.

PROCESS INPUTS - PAUSE-ENABLED CAPABILITY (SR3146.9.5.10.1)

The system provides the Pause-enabled capability, which supports the following attributes:

- Read-only:
 - Identifier
 - Action
- Editable:
 - Exception enabled
 - Risk assessment
 - Exception text

PROCESS INPUTS - DETACHABLE CAPABILITY - UNIT PROCEDURE (SR3146.9.5.10.2)

The system provides the Detachable capability, which supports the following attributes:

- Read-only:
 - Identifier
 - Action
- Editable:
 - Exception enabled
 - Risk assessment
 - Exception text

This capability allows an operator to detach a running unit procedure from its work center and station (see **Detaching a Unit Procedure (SR1089.8.4+)** in "Functional Requirement Specification Execution Framework" [A2] (page 147)).

This capability cannot be used in Dispense unit procedures.

PROCESS INPUTS - ABORT-AND-REACTIVATE-ENABLED CAPABILITY (SR3146.9.5.10.3)

The system provides the Abort-and reactivate-enabled capability, which supports the following attributes:

- Read-only:
 - Identifier

- Action
- Editable:
 - Exception enabled
 - Risk assessment
 - Exception text

This capability allows an operator to abort and reactivate a running unit procedure (see **Reactivating a Unit Procedure (SR1089.8.5)** in "Functional Requirement Specification Execution Framework" [A2] (page [147](#))).

This capability cannot be used in Dispense unit procedures and master workflows.

Process Outputs (SR3146.9.6+)

According to S88, process outputs are supported as part of the formula.

The attributes of process output parameters are maintained in the **Parameter Panel (SR3146.9.2.7)** UI element (page [126](#)).

PROCESS OUTPUTS - MATERIAL (SR3146.9.6.1)

Materials in the sense of product are supported as process output parameters.

The minimum and maximum numbers of process outputs are taken over from the phase building block.

PROCESS OUTPUTS - MATERIAL ATTRIBUTES (SR3146.9.6.1.1)

➤ For recent changes, see revision history (page [149](#)).

The material output parameters support the following attributes:

- Basic data
 - Lock
 - Material identifier (read-only)
 - [Material] Short description (read-only)
 - Planned quantity mode (None, As defined)
 - Planned quantity (with UoM)
 - Usage type (e.g. Input, Output) (read-only)
 - Usage name
 - MFC-relevant (Yes/No), default = No

- Position
- Hierarchy level (read-only)
- Path (read-only)
- Weighing:
 - Default weighing method
 - Allowed weighing methods
 - Tolerances (lower and upper, absolute and relative)
 - Weighing material type
(not relevant to material outputs and ignored during execution)
 - Planned potency
(not relevant to material outputs and ignored during execution)
 - Fixed quantity (Yes/No)
(not relevant to material outputs and ignored during execution)
 - Target weight
 - Target weight tolerances (lower and upper, absolute and relative)
- Packaging Levels (editable with the **Packaging Level Data (SR3146.9.2.17.3)** editor (page 140)):
 - L0 (Meaning, Contained number, Inventory level, Hide during execution)
(the Contained number attribute is not relevant to the L0 level)
 - ...
 - L5 (Meaning, Contained number, Inventory level, Hide during execution)

Materials can be flagged as MFC-relevant and thus be available for MFC management (page 37).

Process Parameters (SR3146.9.7+)

According to S88, process parameters are supported as part of the formula.

- Process parameters and their attributes are taken over from the phase building block and parameter class, respectively.
- For phases located within the same operation and across operations, the output of a phase can be used as input of a parameter attribute of another phase.
- Inputs can be edited by an **Expression editor (SR3146.9.9.4)** UI element (page 131).

- The phase that provides the output must have been processed before the process step that uses the output as input. Phases and operations may or may not be passed (e.g. in simultaneous branches). If they are not passed, the output is Null. This needs to be reflected in the referencing expression.

PROCESS PARAMETERS - CLASSIFICATION (SR3146.9.7.1)

For process parameters, the following classifications are available:

- Critical quality attribute
- Critical process parameter
- Key quality attribute
- Key process parameter
- General quality attribute
- General process parameter

DYNAMIC CREATION OF PROCESS PARAMETERS (SR3146.9.7.4)

The **Parameter Panel (SR3146.9.2.7)** UI element (page [126](#)) for process parameters supports the dynamic creation of process parameters for specific phases during recipe and workflow authoring. The availability of a parameter is controlled by the current phase.

When a new parameter bundle is created, the recipe or workflow author defines the bundle identifier.

Process parameters that have been created in this way can also be deleted within the Parameter Panel.

Bundle Process Parameters (SR3146.9.7.4.1)

For the master process parameter of a bundle, its internal identifier is populated from the bundle identifier.

For all other process parameters of the bundle, their internal identifier is a concatenation of the bundle identifier and the process parameter name.

Bundle Output Variables (SR3146.9.7.4.2)

For all output variables of the same bundle, the output variable identifier is a concatenation of the bundle identifier and the output variable name.

The naming convention applies to PharmaSuite product phases.

Locking of Parameters (SR3146.9.13+)

Locking and unlocking of parameters is supported.

LOCKED PARAMETERS (SR3146.9.13.1)

Parameters can be locked and unlocked as long as the related master recipe or master workflow is not in the **Valid (SR3146.6.4)** status (page 43) or the custom BB is not in the **Approved (SR3146.10.1)** status (page 75).

The attributes of a locked parameter are read-only for this specific occurrence. They can again be updated when the parameter is unlocked.

FROZEN PARAMETERS (SR3146.9.13.2)

A frozen parameter of a master recipe, master workflow, or custom BB is permanently locked if the related source custom BB was in the **Approved (SR3146.10.1)** status (page 75) at the time when it was copied for usage. Thus, the related attributes cannot be updated anymore.

Locking of Transitions (SR3146.9.14+)

Locking and unlocking of transitions is supported.

LOCKED TRANSITIONS (SR3146.9.14.1)

Transitions can be locked and unlocked as long as the related master recipe or master workflow is not in the **Verification** or **Valid (SR3146.6.4)** status (page 43) or the custom BB is not in the **Approved (SR3146.10.1)** status (page 75).

The attributes of a locked transition are read-only for this specific occurrence. They can again be updated when the transition is unlocked.

FROZEN TRANSITIONS (SR3146.9.14.2)

A frozen transition of a master recipe, master workflow, or custom BB is permanently locked if the related source custom BB was in the **Approved (SR3146.10.1)** status (page 75) at the time when it was copied for usage. Thus, the related attributes cannot be updated anymore.

Workbench Object Operations (SR3146.9.10+)

The Designer workbench supports specific object operations for:

- master recipe (SR3146.9.10.1) objects (page 28),
- master workflows (SR3146.9.10.9) objects (page 29),
- element (SR3146.9.10.6) objects (page 31),
- custom building block (SR3146.9.10.7) objects (page 74), and
- change request (SR3146.12+) objects (page 83).

Operations for links and transitions (page 32) are also supported.

Workbench Object Operations - Master Recipe (SR3146.9.10.1)

- For recent changes, see revision history (page 149).

The following operations are supported for master recipes:

- Create new (product material-related) master recipe (page 29)
- Open existing (see **Open (SR3146.9.2.20+)** dialog (page 122))
- View locked master recipe (page 29)
- Save (current or all) master recipe (page 29)
- Copy (Save as...)

For master recipes that are protected by an access privilege: The access privilege of the original master recipe is passed on to the new master recipe. In the default configuration, the access privilege cannot be changed.

- Rename
- Delete
- Create building block from selected object
For master recipes that are protected by an access privilege: The access privilege of the original master recipe is passed on to the new building block. In the default configuration, the access privilege cannot be changed.
- Insert existing procedure BB
- Replace Material (Product) (SR3146.9.10.11) (page 36)
- Status-related operations (see Version Control (SR3146.4.5) (page 141))
- Report
 - Master Recipe Report - Batch (SR3146.11+) (page 95)
 - Master Recipe Report - Device (SR3146.14+) (page 103)
- Statistics (SR3146.9.2.18) (page 140)

- Close (current or all)

Master recipes can either be saved under a new identifier or under the same identifier and a new version number.

WORKBENCH OBJECT OPERATIONS - NEW MASTER RECIPE DIALOG (SR3146.9.10.1.2)

- For recent changes, see revision history (page 149).

When a new master recipe is created, first, the recipe author defines the product material of the master recipe. The list of available materials shall be restricted by the access privileges of the logged-in user. The material identifier is populated as master recipe identifier. The master recipe identifier is editable.

The system shall allow to select an access privilege for the protection of the master recipe from unauthorized access for users with access privileges that allow to maintain protected recipes. The list of available access privileges shall be restricted to the access privileges of the logged-in user.

If an access privilege is defined for the material, this access privilege is passed on to the new master recipe. In the default configuration, the access privilege cannot be changed.

For master recipes based on an ERP BOM, the ERP BOM is selected from the **library for materials (ERP BOM) (SR3146.9.2.4.2)** UI element (page 123). The **material process input (SR3146.9.5.1)** data (page 16) defined within Recipe Designer must match the materials, positions, and quantities set down in the ERP BOM. The ERP BOM items are available via **Setlist (SR3146.9.2.9.1)** UI element (page 125).

WORKBENCH OBJECT OPERATIONS - PREPARE (DEEP SAVE) (SR3146.9.10.1.3)

Master recipes that have no error or warning message left in the **Messages (SR3146.9.2.16+)** window (page 128) can be prepared for **status changes (SR3146.6.4)** by means of the version control (page 43).

WORKBENCH OBJECT OPERATIONS - VIEW LOCKED MASTER RECIPE (SR3146.9.10.1.4)

If a master recipe is locked by another user, the system allows to view the master recipe without the possibility to change it.

The master recipe can be copied (save as). Operations that do not change a master recipe are also available (e.g. status history, print report, graph pagination, statistics).

Workbench Object Operations - Master Workflow (SR3146.9.10.9)

- For recent changes, see revision history (page 149).

The following operations are supported for master workflows:

- Create new master workflow (page 30)

- Open existing (see **Open (SR3146.9.2.20+)** dialog (page 122))
- View locked master workflow (page 31)
- Save (current or all) master workflow (page 30)
- Copy (Save as...)
For master workflows that are protected by an access privilege: The access privilege of the original master workflow is passed on the new master workflow. In the default configuration, the access privilege cannot be changed.
- Rename
- Delete
- Create building block from selected object
For master workflows that are protected by an access privilege: The access privilege of the original master workflow is passed on to the new building block. In the default configuration, the access privilege cannot be changed.
- Insert existing procedure BB
- Status-related operations (see Version Control (SR3146.4.5) (page 141))
- Master Workflow Report (SR3146.13+) (page 111)
- Statistics (SR3146.9.2.18) (page 140)
- Close (current or all)

Master workflows can either be saved under a new identifier or under the same identifier and a new version number.

WORKBENCH OBJECT OPERATIONS - NEW MASTER WORKFLOW DIALOG (SR3146.9.10.9.1)

- For recent changes, see revision history (page 149).

When a new master workflow is created, first, the workflow author defines the master workflow identifier.

The system shall allow to select an access privilege for the protection of the master workflow from unauthorized access for users with access privileges that allow to maintain protected workflows. The list of available access privileges shall be restricted to the access privileges of the logged-in user.

WORKBENCH OBJECT OPERATIONS - PREPARE (DEEP SAVE) (SR3146.9.10.9.2)

Master workflows that have no error or warning message left in the **Messages (SR3146.9.2.16+)** window (page 128) can be prepared for **status changes (SR3146.6.4)** by means of the version control (page 43).

WORKBENCH OBJECT OPERATIONS - VIEW LOCKED MASTER WORKFLOW (SR3146.9.10.9.3)

If a master workflow is locked by another user, the system allows to view the master workflow without the possibility to change it.

The master workflow can be copied (save as). Operations that do not change a master workflow are also available (e.g. status history, print report, graph pagination, statistics).

Workbench Object Operations - Element (SR3146.9.10.6)

➤ For recent changes, see revision history (page 149).

The following operations are supported for elements that are derived from building blocks (e.g. recipe element within a recipe, workflow element within a workflow, BB element within a custom BB). They apply to elements on all levels (procedure, unit procedure, operation, phase).

- Create new (e.g. dummy phase)
- Insert element from an existing BB via the **Setlist (SR3146.9.2.9)** UI element (page 124)
- Create BB from selected element
For elements that are part of a master recipe, master workflow, or custom building block that is that is protected by an access privilege: The access privilege of the header object is passed on to the new building block. In the default configuration, the access privilege cannot be changed.
- Replace (overrides parameters) via the **Setlist (SR3146.9.2.9)** UI element (page 124)
- Smart replace via the **Setlist (SR3146.9.2.9)** UI element (page 124) (see **Workbench Object Operations - Smart Replace (SR3146.9.10.10)** operation (page 32))
- Rename
- Delete selected elements

- Connect elements manually. However, Recipe and Workflow Designer inserts possible or required transitions (page 14) and required selection/simultaneous branches (page 14) automatically.
- Insert branch with exclusive alternative steps (XOR). A selection branch (page 14) and required transitions (page 14) connect the elements.
- Insert loop. A loop, required selection branches, and a transition connect the elements.

- Insert branch with parallel steps (AND). A simultaneous branch (page 14) connects the elements.
- If a BB is created from a selected element, the new building block is available for selection from the **Universe (SR3146.9.2.4+)** UI element (page 122) and the **Setlist (SR3146.9.2.9)** UI element (page 124) and is opened in a new tab. Available revisions of building blocks with the same identifier are listed for information.
- Insert sequence. A link (page 14) and a transition (page 14) connect the elements.
- Select all elements and transitions along with or without links.
- Find elements by searching for a specific string.

Workbench Object Operations - Link and Transition

The following operations are supported for links (page 14) and transitions (page 14):

- Delete links and transitions selected in a graph.
- Select all links and transitions along with all elements in the current Graph Window.
- Select all transitions along with all elements in the current Graph Window.
- Select all links in the current Graph Window.

Workbench Object Operations - Smart Replace (SR3146.9.10.10)

With the **Smart replace** operation, custom building blocks (cBBs) and elements of recipes, workflows, and cBBs are replaced by another cBB from the Setlist. The operation differs from the **Replace** operation in its parameter and transition handling. It merges the parameter data and transitions of the replacing (new) building block into the parameter data and transitions of the building block to be replaced (old building block).

If the new cBB is in the **Draft** or **Approved** status (page 75), the parameters are merged depending on their type, **locking state** (SR3146.9.13.1) (page 27), and existence; transitions are merged depending on their **locking state** (SR3146.9.14.1) (page 27) and existence.

- Merge strategy for **process parameter** (SR3146.9.7+) data (page 25) and their attributes

New cBB	Old cBB	Resulting building block after smart replace
Parameter is permanently locked (Frozen).	Parameter is permanently locked (Frozen).	The new parameter is taken over , it is still permanently locked (Frozen). The parameter of the old cBB is ignored.
	Parameter is locked for this specific occurrence on the BB, but could be unlocked.	
	Parameter is not locked .	
Parameter is locked for this specific occurrence on the BB, but could be unlocked.	Parameter is permanently locked (Frozen).	The new parameter is taken over , it is locked for this specific occurrence on the BB, but could be unlocked. The parameter of the old cBB is ignored. This combination is not relevant to change request (SR3146.12+) objects (page 83).
	Parameter is locked for this specific occurrence on the BB, but could be unlocked.	
	Parameter is not locked .	
Parameter is not locked .	Parameter is permanently locked (Frozen).	The previous parameter is taken over , but is not locked . The parameter of the new cBB is ignored.
	Parameter is locked for this specific occurrence on the BB, but could be unlocked.	The previous parameter is taken over , it is still locked for this specific occurrence on the BB, but could be unlocked. The parameter of the new cBB is ignored.
	Parameter is not locked .	The previous parameter is taken over , it is still not locked . The parameter of the new cBB is ignored.
Parameter in any locking state.	No parameter present.	The new parameter is taken over in its locking state.
No parameter present.	Parameter in any locking state.	The previous parameter does no longer exist.

- Merge strategy for the following parameter-specific attributes:
 - The attributes are unique keys.
 - Usage attribute of **Signature privilege** (SR3146.9.5.3) parameters (page 17)
 - **BOE position** attribute of **Equipment requirement (class)** (SR3146.9.5.5) parameters (page 19)

- **Work center identifier** attribute of **Work center (SR3146.9.5.4)** parameters (page 18)
- **Station identifier** attribute of **Station (SR3146.9.5.8)** parameters (page 18)
- **Capability identifier** attribute of **Capability (SR3146.9.5.9)** parameters (page 20)

New cBB	Old cBB	Resulting building block after smart replace
Parameter is permanently locked (Frozen).	Parameter is permanently locked (Frozen).	The new parameter is taken over , it is still permanently locked (Frozen). The parameter of the old cBB is ignored.
	Parameter is locked for this specific occurrence on the BB, but could be unlocked.	
	Parameter is not locked .	
	No parameter present.	The new parameter is taken over , it is still permanently locked (Frozen).
Parameter is locked for this specific occurrence on the BB, but could be unlocked.	Parameter is permanently locked (Frozen).	The new parameter is taken over , it is still locked for this specific occurrence on the BB, but could be unlocked. The parameter of the old cBB is ignored. This combination is not relevant to change request (SR3146.12+) objects (page 83).
	Parameter is locked for this specific occurrence on the BB, but could be unlocked.	
	Parameter is not locked .	
	No parameter present.	The new parameter is taken over , it is still locked for this specific occurrence on the BB, but could be unlocked. This combination is not relevant to change request (SR3146.12+) objects (page 83).
Parameter is not locked .	Parameter is permanently locked (Frozen).	The previous parameter is taken over , but is not locked . The parameter of the new cBB is ignored.
	Parameter is locked for this specific occurrence on the BB, but could be unlocked.	The previous parameter is taken over , it is locked for this specific occurrence on the BB, but could be unlocked. The parameter of the new cBB is ignored.
	Parameter is not locked .	The previous parameter is taken over , it is still not locked . The parameter of the new cBB is ignored.
	No parameter present.	There is still no parameter present. The parameter of the new cBB is ignored.

New cBB	Old cBB	Resulting building block after smart replace
No parameter present.	Parameter in any locking state.	The previous parameter is still locked (if locked before) and not locked in all other cases.

- Material parameters (**process input (SR3146.9.5.1)** data (page 16), **process output (SR3146.9.6.1)** data (page 24))
They are completely ignored during the **Smart replace** operation.
- Merge strategy for transitions (page 20) and their attributes
 - The identifier is the unique key within a procedural element.

New cBB	Old cBB	Resulting building block after smart replace
Transition is permanently locked (Frozen).	Transition is permanently locked (Frozen).	The new transition is taken over , it is still permanently locked (Frozen). The transition of the old cBB is ignored.
	Transition is locked for this specific occurrence on the BB, but could be unlocked.	
	Transition is not locked .	
Transition is locked for this specific occurrence on the BB, but could be unlocked.	Transition is permanently locked (Frozen).	The new transition is taken over , it is locked for this specific occurrence on the BB, but could be unlocked. The transition of the old cBB is ignored. This combination is not relevant to change request (SR3146.12+) objects (page 83).
	Transition is locked for this specific occurrence on the BB, but could be unlocked.	
	Transition is not locked .	
Transition is not locked .	Transition is permanently locked (Frozen).	The previous transition is taken over , but is not locked . The transition of the new cBB is ignored.
	Transition is locked for this specific occurrence on the BB, but could be unlocked.	The previous transition is taken over , it is still locked for this specific occurrence on the BB, but could be unlocked. The transition of the new cBB is ignored.
	Transition is not locked .	The previous transition is taken over , it is still not locked . The transition of the new cBB is ignored.
Transition in any locking state.	No transition present.	The new transition is taken over in its locking state.
No transition present.	Transition in any locking state.	The previous transition does no longer exist.

Applies to all parameters and transitions:

- If a building block on a higher level (e.g. operation level) is replaced using the **Smart Replace** operation and the new cBB contains a new recipe element, workflow element, BB element, or transition that did not exist in the previous version, this new element or transition is taken over including all of its parameter and transition configurations.

Applies to material parameters:

- If a building block is replaced using the regular **Replace** operation or is added to a **procedural structure (SR3146.9.4)** hierarchy (page 13), all of its already assigned material parameters are still taken over. They are only ignored in case the **Smart replace** operation is used.

Workbench Object Operations - Replace Material (Product) (SR3146.9.10.11)

- For recent changes, see revision history (page 149).
- Does not apply to Workflow Designer.

With the **Replace material (product)** operation, the material (product) of an existing master recipe, including its related ERP BOM (if applicable), can be replaced.

The system shall only allow to select ERP BOMs and materials whose access privilege matches the access privilege of the logged-in user.

If an access privilege is defined for the new ERP BOM or material that differs from the access privilege of the master recipe, the new access privilege is passed on to the master recipe.

Before the new material can be selected, the system requests a confirmation by the recipe author since the **Replace material (product)** operation can lead to a significant number of changes in the master recipe data.

Along with replacing the material (product), the system automatically updates the following data:

- Material (product) of the master recipe's header attributes
- Material of the final MFC output
- Material of the MFC-relevant material input and output parameters (that used the previous material)
- Packaging levels (Meaning, Contained number) of the material parameters (L0, ..., L5) (that used the previous material)
- Only if the new material (product) is related to an ERP BOM:
 - ERP BOM items in the Setlist
 - Planned quantity of the master recipe

- Packaging levels of the master recipe (L6, ..., L9)
- Material inputs by BOM position (material identifier, material short description, planned quantity, and fixed quantity)
- In case linked MFC inputs/outputs, according to the MFC definition, refer to the same material:
The system also updates the material identifier and planned quantity of the related material outputs (same unit procedure, e.g. Dispense, planned quantity only for the first linked output) and the material identifier of the subsequent material inputs (subsequent unit procedure, e.g. Mixing).

In case there are remaining material-related inconsistencies that cannot be resolved automatically, the system displays the corresponding messages in the **Messages (SR3146.9.2.16+)** window (page 128). For details, see **Validator Consistency Checks (SR3146.9.4.3+)**, Checks that Apply to Master Recipes Based on an ERP BOM (page 54) and Checks that Apply to Material Parameters of a Phase (page 56).

Details of the performed changes, are visible in the **Comparison (SR3146.9.2.19+)** window (page 130), based on a defined comparison baseline.

Material Flow Control (MFC) (SR3146.9.8)

A graphical tool allows to model the material flow of a master recipe or master workflow through a procedure within the Recipe and Workflow Designer workbench. The MFC-relevant material parameters (SR3146.9.5.1 (page 16), SR3146.9.6.1 (page 24)) provide the basis for the material flow.

Graphical View (SR3146.9.8.1)

- For recent changes, see revision history (page 149).

The graphical representation displays the following components:

- Unit procedures (of the procedure)
- MFC-relevant material inputs
- MFC-relevant material outputs

- Unit procedures hold all MFC-relevant material nodes defined for their respective phases.
Before a merge, there are no connections available between unit procedures.
- Direct MFC inputs are derived from the material input parameters of the procedure's phases.
A direct MFC input is displayed as triangle with an outgoing connector.
The number it shows corresponds to the numeric suffix of the MFC item identifier visible in the table view of the MFC data.

- Does not apply to Recipe Designer - Device and Workflow Designer.
Direct MFC outputs are derived from the material output parameters of the procedure's phases.
A direct MFC output is displayed as triangle with an incoming connector.
The number it shows corresponds to the numeric suffix of the MFC item identifier visible in the table view of the MFC data.
- Does not apply to Recipe Designer - Device and Workflow Designer.
Switch nodes are created when direct MFC inputs and outputs are merged. A direct MFC input is replaced by an input switch and a direct MFC output is replaced by an output switch.
A switch node is displayed as a small circle between unit procedures.
- A double-click of a direct MFC input, direct MFC output, input switch, output switch, confluence, or connector shows the related material parameter(s) within the **Parameter Panel (SR3146.9.2.7)** UI element (page 126).
- The number an MFC input or output shows, corresponds to the numeric suffix of the MFC item identifier, which the system generates and assigns automatically.
Before a merge, there are no connections available between inputs and outputs.
- Does not apply to Recipe Designer - Device and Workflow Designer.
A confluence node is displayed as circle and represents the point in the material flow when an input is processed and thus converted into an output or merged with other inputs to form an output. A confluence can receive one or more inputs, but can only issue one output.
The confluence node of an MFC input that has been merged with the final output is displayed as a triangle with an incoming connector and no number.
- Does not apply to Recipe Designer - Device and Workflow Designer.
Connectors represent the material flows between the direct MFC input, direct MFC output, confluence nodes, and switch nodes. Connectors between switch nodes are dashed lines if they represent alternative transfers.

Table View (SR3146.9.8.1.1)

The table view of the MFC data provides following features:

- Representation of MFC-related attributes, plus position and planned quantity for ingoing, outgoing, and transfer items
- If the planned quantity mode is **None** or **As produced**, the mode is displayed instead of any planned quantity.
- Synchronization of highlighted MFC items within the table and the graph
- A double-click shows the related material parameter(s) within the **Parameter Panel (SR3146.9.2.7)** UI element (page 126).

Does not apply to Recipe Designer - Device and Workflow Designer.
The number of a transfer item consists of the numbers of the output and the input items from which it was merged.

Tooltips (SR3146.9.8.2)

- For recent changes, see revision history (page [149](#)).

Tooltips are available for

- direct MFC inputs,
- direct MFC outputs,
Does not apply to Recipe Designer - Device and Workflow Designer.
- confluences, switches, and connectors (transfers).
Does not apply to Recipe Designer - Device and Workflow Designer.

They present the following attributes:

- Material identifier
- [Material] Short description,
- Position (if available)
- Planned quantity (if available)
- If the planned quantity mode is **None** or **As produced**, the mode is displayed instead of any planned quantity.
- Path information (not for connectors)

Automatic Merge - Default MFC Definition (SR3146.9.8.3)

- For recent changes, see revision history (page [149](#)).
- Does not apply to Recipe Designer - Device and Workflow Designer.

The recipe author can trigger the automatic merge of direct MFC inputs, direct MFC outputs, input switches, and output switches at any time during the MFC definition process.

Already merged inputs and outputs are not affected.

Nodes are merged when the system can unambiguously determine their merge targets. When the system identifies two or more merge targets, it does not perform a merge.

DEFAULT MFC DEFINITION WITHIN A UNIT PROCEDURE (SR3146.9.8.4)

- Does not apply to Recipe Designer - Device and Workflow Designer.

The following rules apply to the automatic merge of MFC inputs and outputs within a unit procedure:

- The inputs and outputs of a Dispense unit procedure are merged 1:1, if they have the same material assigned.
- The inputs and outputs of a non-Dispense unit procedure are merged n:1, if only one output is defined. Different materials can be assigned to the inputs and the output.
- In any other case, the inputs and outputs are not merged.
Already merged inputs and outputs are not affected.

DEFAULT MFC DEFINITION BETWEEN UNIT PROCEDURES (SR3146.9.8.5)

- For recent changes, see revision history (page [149](#)).
- Does not apply to Recipe Designer - Device and Workflow Designer.

The following rules apply to the automatic merge of MFC inputs and outputs between unit procedures:

- The outputs and inputs are merged 1:1, if there are unique output-input relationships available based on the material assignments and the SFC-related modeling of unit procedures.
In case of a selection branch (XOR), the merge results in alternative transfers.
- In any other case, the inputs and outputs are not merged.
Already merged inputs and outputs are not affected.

Manual Merge

MANUAL MERGE WITHIN A UNIT PROCEDURE (SR3146.9.8.6)

- Does not apply to Recipe Designer - Device and Workflow Designer.

The following rules apply to the manual merge of MFC inputs and outputs within a unit procedure:

- Outputs cannot be merged.
- For Dispense unit procedures:
 - Inputs cannot be merged.
 - Inputs and outputs can be merged 1:1, even if they have different materials assigned.
 - When one input (output) is selected, all outputs (inputs) that have the same material assigned are highlighted.

- For non-Dispense unit procedures:
 - Multiple inputs can be merged to one output n:1, even if they have different materials assigned.
 - When one input is selected, all inputs and outputs with which it can potentially be merged are highlighted.
 - When one output is selected, all inputs with which it can potentially be merged are highlighted.

The system can be enabled to support merging multiple inputs to one output not only for non-Dispense, but also for Dispense unit procedures by means of a configuration key (see Technical Manual "Configuration and Extension").

MANUAL MERGE BETWEEN UNIT PROCEDURES (SR3146.9.8.7)

- For recent changes, see revision history (page 149).
- Does not apply to Recipe Designer - Device and Workflow Designer.

The following rules apply to the manual merge of direct MFC inputs, direct MFC outputs, input switches, and output switches between unit procedures:

- Inputs and outputs can only be merged
 - if they have the same material assigned and
 - if the input belongs to a succeeding unit procedure according to the SFC model or if unit procedures are modeled in parallel and
- In case of simultaneous branches (no XOR), inputs and outputs can only be merged 1:1.
- In case of selection branches (XOR), one input can be merged into multiple alternative outputs and one output can be merged into multiple alternative inputs. In both cases alternative transfers are created.
- When an input (output) is selected, all outputs (inputs) that have the same material assigned and are suitable according to the MFC structure are highlighted.

MANUAL MERGE OF MFC INPUTS WITH THE FINAL OUTPUT (SR3146.9.8.13)

- Does not apply to Recipe Designer - Device and Workflow Designer.

The following rules apply to the manual merge of MFC inputs with the final output:

- Merge of inputs with the final output is not allowed for Dispense unit procedures (see **Validator - MFC – Final output for Dispense unit procedures (SR3146.9.4.3.118)** message (page 61)).
- Any confluence of a direct MFC input can be merged with the final output of the master recipe.

- Multi-selection is supported.

Split (SR3146.9.8.8)

- For recent changes, see revision history (page 149).
- Does not apply to Recipe Designer - Device and Workflow Designer.

Inputs and outputs that are presented in a merged status can be split.

Multi-selection of MFC inputs merged with the final output is supported.

- Split connector: dissolves the connection between two nodes.
- Split direct MFC input or output: removes input or output from its current confluence and recreates its initial (single) confluence node.
- Split confluence: dissolves the entire confluence and recreates all initial (single) confluence nodes of its inputs and its outputs.
- Split output/input switch: dissolves all its outgoing/incoming connectors and replaces the switch with its initial output/input.

Reset (SR3146.9.8.9)

- Does not apply to Recipe Designer - Device and Workflow Designer.

The MFC definition can be reset to its initial, unmerged state by dissolving all merges.

Automated Update of MFC Representation (SR3146.9.8.10)

MFC-relevant changes of the master recipe and master workflow result in an automatic update of the graph and the tabular form (e.g. MFC-related material parameters, recipe structure, workflow structure).

MFC data and its representation that is not affected by the change of the master recipe or master workflow remains unchanged.

MFC Definition Linked to Procedure Level (SR3146.9.8.12)

The MFC definition of a master recipe or master workflow is linked to the procedure level, so that MFC definitions can also be modeled as part of a procedure-level custom building block and can be reused along with this custom building block.

Information Flow (SR3146.9.9+)

The information flow supports references of parameters across phases, operations, unit procedures, and supports transitions. The references are maintained in the **Expression editor (SR3146.9.9.4)** UI element (page 131).

Output Mechanism (SR3146.9.9.1)

The output mechanism allows phases to have specific output variables defined.

Referencing (SR3146.9.9.2)

Upon the definition of parameter attributes and transitions, output variables can be referenced.

Only those output variables are available for referencing whose phases can be potentially completed before the execution of the current phase is started, according to the SFC model.

References are supported within an operation and across operations and unit procedures.

NO DISTRIBUTED EXECUTION (SR3146.9.9.2.1)

Output variables that belong to operations marked for Distributed Execution cannot be referenced.

Validator (SR3146.9.9.3)

The automatic validation of references of output variables covers:

- Broken links
- Wrong syntax of expressions/calculations
- Semantic reference errors: referenced phase is not completed when referencing phase becomes active or referenced phase belongs to Dispense operation.

Standard Statuses for Version Control (SR3146.6.4)

The **version control (SR3146.4.5)** function (page 141) supports the following statuses:

- **Initial:** Item is new and has not been saved yet.
- **Edit:** Item is saved and can be modified.
- **Verification:** Item is ready for review. If updates are required, a status change back to **Edit** has to be performed.
- **Scheduled:** Item is valid but its period of validity is in the future related to the current global time.
This status allows fast filtering for specific items and the creation of a consistent documentation.
- **Valid:** Item is valid. All of its attributes that are under version control cannot be changed.
- **Archived:** Item is no longer valid. It can no longer be used.
- **Obsolete:** Item was never valid and should not be used any more.

The following important rules apply:

- At a specific time, only one version of an item can be **Valid**.
- To extend the period of validity, a status change from **Valid** to **Valid** must be performed.

Default Version Graph for Master Recipes (SR3146.6.4.1)

The **version control (SR3146.4.5)** function (page 141) supports the following default transitions for master recipes:

Transition (ID - From » To)	Signature (Access privilege, according to FSM)	Additional information
1.1 - Initial » Edit	---	Automatic transition.
1.2 - Verification » Edit	Vers_Trans_Verification-Edit	Default user groups: Recipe Author, Qualified Person, Supervisor (Shop Floor)
2.1 - Edit » Verification	Vers_Trans_Edit-Verification	Default user groups: Recipe Author, Qualified Person, Supervisor (Shop Floor)
3.1 - Verification » Valid	Vers_Trans_Verification-Valid	Default user group: Qualified Person
3.2 - Valid » Valid	Vers_Trans_Valid-Valid	E.g. if period of validity has to be changed. Default user group: Qualified Person
3.3 - Scheduled » Valid	Vers_Trans_Scheduled-Valid	Automatic transition when the current date reaches the start date of the period of validity. Default user group: Qualified Person
4.1 - Valid » Archived	Vers_Trans_Valid-Archived	Automatic transition when the end date of the period of validity is in the past. Automatic transition when the status of another version of a master recipe is set to Valid (except if there is a version in the Scheduled status). Manual transition when a status transition to Archived is executed by the recipe author. Default user group: Qualified Person
5.1 - Verification » Scheduled	Vers_Trans_Verification-Scheduled	If the scheduled period overlaps with a period of validity, the end date of the period of validity is moved to a date directly before the

Transition (ID - From » To)	Signature (Access privilege, according to FSM)	Additional information
		start date of the scheduled period. Default user group: Qualified Person
5.2 - Scheduled » Scheduled	Vers_Trans_Scheduled-Scheduled	E.g. if period of validity has to be changed. Default user group: Qualified Person
6.1 - Edit » Obsolete	Vers_Trans_Edit-Obsolete	Default user groups: Recipe Author, Qualified Person, Supervisor (Shop Floor)
6.2 - Verification » Obsolete	Vers_Trans_Verification-Obsolete	Default user groups: Recipe Author, Qualified Person, Supervisor (Shop Floor)
6.3 - Scheduled » Obsolete	Vers_Trans_Scheduled-Obsolete	Default user group: Qualified Person

The diagram illustrates the **MESMasterRecipeVersionGraph** flexible state model.

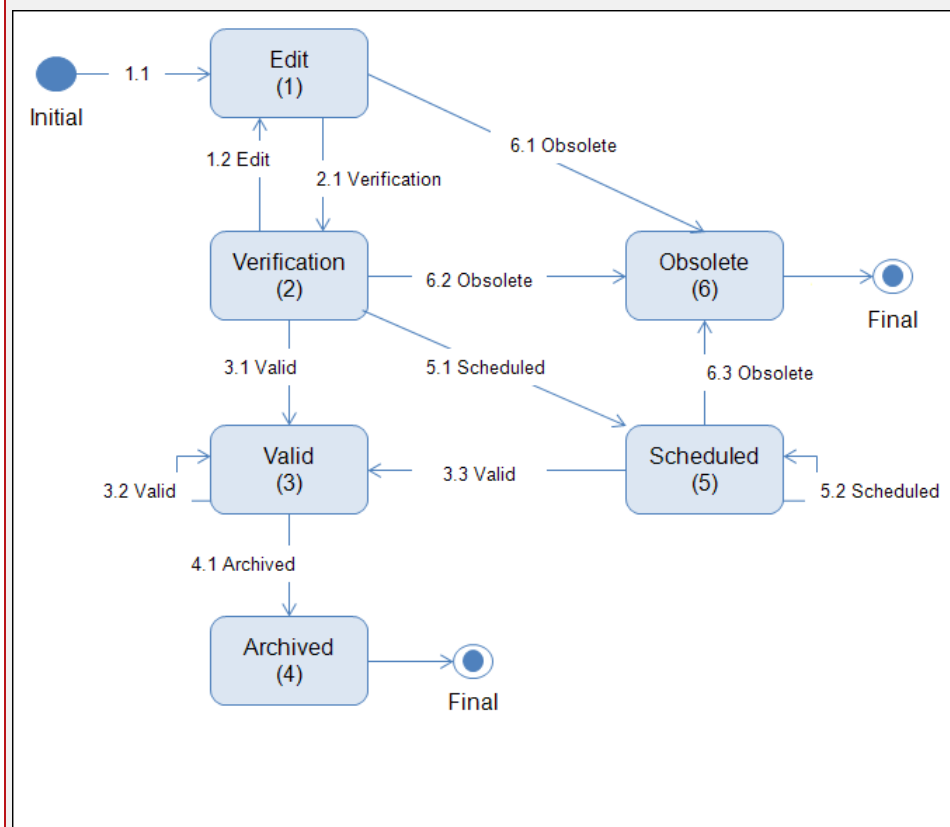


Figure 3: Master recipe - status graph

Default Version Graph for Master Workflows (SR3146.6.4.4)

The **version control (SR3146.4.5)** function (page 141) supports the following default transitions for master workflows:

Transition (ID - From » To)	Signature (Access privilege, according to FSM)	Additional information
1.1 - Initial » Edit	---	Automatic transition.
1.2 - Verification » Edit	Vers_Trans_WF_Verification-Edit	Default user groups: Workflow Author, Qualified Person, Supervisor (Shop Floor)
2.1 - Edit » Verification	Vers_Trans_WF_Edit-Verification	Default user groups: Workflow Author, Qualified Person, Supervisor (Shop Floor)
3.1 - Verification » Valid	Vers_Trans_WF_Verification-Valid	Default user groups: Workflow Author, Qualified Person
3.2 - Valid » Valid	Vers_Trans_WF_Valid-Valid	E.g. if period of validity has to be changed. Default user groups: Workflow Author, Qualified Person
3.3 - Scheduled » Valid	Vers_Trans_WF_Scheduled-Valid	Automatic transition when the current date reaches the start date of the period of validity. Default user groups: Workflow Author, Qualified Person
4.1 - Valid » Archived	Vers_Trans_WF_Valid-Archived	Automatic transition when the end date of the period of validity is in the past. Automatic transition when the status of another version of a master workflow is set to Valid (except if there is a version in the Scheduled status). Manual transition when a status transition to Archived is executed by the workflow author. Default user groups: Workflow Author, Qualified Person
5.1 - Verification » Scheduled	Vers_Trans_WF_Verification-Scheduled	If the scheduled period overlaps with a period of validity, the end date of the period of validity is moved to a date directly before the start date of the scheduled period. Default user groups: Workflow

Transition (ID - From » To)	Signature (Access privilege, according to FSM)	Additional information
		Author, Qualified Person
5.2 - Scheduled » Scheduled	Vers_Trans_WF_Scheduled-Scheduled	E.g. if period of validity has to be changed. Default user groups: Workflow Author, Qualified Person
6.1 - Edit » Obsolete	Vers_Trans_WF_Edit-Obsolete	Default user groups: Workflow Author, Qualified Person, Supervisor (Shop Floor)
6.2 - Verification » Obsolete	Vers_Trans_WF_Verification-Obsolete	Default user groups: Workflow Author, Qualified Person, Supervisor (Shop Floor)
6.3 - Scheduled » Obsolete	Vers_Trans_WF_Scheduled-Obsolete	Default user groups: Workflow Author, Qualified Person

The diagram illustrates the **MESMasterWorkflowVersionGraph** flexible state model.

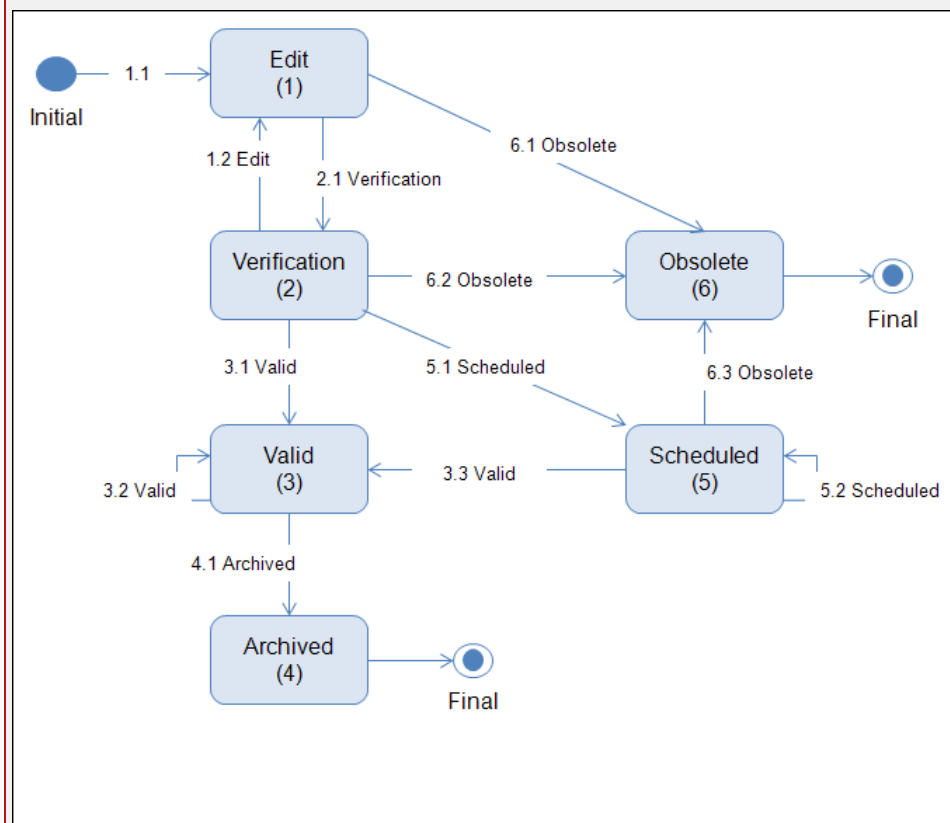


Figure 4: Master workflow - status graph

Validator Consistency Checks (SR3146.9.4.3+)

The following messages are available to inform the recipe or workflow author about issues found by the validator consistency check. They are displayed in the **Messages** (SR3146.9.2.16+) window (page [128](#)).

Checks that Apply to SFC Modeling

Validator - Dummy elements (SR3146.9.4.3.1)

Message category: Error

Issue: Dummy elements are used within a master recipe, master workflow, or (custom) building block.

Validator - Loop (SR3146.9.4.3.3)

➤ For recent changes, see revision history (page [149](#)).

Message category: Error

Issue: Loops are built between unit procedures of a master recipe or (custom) building block.

Validator - Outgoing link (SR3146.9.4.3.12)

Message category: Error

Issue: There is no outgoing link defined for a procedural element.

Validator - Start and end step (SR3146.9.4.3.13)

Message category: Error

Issue: Start and end step are not on the same branch.

Validator - Simultaneous branch (SR3146.9.4.3.14)

Message category: Error

Issue: The end of a simultaneous branch is not merged as a simultaneous branch.

Validator - Selection branch (SR3146.9.4.3.15)

Message category: Error

Issue: The end of a selection branch is not merged as a selection branch.

Validator - Simultaneous branch - Unsuitable components (SR3146.9.4.3.17)

Message category: Error

Issue: The end of a simultaneous branch connects unsuitable components (e.g. simultaneous branch includes open ends).

Validator - Errors in preceding components (SR3146.9.4.3.18)

Message category: Warning

Issue: Elements could not be validated due to errors in the preceding components.

Validator - Start step path (SR3146.9.4.3.19)

Message category: Error

Issue: Elements cannot be reached from the start step.

Validator - End step path (SR3146.9.4.3.20)

Message category: Error

Issue: Elements cannot reach the end step.

Checks that Apply to Transition Conditions

Validator - Transition condition is missing (SR3146.9.4.3.27)

Message category: Error

Issue: There is no condition expression defined for a transition of a selection branch.

Validator - Transition condition - Syntax (SR3146.9.4.3.31)

Message category: Error

Issue: There is a syntax error in the expression of a transition condition.

Validator - Transition condition - Broken link (SR3146.9.4.3.37)

Message category: Depends on the context (see child requirements)

Issue: The expression of a transition condition includes a broken link that needs to be resolved.

Validator - Transition condition - Broken link - Warning (SR3146.9.4.3.37.1)

Context: Custom building block and transition is **not** flagged for **locking**

Message category: Warning

Validator - Transition condition - Broken link - Error (SR3146.9.4.3.37.2)

Context: Custom building block and transition is flagged for **locking**, master recipe, or master workflow

Message category: Error

Validator - Transition condition - Invalid reference (SR3146.9.4.3.38)

Message category: Error

Issue: The phase referenced by a transition condition is invalid.

Checks that Apply to Expressions within Parameter Attributes

Validator - Phase parameter - Syntax (SR3146.9.4.3.21)

Message category: Error

Issue: There is a syntax error in the expression of a parameter attribute.

Validator - Phase parameter - Broken link (SR3146.9.4.3.25)

Message category: Depends on the context (see child requirements)

Issue: The expression of a parameter attribute includes a broken link that needs to be resolved.

Validator - Phase parameter - Broken link - Warning (SR3146.9.4.3.25.1)

Context: Custom building block and parameter is **not** flagged for **locking**

Message category: Warning

Validator - Phase parameter - Broken link - Error (SR3146.9.4.3.25.2)

Context: Custom building block and parameter is flagged for **locking**, master recipe, or master workflow

Message category: Error

Validator - Phase parameter - Invalid reference (SR3146.9.4.3.36)

Message category: Error

Issue: The phase referenced by a parameter attribute is invalid.

Checks that Apply to MFC Modeling of a Procedure within Recipe and Workflow Designer

Validator - MFC - Final output does not match material (product) (SR3146.9.4.3.61)

- Does not apply to Recipe Designer - Device and Workflow Designer.

Message category: Error

Issue: The final output material of an MFC definition does not match the material (product) defined within the master recipe header.

Validator - MFC - Multiple outputs for non-Dispense unit procedures (SR3146.9.4.3.62)

- Does not apply to Recipe Designer - Device and Workflow Designer.

Message category: Depends on the context (see child requirements)

Issue: For non-standalone Dispense unit procedures, there are multiple direct MFC outputs defined or the one defined output is not the final output of the recipe.

Validator - MFC - Multiple outputs for non-Dispense unit procedures - Warning (SR3146.9.4.3.62.1)

Context: Custom building block

Message category: Warning

Validator - MFC - Multiple outputs for non-Dispense unit procedures - Error (SR3146.9.4.3.62.2)

Context: Master recipe

Message category: Error

Validator - MFC - No output defined for input (SR3146.9.4.3.63)

- For recent changes, see revision history (page [149](#)).
- Does not apply to Recipe Designer - Device and Workflow Designer.

Message category: Depends on the context (see child requirements)

Issue: There is no target output defined for a direct MFC input or an input switch.

Validator - MFC - No output defined for input - Warning (SR3146.9.4.3.63.1)

Context: Custom building block

Message category: Warning

Validator - MFC - No output defined for input - Error (SR3146.9.4.3.63.2)

Context: Master recipe

Message category: Error

Validator - MFC - No input defined for output (SR3146.9.4.3.64)

- For recent changes, see revision history (page [149](#)).
- Does not apply to Recipe Designer - Device and Workflow Designer.

Message category: Depends on the context (see child requirements)

Issue: There is no predecessor input defined for a direct MFC output or an output switch.

Validator - MFC - No input defined for output - Warning (SR3146.9.4.3.64.1)

Context: Custom building block

Message category: Warning

Validator - MFC - No input defined for output - Error (SR3146.9.4.3.64.2)

Context: Master recipe

Message category: Error

Validator - MFC - Output defined for device-specific recipe or workflow (SR3146.9.4.3.110)

- Does not apply to Recipe Designer - Batch.

Message category: Error

Issue: An output has been defined for a device-specific master recipe or a master workflow.

Validator - MFC - No direct output defined (SR3146.9.4.3.117)

- For recent changes, see revision history (page 149).

- Does not apply to Recipe Designer - Device and Workflow Designer.

Message category: Depends on the context (see child requirements)

Issue: There is no direct output defined or at least, for one execution path, no direct output is defined.

Validator - MFC - No direct output defined - Warning (SR3146.9.4.3.117.1)

- For recent changes, see revision history (page 149).

Context: Custom building block

Message category: Warning

Validator - MFC - No direct output defined - Error (SR3146.9.4.3.117.2)

- For recent changes, see revision history (page 149).

Context: Master recipe

Message category: Error

Validator - MFC - Attributes of material parameters do not match (SR3146.9.4.3.136)

- For recent changes, see revision history (page 149).

- Does not apply to Recipe Designer - Device and Workflow Designer.

Message category: Error

Issue: The attributes of material parameters with the same MFC position are different.

Validator - MFC - Outgoing or incoming connector is missing (SR3146.9.4.3.138)

- For recent changes, see revision history (page 149).

- Does not apply to Recipe Designer - Device and Workflow Designer.

Message category: Error

Issue: For at least one execution path, there is no outgoing connector defined for a switch output or no incoming connector defined for a switch input.

Validator - MFC - More than one outgoing or incoming connector on one execution path (SR3146.9.4.3.139)

- For recent changes, see revision history (page 149).
- Does not apply to Recipe Designer - Device and Workflow Designer.

Message category: Error

Issue: For at least one execution path, a switch output has more than one outgoing connector or a switch input has more than one incoming connector.

Checks that Apply to Master Recipes in General

Validator - Recipe - Reference quantity not defined (SR3146.9.4.3.67)

- Does not apply to Recipe Designer - Device and Workflow Designer.

Message category: Error

Issue: The planned quantity of a master recipe or the UoM of the planned quantity is missing.

Per default configuration, the check is disabled for master recipes of the **Cost center** usage type.

Validator - Recipe - Min/Max quantities without UoM (SR3146.9.4.3.111)

- Does not apply to Workflow Designer.

Message category: Error

Issue: The minimum or maximum quantity-related unit of measure of a master recipe is missing or the unit of measure cannot be converted with the available conversion factors.

Validator - Master recipe and material parameter - Conflicting inventory level (SR3146.9.4.3.52)

- For recent changes, see revision history (page 149).

Message category: Error

Issue: There are conflicts in the inventory level definition of the packaging levels of the master recipe and its material output parameters. Applies only to material output parameters that refer to the material (product) of the master recipe.

Validator - Packaging level of master recipe (SR3146.9.4.3.107)

- For recent changes, see revision history (page 149).

Message category: Error

Issue: For the packaging level defined in the master recipe header there is a contained number defined but no meaning.

Checks that Apply to Master Recipes Based on an ERP BOM

Validator - Recipe - Planned quantity does not fit ERP BOM (SR3146.9.4.3.57)

- Does not apply to Workflow Designer.

Message category: Error

Issue: The planned quantity of the master recipe does not match the base quantity of the related ERP BOM.

Validator - Recipe - Planned quantity does not fit ERP BOM - Resolve (SR3146.9.4.3.57.1)

The inconsistency can be resolved with the **Messages - Resolve issue (SR3146.9.2.16.4)** function (page 129). The system automatically updates the following data:

- Planned quantity of the master recipe

The automatic resolution of the ERP BOM discrepancy is also included in the **Workbench Object Operations - Replace Material (Product) (SR3146.9.10.11)** operation (page 36).

Validator - Recipe - Packaging level does not fit ERP BOM (SR3146.9.4.3.58)

- For recent changes, see revision history (page 149).

- Does not apply to Workflow Designer.

Message category: Error

Issue: The packaging level (content) of the master recipe does not match the packaging level (content) of the related ERP BOM.

Validator - Recipe - Packaging level does not fit ERP BOM - Resolve (SR3146.9.4.3.58.1)

- For recent changes, see revision history (page 149).

The inconsistencies can be resolved with the **Messages - Resolve issue (SR3146.9.2.16.4)** function (page 129). The system automatically updates the following data:

- Packaging level of the master recipe
- Packaging level content of the master recipe

The automatic resolution of the ERP BOM discrepancies is also included in the **Workbench Object Operations - Replace Material (Product) (SR3146.9.10.11)** operation (page 36).

Validator - Phase parameter - Missing ERP BOM position (SR3146.9.4.3.50)

- Does not apply to Workflow Designer.

Message category: Error

Issue: In a master recipe that is based on an ERP BOM, there are positions that represent a direct MFC input missing in the material parameters.

Validator - Phase parameter - Multiple ERP BOM positions (SR3146.9.4.3.51)

- For recent changes, see revision history (page 149).
- Does not apply to Workflow Designer.

Message category: Error

Issue: There are positions defined multiple times in the material parameters of a master recipe or procedure building block that is based on an ERP BOM. The positions are not unique throughout an entire execution path.

Validator - Phase parameter - Non-ERP BOM position (SR3146.9.4.3.55)

- Does not apply to Workflow Designer.

Message category: Error

Issue: In a master recipe that is based on an ERP BOM, there are material parameter positions that represent direct MFC inputs, which are not part of the ERP BOM.

Validator - Phase parameter - ERP BOM attribute discrepancy (SR3146.9.4.3.56)

- Does not apply to Workflow Designer.

Message category: Error

Issue: At least one of the common attributes of a position of a material parameter and the master recipe-related ERP BOM item do not match.

Common attributes are material identifier and short description, quantity with UoM, Fixed quantity.

Validator - Phase parameter - ERP BOM attribute discrepancy - Resolve (SR3146.9.4.3.56.1)

The inconsistencies can be resolved with the **Messages - Resolve issue (SR3146.9.2.16.4)** function (page 129). The system automatically updates the following data:

- Material inputs by BOM position (material identifier, material short description, planned quantity, and fixed quantity)
- In case linked MFC inputs/outputs, according to the MFC definition, refer to the same material:
The system also updates the material identifier and planned quantity of the related material outputs (same unit procedure, e.g. Dispense, planned quantity only for

the first linked output) and the material identifier of the subsequent material inputs (subsequent unit procedure, e.g. Mixing).

The automatic resolution of the ERP BOM discrepancies is also included in the **Workbench Object Operations - Replace Material (Product) (SR3146.9.10.11)** operation (page 36).

Validator - MFC (final output) - Planned quantity does not fit master recipe (SR3146.9.4.3.133)

➤ Does not apply to Workflow Designer.

Message category: Information

Issue: The planned quantity of the final output does not match the planned quantity to be produced by the master recipe.

Checks that Apply to Material Parameters of a Phase

Validator - Minimum required material parameters (SR3146.9.4.3.10)

Message category: Depends on the context (see child requirements)

Issue: The number of minimum required material input or output parameters of a phase is not met.

Validator - Minimum required material parameters - Warning (SR3146.9.4.3.10.1)

Context: Custom building block

Message category: Warning

Validator - Minimum required material parameters - Error (SR3146.9.4.3.10.2)

Context: Master recipe and master workflow

Message category: Error

Per default configuration, the check is disabled for master recipes of the **Cost center** usage type.

Validator - MFC - Input without planned quantity (SR3146.9.4.3.66)

Message category: Depends on the context (see child requirements)

Issue: The planned quantity of an MFC-related input parameter that represents a direct input is missing.

Validator - MFC - Input without planned quantity - Warning (SR3146.9.4.3.66.1)

Context: Custom building block, on procedure level only

Message category: Warning

Validator - MFC - Input without planned quantity - Error (SR3146.9.4.3.66.2)

Context: Master recipe and master workflow

Message category: Error

Validator - MFC - Uniqueness and format of MFC position (SR3146.9.4.3.68)

- For recent changes, see revision history (page 149).

Message category: Error

Issue:

- The MFC positions assigned to MFC-related parameters (input and output) are not unique throughout an entire execution path or
- for filler materials, the assigned MFC position is not numeric.
The filler-related check only applies in the context of batch-related recipes and custom building blocks on procedure level.

Validator - Material - Incompatible UoM (SR3146.9.4.3.69)

- For recent changes, see revision history (page 149).

Message category: Error

Issue: The UoM assigned to an attribute of a material parameter is not supported or incompatible with the UoM of the material. The following attributes are affected:

- Planned quantity, Target weight
- Non-weighable UoM in case of a Dispense operation
Does not apply to Recipe Designer - Device and Workflow Designer.
- UoM cannot be converted in case of an input or output material parameter of a Dispense operation
Does not apply to Recipe Designer - Device and Workflow Designer.
- Absolute tolerances (of planned quantity and target weight)
- Relative tolerances (of planned quantity and target weight)

Validator - Material - UoM of planned potency (SR3146.9.4.3.70)

Message category: Error

Issue: The UoM of the planned potency of a material parameter is not suitable.

Validator - Material - One tolerance is missing (SR3146.9.4.3.71)

- For recent changes, see revision history (page 149).

Message category: Error

Issue: Only one tolerance of a material parameter is defined and either the lower or the upper tolerance is missing. If a target weight has been defined, either its lower or upper tolerance is missing.

Validator - Material - Negative values (SR3146.9.4.3.77)

- For recent changes, see revision history (page 149).

Message category: Error

Issue: The defined quantity or tolerance for a material parameter attribute is negative. If a target weight has been defined, its quantity or tolerance is negative.

Validator - MFC - MFC position is missing (SR3146.9.4.3.78)

Message category: Depends on the context (see child requirements)

Issue: There is no MFC position assigned to an MFC-related input parameter that represents a direct input.

Validator - MFC - MFC position is missing - Warning (SR3146.9.4.3.78.1)

Context: Custom building block, on procedure level only

Message category: Warning

Validator - MFC - MFC position is missing - Error (SR3146.9.4.3.78.2)

Context: Master recipe and master workflow

Message category: Error

Validator - BOM position not used (SR3146.9.4.3.137)

- For recent changes, see revision history (page 149).

Message category: Error

Issue: For a particular BOM position, there is no corresponding direct MFC input on at least one execution path.

Validator - MFC - Unallowed planned quantity mode (SR3146.9.4.3.113)

Message category: Depends on the context (see child requirements)

Issue: The planned quantity mode setting of an MFC-relevant material parameter is not allowed.

- The **As produced** planned quantity mode is only allowed for ingoing MFC transfer items.
- The **None** planned quantity mode is only allowed for materials of the Auxiliary substance weighing material type (or if no weighing material type is set).

Validator - MFC - Unallowed planned quantity mode - Warning (SR3146.9.4.3.113.1)

Context for the **As produced** planned quantity mode: Custom building block, on procedure level only

Message category: Warning

Validator - MFC - Unallowed planned quantity mode - Error (SR3146.9.4.3.113.2)

Context for the **As produced** planned quantity mode: Master recipe and master workflow

Context for the **None** planned quantity mode: Custom building block, master recipe, and master workflow

Message category: Error

Validator - MFC - Fixed quantity not allowed (SR3146.9.4.3.114)

Message category: Information

Issue: An MFC-relevant output parameter is set to have a fixed quantity.

Validator - MFC - Fixed quantity not allowed for Active or Compensator (SR3146.9.4.3.128)

Message category: Information

Issue: An MFC-relevant material input parameter whose weighing material type is **Active** or **Compensator** is set to have a fixed quantity.

Validator - MFC - Planned quantities of input and output parameters do not match (SR3146.9.4.3.116)

Message category: Information

Issue: The planned quantity of an MFC-relevant input parameter does not match the planned quantity of its output parameter.

Validator - MFC - Planned quantity mode is As produced - Tolerances are missing (SR3146.9.4.3.119)

Message category: Error

Issue: There are no tolerances defined for an MFC-relevant parameter with **As produced** as planned quantity mode.

Validator - Material parameter - Conflicting inventory level (SR3146.9.4.3.53)

- For recent changes, see revision history (page [149](#)).

Message category: Error

Issue: There are conflicts in the inventory level definition of the packaging levels.

Applies only to material parameters that do not refer to the material (product) of the master recipe.

Validator - Packaging level of material parameter (SR3146.9.4.3.108)

- For recent changes, see revision history (page [149](#)).

Message category: Error

Issue: A packaging level of a material parameter has a contained number defined but no meaning.

Validator - Material parameter - Packaging level does not fit material's master data (SR3146.9.4.3.109)

- For recent changes, see revision history (page 149).

Message category: Error

Issue: The packaging level (content) of the material parameter does not match the packaging level (content) of the material's master data.

Validator - Material parameter - Packaging level does not fit material's master data - Resolve (SR3146.9.4.3.109.1)

- For recent changes, see revision history (page 149).

The inconsistencies can be resolved with the **Messages - Resolve issue (SR3146.9.2.16.4)** function (page 129). The system automatically updates the following data:

- Packaging level of the material parameter
- Packaging level content of the material parameter

The automatic resolution of the ERP BOM discrepancies is also included in the **Workbench Object Operations - Replace Material (Product) (SR3146.9.10.11)** operation (page 36).

Checks that Apply to Specific Dispense Building Blocks

Validator - Dispense operation (SR3146.9.4.3.35)

Message category: Error

Issue: The unit procedure of a Dispense operation holds more than one operation.

Validator - MFC - Transfer input for Dispense unit procedures (SR3146.9.4.3.65)

- Does not apply to Recipe Designer - Device and Workflow Designer.

Message category: Depends on the context (see child requirements)

Issue: An MFC transfer is defined as an input material for a Dispense unit procedure.

Validator - MFC - Transfer input for Dispense unit procedures - Warning (SR3146.9.4.3.65.1)

Context: Custom building block

Message category: Warning

Validator - MFC - Transfer input for Dispense unit procedures - Error (SR3146.9.4.3.65.2)

Context: Master recipe

Message category: Error

Validator - MFC - Final output for Dispense unit procedures (SR3146.9.4.3.118)

- Does not apply to Recipe Designer - Device and Workflow Designer.

Message category: Depends on the context (see child requirements)

Issue: One of the materials of a Dispense unit procedure is merged with the final output.

Validator - MFC - Final output for Dispense unit procedures - Warning (SR3146.9.4.3.118.1)

Context: Custom building block

Message category: Warning

Validator - MFC - Final output for Dispense unit procedures - Error (SR3146.9.4.3.118.2)

Context: Master recipe

Message category: Error

Validator - MFC - None planned quantity mode for Dispense unit procedure (SR3146.9.4.3.115)

Message category: Error

Issue: An MFC-relevant material parameter with **None** as planned quantity mode is defined for a Dispense unit procedure.

Validator - D_Material - Default weighing method not defined (SR3146.9.4.3.72)

- Does not apply to Recipe Designer - Device and Workflow Designer.

Message category: Error

Issue: The default weighing method is not defined for a material parameter of a Dispense operation.

Validator - D_Material - Default weighing method not allowed (SR3146.9.4.3.132)

- Does not apply to Recipe Designer - Device and Workflow Designer.

Message category: Error

Issue: The default weighing method of a material input parameter of a Dispense operation is not among its allowed weighing methods.

Validator - D_Material - Allowed weighing methods not defined (SR3146.9.4.3.73)

- Does not apply to Recipe Designer - Device and Workflow Designer.

Message category: Error

Issue: The allowed weighing methods are not defined for a material parameter of a Dispense operation.

Validator - D_Material - Tolerances not defined (SR3146.9.4.3.74)

- For recent changes, see revision history (page 149).
- Does not apply to Recipe Designer - Device and Workflow Designer.

Message category: Error

Issue:

- Dispense: The tolerances are not defined for a material parameter of a Dispense operation.
- Output weighing: If a target weight has been defined, its tolerances are not defined.

Validator - D_Material - Weighing material type not defined (SR3146.9.4.3.75)

- Does not apply to Recipe Designer - Device and Workflow Designer.

Message category: Error

Issue: The weighing material type is not defined for a material parameter of a Dispense operation.

Validator - D_Material - Planned potency not defined (SR3146.9.4.3.76)

- Does not apply to Recipe Designer - Device and Workflow Designer.

Message category: Error

Issue: The planned potency is not defined for an active material of a Dispense operation.

Checks that Apply to Inline Weighing

Validator - MFC - Unallowed weighing method (SR3146.9.4.3.120)

Message category: Error

Issue: For the **Identify material** phase of an Inline Weighing operation, **Pallet** weighing is not allowed for MFC-relevant material parameters.

Validator - MFC - Planned quantity mode is As defined - Tolerances are missing (SR3146.9.4.3.121)

Message category: Error

Issue: For the **Identify material** phase of an Inline Weighing operation, there are no tolerances defined for an MFC-relevant material parameter with **As defined** as planned quantity mode.

Validator - MFC - Planned quantity mode is As defined - Planned quantity is missing (SR3146.9.4.3.122)

Message category: Error

Issue: For the **Identify material** phase of an Inline Weighing operation, there is no planned quantity defined for an MFC-relevant material parameter with **As defined** as planned quantity mode.

Validator - MFC - Weighing material type not defined (SR3146.9.4.3.123)

Message category: Error

Issue: For the **Identify material** phase of an Inline Weighing operation, there is no weighing material type defined for an MFC-relevant material parameter.

Validator - MFC - MFC position is missing (SR3146.9.4.3.124)

Message category: Error

Issue: For the **Identify material** phase of an Inline Weighing operation, there is no MFC position assigned to an MFC-related input material parameter. This applies to ingoing MFC transfer item material parameters only.

Validator - MFC - Ingoing MFC transfer item - Unallowed weighing material type (Active) (SR3146.9.4.3.125)

Message category: Error

Issue: For the **Identify material** phase of an Inline Weighing operation, the **Active** weighing material type is not allowed for an ingoing MFC transfer item.

See also **Unallowed weighing material type (Compensator, Filler)** (SR3146.9.4.3.134) message (page 63).

Validator - MFC - Ingoing MFC transfer item - Unallowed weighing material type (Compensator, Filler) (SR3146.9.4.3.134)

Message category: Error

Issue: For the **Identify material** phase of an Inline Weighing operation, the **Compensator** and **Filler** weighing material types are not allowed for an ingoing MFC transfer item if they were already set to a weighing material type other than **Auxiliary** before in at least one preceding material position.

See also **Unallowed weighing material type (Active)** (SR3146.9.4.3.125) message (page 63).

Validator - MFC - Ingoing MFC transfer item - Unallowed planned quantity mode (SR3146.9.4.3.135)

Message category: Error

Issue: For the **Identify material** phase of an Inline Weighing operation, the **As produced** planned quantity mode is not allowed for an ingoing MFC transfer item of the **Compensator** and **Filler** weighing material types.

Validator - D_Material - Planned potency not defined (SR3146.9.4.3.127)

➤ Does not apply to Recipe Designer - Device and Workflow Designer.

Message category: Error

Issue: The planned potency is not defined for an active material of an Inline Weighing operation.

Checks that Apply to Output Weighing

Validator - MFC - Planned quantity mode is As defined - Tolerances are missing (SR3146.9.4.3.129)

Message category: Error

Issue: For the **Weigh** phase of an Output Weighing operation, there are no tolerances defined for an MFC-relevant material parameter with **As defined** as planned quantity mode.

Validator - MFC - Planned quantity mode is As defined - Planned quantity is missing (SR3146.9.4.3.130)

Message category: Error

Issue: For the **Weigh** phase of an Output Weighing operation, there is no planned quantity defined for an MFC-relevant material parameter with **As defined** as planned quantity mode.

Validator - MFC - Planned quantity mode is As defined - Allowed weighing method is missing (SR3146.9.4.3.131)

Message category: Error

Issue: For the **Weigh** phase of an Output Weighing operation, there is no allowed weighing method defined for an MFC-relevant material parameter with **As defined** as planned quantity mode.

Additional checks for Output Weighing

- Validator - D_Material - Tolerances Not Defined (SR3146.9.4.3.74) check (page [62](#))

Checks that Apply to Non-production Master Recipes in General

Validator - Non-production MR - More than one unit procedure or operation (SR3146.9.4.3.100)

Message category: Error

Issue: The procedural structure of a non-production master recipe contains more than one unit procedure or operation.

Checks that are disabled for master recipes of the Cost center usage type

- Validator - Minimum required material parameters - Error (SR3146.9.4.3.10.2) check (page 56)
- Validator - Recipe - Reference quantity not defined (SR3146.9.4.3.67) check (page 53)

Checks that Apply to Master Workflows in General

Validator - Workflow - Phase is incompatible (SR3146.9.4.3.80)

Message category: Error

Issue: The inserted custom building block contains phases that are not compatible with master workflows.

Validator - Workflow - Workflow type not defined (SR3146.9.4.3.81)

Message category: Error

Issue: The workflow type is missing that defines to which group the one-click startable workflow belongs in the Cockpit of the Production Execution Client.

Validator - Workflow - Only appendable if production-relevant (SR3146.9.4.3.112)

Message category: Error

Issue: A master workflow is not production-relevant, but is set to be appendable during processing.

Checks that Apply to Signature Privileges of a Phase

Validator - Usage type (SR3146.9.4.3.8)

Message category: Error

Issue: The same usage type is defined for multiple signature privileges of a single phase.

Validator - Exception signature supported (SR3146.9.4.3.9)

Message category: Error

Issue: An exception signature is defined, although exceptions are not supported by the phase.

Validator - Double signature for sequential phase completion (SR3146.9.4.3.5)

Message category: Error

Issue: A double signature for a sequential phase completion signature has been defined.

Validator - (Sequential) phase completion signatures (SR3146.9.4.3.6)

Message category: Error

Issue: One phase completion signature and one sequential phase completion signature have been assigned to a phase.

Checks that Apply to Equipment Requirements of a Phase

Validator - Property definition is inconsistent (SR3146.9.4.3.82)

Message category: Error

Issue: The content definition of a property type does not match the content definition of the respective class property.

Validator - Minimum required equipment parameters (SR3146.9.4.3.83)

Message category: Depends on the context (see child requirements)

Issue: The number of minimum required equipment requirement parameters of a phase is not met.

Validator - Minimum required equipment parameters - Warning (SR3146.9.4.3.83.1)

Context: Custom building block

Message category: Warning

Validator - Minimum required equipment parameters - Error (SR3146.9.4.3.83.2)

Context: Master recipe and master workflow

Message category: Error

Validator - Equipment class requirements are not unique (SR3146.9.4.3.85)

Message category: Error

Issue: The equipment requirements on class level are not unique. One equipment class has been assigned multiple times.

Validator - Equipment requirement - BOE is not unique (SR3146.9.4.3.88)

Message category: Error

Issue: The BOE position assigned to an equipment requirement parameter is not unique throughout the master recipe, master workflow, or building block.

Validator - Equipment requirement identifier is not unique (SR3146.9.4.3.86)

Message category: Error

Issue: The equipment requirement identifier is not unique within a class.

Validator - Equipment requirement rule - Identifier is missing (SR3146.9.4.3.87)

Message category: Error

Issue: There is no identifier defined for an equipment requirement rule.

Validator - Equipment requirement rule - Syntax error (SR3146.9.4.3.89)

Message category: Error

Issue: There is a syntax error in the expression of an equipment requirement rule.

Validator - Equipment requirement rule - Broken link (SR3146.9.4.3.90)

Message category: Depends on the context (see child requirements)

Issue: The expression of an equipment requirement rule includes a broken link that needs to be resolved.

Validator - Equipment requirement rule - Broken link - Warning (SR3146.9.4.3.90.1)

Context: Custom building block and parameter is **not** flagged for **locking**

Message category: Warning

Validator - Equipment requirement rule Broken link - Error (SR3146.9.4.3.90.2)

Context: Custom building block and parameter is flagged for **locking**, master recipe, or master workflow

Message category: Error

Validator - Equipment requirement rule - Invalid reference (SR3146.9.4.3.91)

Message category: Error

Issue: The phase referenced by an equipment requirement rule is invalid.

Validator - Equipment requirement rule - Rule expression is missing (SR3146.9.4.3.92)

Message category: Error

Issue: There is no rule expression defined for an equipment requirement rule.

Checks that Apply to Work Center Assignments

Station assignment does not match a given work center assignment (SR3146.9.4.3.40)

Message category: Error

Issue: Per unit procedure, at least one station assigned on operation level does not match the work centers assigned on unit procedure level. Therefore, work centers have to be added or stations have to be removed.

Station assignment is missing for a given work center assignment (SR3146.9.4.3.42)

Message category: Error

Issue: At least one station is already assigned to an operation. However, there are more work centers assigned on unit procedure level and the dispatching of this operation does not yet reflect the stations of all of the work centers.

Unusable station (SR3146.9.4.3.43)

Message category: Error

Issue: At least one station of this operation cannot be used because it does not belong to any work center.

Unusable work center (SR3146.9.4.3.44)

Message category: Error

Issue: At least one work center of this unit procedure cannot be used because it has no station defined.

Checks that Apply to Capability Assignments in General

Usage of capability not allowed for workflows (SR3146.9.4.3.79)

Message category: Error

Issue: The operation/unit procedure contains the Abort-and-reactivate-enabled capability. But the capability is not permitted for use in workflows.

Usage of capability not allowed for Dispense (SR3146.9.4.3.45)

Message category: Error

Issue: The operation/unit procedure contains at least one of the following capabilities: Abort-and-reactivate-enabled, Auto-startable, Detachable, Escalation-enabled, Event-triggered, Server-run, or Trigger-enabled. These capabilities, however, are not permitted for use in Dispense operations or unit procedures.

Usage of capability not allowed for simultaneous branches (SR3146.9.4.3.106)

Message category: Error

Issue: The operation contains the Auto-startable capability. But this capability is not permitted for use in operations located on a simultaneous branch.

Capability assignment does not match usage type of privilege parameter (SR3146.9.4.3.93)

Message category: Depends on the context (see child requirements)

Issue: The operation or unit procedure contains a capability, but has no signature privilege of the same usage type to match the capability's risk assessment.

Capability assignment does not match usage type of privilege parameter - Warning (SR3146.9.4.3.93.1)

Context: Custom building block

Message category: Warning

Capability assignment does not match usage type of privilege parameter - Error (SR3146.9.4.3.93.2)

Context: Master recipe and master workflow

Message category: Error

Incompatible capability assignment (SR3146.9.4.3.46)

Message category: Error

Issue: Capabilities that are assigned to the same operation are not compatible. This applies to the following combinations:

Capability	Auto-startable	Server-run
Auto-startable	N/A	Incompatible*
Detachable	Compatible	Incompatible
Escalation-enabled	Incompatible	Incompatible
Event-triggered	Incompatible	Incompatible
Trigger-enabled	Compatible	Incompatible

*Server-run operations will start automatically by design. The Auto-startable capability only applies as an option to operations that run on a client.

Incompatible capability assignment between different levels (SR3146.9.4.3.105)

Message category: Error

Issue: Capabilities that are assigned to a unit procedure are not compatible with the capabilities of its operations. This applies to the following combinations:

Capability	Detachable (unit procedure)
Server-run (operation)	Incompatible
Event-triggered (operation)	Incompatible

Capability assignment is missing (SR3146.9.4.3.47)

Message category: Error

Issue: A capability has been assigned, which also requires another capability. This applies to the Trigger-enabled and Escalation-enabled capabilities, which are only allowed in combination with the Event-triggered capability.

Checks that Apply to Escalation-enabled Capability Assignments

Escalation validator - Duration is missing for enabled exceptions (SR3146.9.4.3.94)

Message category: Error

Issue: The Escalation-enabled capability has been assigned and exceptions are enabled, but not all related durations are defined.

Escalation validator - Durations are conflicting (SR3146.9.4.3.95)

Message category: Error

Issue: The Escalation-enabled capability has been assigned, but the durations defined for its actions are conflicting.

Escalation validator - No duration defined (SR3146.9.4.3.96)

Message category: Depends on the context (see child requirements)

Issue: The Escalation-enabled capability has been assigned, but no duration is defined.

Escalation validator - No duration defined - Warning (SR3146.9.4.3.96.1)

Context: Custom building block

Message category: Warning

Escalation validator - No duration defined - Error (SR3146.9.4.3.96.2)

Context: Master recipe and master workflow

Message category: Error

Checks that Apply to Trigger-enabled Capability Assignments

Trigger validator - No valid trigger reference defined (SR3146.9.4.3.97)

Message category: Depends on the context (see child requirements)

Issue: A Trigger-enabled capability contains triggers with a missing or invalid trigger phase reference.

Trigger validator - No valid trigger reference defined - Warning (SR3146.9.4.3.97.1)

Context: Custom building block

Message category: Warning

Trigger validator - No valid trigger reference defined - Error (SR3146.9.4.3.97.2)

Context: Master recipe and master workflow

Message category: Error

Trigger validator - Trigger identifiers are not unique (SR3146.9.4.3.98)

Message category: Error

Issue: A Trigger-enabled capability contains triggers whose identifiers are not unique.

Trigger validator - Triggers are missing (SR3146.9.4.3.99)

Message category: Error

Issue: A Trigger-enabled capability does not contain any triggers.

Trigger phase without event-triggered operation (SR3146.9.4.3.101)

Message category: Information

Issue: A trigger phase is used in the recipe or workflow that is not referenced from an event-triggered operation.

Checks that Apply to Server-run Capability Assignments

Server-run validator - Conflicting station assignment (SR3146.9.4.3.102)

Message category: Error

Issue: The Server-run capability has been assigned, but stations are also assigned.

Server-run validator - Phase cannot be run on a server (SR3146.9.4.3.103)

Message category: Error

Issue: The Server-run capability has been assigned, but the phase cannot be run on a server.

Server-run validator - Phase needs to be run on a server (SR3146.9.4.3.104)

Message category: Error

Issue: The Server-run capability has not been assigned, but the phase needs to be run on a server.

Server-run validator - No phase completion signature or mandatory comment allowed (SR3146.9.4.3.7)

Message category: Error

Issue: A (sequential) phase completion signature or a signature with a mandatory comment has been assigned to a phase that needs to be run on a server.

Management of Building Blocks (SR3146.10+)

The management of building blocks differs from the management of master recipes and master workflows in the following areas:

Header Attributes (SR3146.10.4)

➤ For recent changes, see revision history (page [149](#)).

The building block object supports the following header attributes:

- Basic data:
 - Identifier
 - Revision
 - Short description
 - Comparison baseline
 - Access privilege
 - Category 1 - 3
 - Description
- Status data:
 - Status
 - Approved on
 - Approved by
 - Archived on
 - Archived by

The attributes are maintained in the **header component-specific (SR3146.9.2.15)** property window (page [128](#)).

Workbench Object Operations (SR3146.9.10.7)

➤ For recent changes, see revision history (page 149).

The following operations are supported for custom building blocks. They apply to elements on all levels (procedure, unit procedure, operation, phase).

- Create new building block (page 75)
- Open existing (see **Open (SR3146.9.2.20+)** dialog (page 122))
- View locked building block (page 75)
- Save (current or all)
- Save as new custom BB
For custom building blocks that are protected by an access privilege: The access privilege of the original custom building block is passed on to the new custom building block. In the default configuration, the access privilege cannot be changed.
- Rename
- Replace (overrides parameters) via the **Setlist (SR3146.9.2.9)** UI element (page 124)
- Smart replace via the **Setlist (SR3146.9.2.9)** UI element (page 124)
(see **Workbench Object Operations - Smart Replace (SR3146.9.10.10)** operation (page 32))
- Delete
- Status Handling of Building Blocks (SR3146.10.1) (page 75)
- Graph pagination
- Usage List (SR3146.10.3) (page 76)
- Statistics (SR3146.9.2.18) (page 140)
- Close (current or all)

- A building block (procedure, unit procedure, operation, phase) consists of a header and its underlying graphs, except for phases since they represent the lowest level.
Available revisions of building blocks with the same identifier are listed for information.
- While renaming BB: Available revisions of building blocks with the same identifier are listed for information.
- If a BB is copied (save as), the new building block is available for selection from the **Universe (SR3146.9.2.4+)** UI element (page 122) and the **Setlist (SR3146.9.2.9)** UI element (page 124) and is opened in the current tab. Available revisions of building blocks with the same identifier are listed for information.

Workbench Object Operations - New Building Block Dialog (SR3146.9.10.7.2)

- For recent changes, see revision history (page 149).

When a new building block is created, first, the building block author defines the building block identifier and revision.

The system shall allow to select an access privilege for the protection of the building block from unauthorized access for users with access privileges that allow to maintain protected building blocks. The list of available access privileges shall be restricted to the access privileges of the logged-in user.

Workbench Object Operations - View Locked Building Block (SR3146.9.10.7.1)

If a building block is locked by another user, the system allows to view the building block without the possibility to change it.

The building block can be copied (save as). Operations that do not change a building block are also available (e.g. graph pagination, usage list, statistics).

Status Handling of Building Blocks (SR3146.10.1)


Status handling of custom building blocks is supported (**Draft**, **Approved**, **Archived**).

The system automatically locks all transitions without an identifier when a building block is approved.

Building blocks in the **Archived** status are available in the **Universe (SR3146.9.2.4+)** UI element (page 122) (of Recipe and Workflow Designer) and the **Setlist (SR3146.9.2.9)** UI element (page 124), but they cannot be added to a recipe, workflow, or custom building block.

- **Draft:** Initial status of building blocks.
- **Approved:** Building blocks can no longer be modified. If an **Approved** building block is used within a master recipe, master workflow, or a custom building block, parameters or transitions that were locked are frozen and cannot be unlocked. Parameters or transitions that are not locked can still be modified. Underlying graphs are read-only. See "Locking of Parameters (SR3146.9.13+)" (page 27) and "Locking of Transitions (SR3146.9.14+)" (page 27).
- **Archived:** They are still available in the **Universe (SR3146.9.2.4+)** UI element (page 122) (of Recipe and Workflow Designer) and the **Setlist (SR3146.9.2.9)** UI element (page 124), but can only be added to a **change request SR3146.12+** object (page 83) (as old building block) from the Setlist. They are accessible via the **Open (SR3146.9.2.20+)** function (page 122).

Each status change must be confirmed by a single electronic signature.

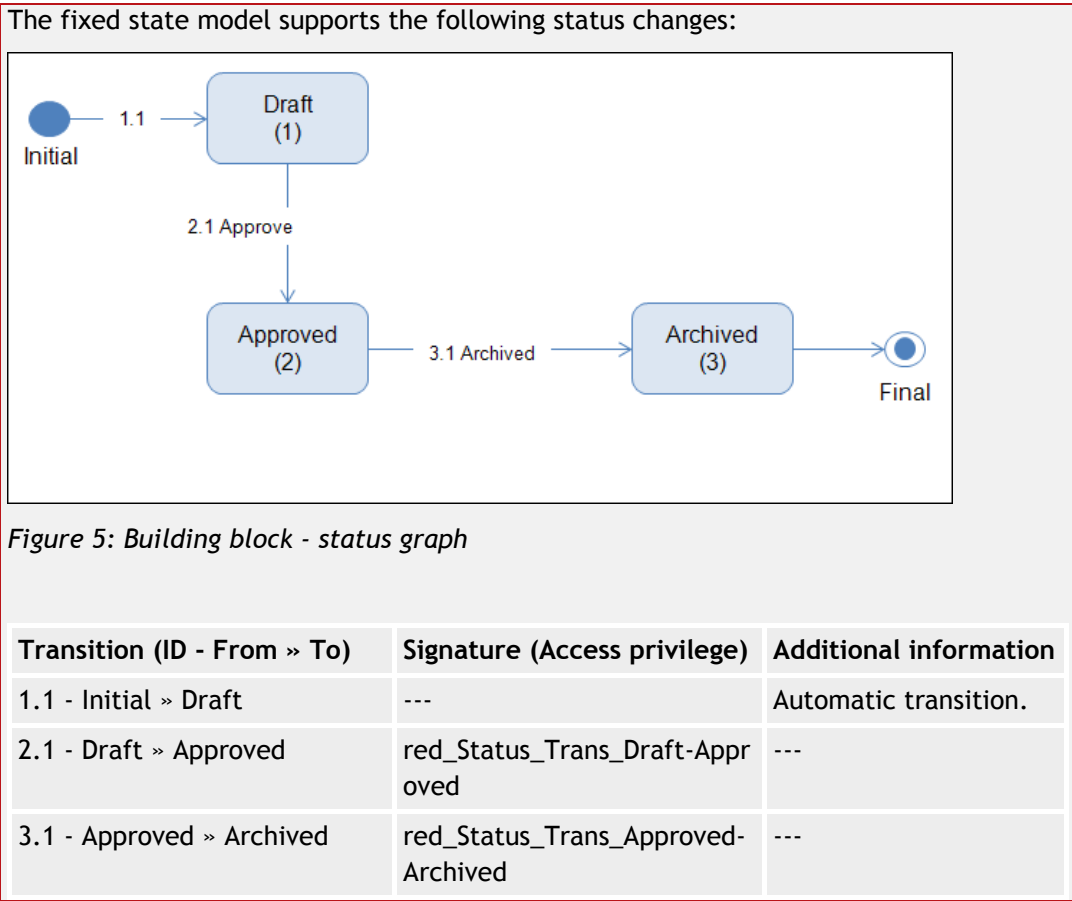

Revisioning (SR3146.10.2)

Revisioning of custom building blocks is supported as a mandatory free-text property.

Unique revisioning is requested when a building block is created (see "Workbench Object Operations (SR3146.9.10.7)" (page 74), "Workbench Object Operations - Element (SR3146.9.10.6)" (page 31)).

Usage List (SR3146.10.3)

- For recent changes, see revision history (page 149).
- The usage list of a custom building block provides all objects that have been designed based on the custom building block in question.
- The following object types are supported by the usage list: custom building blocks, master recipes, and master workflows.
- For objects that are protected by an access privilege: The system shall only provide objects whose access privilege matches the access privilege of the logged-in user.

Content (SR3146.10.3.1)

The usage list allows to track:

- Direct usage of a custom building block
All objects that contain elements that were derived directly from the custom building block in question, without using an intermediate-level-approved custom building block.
Example: A custom phase was built directly into a master recipe.
- Indirect usage of a custom building block
All objects that contain elements that were not derived directly from the custom building block in question, but that were derived from an intermediate-level-approved custom building block.
Example: A custom phase is used to build and approve a custom operation, which then was built directly into a master recipe.

Objects Information (SR3146.10.3.2)

- For recent changes, see revision history (page [149](#)).

Each object that is listed in the usage list provides the following information:

- Object type
- Identifier, version, and status of the object
- Usage (direct, indirect)
- Used building block level
- [Used building block] Identifier, version, and status
 - Direct usage: Identifier of the custom building block for which the usage list has been compiled.
 - Indirect usage: Identifier of the intermediate-level-approved custom building.
- Path information of the recipe, workflow, custom building block elements

Additionally, the usage list provides a list of changed objects whose objects information cannot be refreshed.

Filter (SR3146.10.3.3)

The usage list supports the following filter options:

- Include indirect usages
- Include non-approved artifacts

-
-
- Rockwell Software PharmaSuite® 8.4 - Functional Requirement Specification Recipe and Workflow Management
-
-

Comparing Master Recipes, Master Workflows, and Custom Building Blocks (SR3146.15+)

PharmaSuite supports the approval of master recipes, master workflows, and custom building blocks that are based on an approved version with a comparison mechanism.

The comparison baseline of a master recipe, master workflow, or custom building block defines the object to be compared with the current object.

The baseline is a header attribute of the **master recipe (SR3146.9.3.1)** object (page 9), the **master workflow (SR3146.9.12.1)** object (page 11), and the **building block (SR3146.10.4)** object (page 73).

The comparison result is displayed in the **Comparison (SR3146.9.2.19+)** window (page 130). It informs a reviewer of the recipe, workflow, or custom building block of the differences found by the comparison mechanism. Thus, the reviewer gets a summary of the differences.

Comparison Algorithm (SR3146.15.1+)

The comparison algorithm covers attributes, parameters, procedural elements within their SFC context, MFC-related settings, and the execution-relevant layout of a master recipe, master workflow, or custom building block.

Comparing Header Attributes (SR3146.15.1.1+)

Master recipe-specific header attributes (SR3146.15.1.1.1)

Each attribute is compared with its baseline, except for:

- Identifier
- Version
- Comparison baseline
- Status
- Valid from and Valid until

Master workflow-specific header attributes (SR3146.15.1.1.2)

Each attribute is compared with its baseline, except for:

- Identifier

- Version
- Comparison baseline
- Status
- Valid from and Valid until

Custom building block-specific header attributes (SR3146.15.1.1.3)

Each attribute is compared with its baseline, except for:

- Identifier
- Revision
- Comparison baseline
- Status
- Approved on and Approved by
- Archived on and Archived by

Comparing Element Attributes (SR3146.15.1.2)

Each attribute is compared with its baseline.

Comparing Source Attributes (SR3146.15.1.3)

Only the identifier attribute is compared with its baseline.

The following attributes are not compared:

- Short description
- Usable in
- Trigger
- Server-run
- Pause-aware

Comparing Parameter Attributes (SR3146.15.1.4+)

Material parameter-specific attributes (SR3146.15.1.4.1)

Each attribute is compared with its baseline, except for:

- Material short description

Material Flow Control-specific attributes (SR3146.15.1.4.2)

Each attribute is compared with its baseline.

Equipment requirement parameter-specific attributes (SR3146.15.1.4.3)

Each attribute is compared with its baseline, except for:

- Short description
- Description
- EQM level

Work center assignment parameter-specific attributes (SR3146.15.1.4.4)

Each attribute is compared with its baseline, except for:

- Description
- Storage area

Privilege parameter-specific attributes (SR3146.15.1.4.5)

Each attribute is compared with its baseline, except for:

- Default reason (1st and 2nd)

Capability parameter-specific attributes (SR3146.15.1.4.6)

Each attribute is compared with its baseline.

Process parameter-specific attributes (SR3146.15.1.4.7)

Each attribute is compared with its baseline, except for:

- Content (concatenation of other attributes that are being compared)

Comparing Procedural Elements (SR3146.15.1.5+)

Transition-specific attributes (SR3146.15.1.5.2)

Each attribute is compared with its baseline.

Relation-specific characteristics (SR3146.15.1.5.1)

The following attributes are compared with its baseline:

- Existence of elements
- Successor relations of elements
- Branching type (sequence, selection, simultaneous)

Location-specific characteristics (SR3146.15.1.5.3)

The following attribute is compared with its baseline:

- Position (relative layout) of elements in simultaneous branches

The position of elements in simultaneous branches can affect the preview rendering during execution and the sequence of elements within the batch report.

Defining a Comparison Baseline (SR3146.15.2)

The system allows to define a comparison baseline for master recipes and master workflows in the **Edit** status (page 43) and building blocks in the **Draft** status (page 75).

The following rules apply to comparison baseline objects:

- For master recipes or master workflows, they must be in the **Scheduled**, **Valid**, or **Archived** status (page 43).
- For building blocks, they must be in the **Approved** or **Archived** status (page 75).

Performing a Comparison (SR3146.15.3)

To perform a comparison, a comparison baseline must be assigned. As long as no comparison baseline is assigned, the system displays an appropriate message.

When a comparison has been performed, the system displays the result in the list panel of the **Comparison (SR3146.9.2.19+)** window (page 130).

- If no differences have been detected, the system displays an appropriate message.
- If differences have been detected, the system displays the details.
- If the current object has been changed in the meantime at the current client, the system displays an appropriate message.

Management of Change Requests (SR3146.12+)

A change request consists of its **header attributes (SR3146.12.1)** data (page 83), the **scope definition (SR3146.12.3+)** feature (page 89), and the **action list (SR3146.12.4)** feature (page 90).

The scope of a change request is defined by the old (to be replaced) and new (replacing) custom building blocks and by the master recipes/master workflows and custom building blocks selected from the usage list that shall be subject to changing. The action list shows all master recipes/master workflows and custom building blocks which are to be changed.

During the **execution of a change request (SR3146.12.5)** operation (page 91), the replacement of building blocks follows the same rules as for the **Smart Replace (SR3146.9.10.10)** operation (page 32).

The replacement of used system building blocks is not supported by the change request.

The life cycle of a change request is controlled by the **status management of change requests (SR3146.12.2+)** feature (page 85).

Header Attributes (SR3146.12.1)

➤ For recent changes, see revision history (page 149).

The change request object supports the following header attributes:

- Basic data:
 - Identifier
 - Description
 - Comparison baseline handling (Keep, Update, Delete)
- Recipe-/workflow-related data:
 - Target status (Edit, Verification, All valid statuses)
 - Effective date
- Status data:
 - Status
- Additional data:
 - Structure level
 - New building block

- [New BB] Revision
- [New BB] Access privilege
- Old building block
- [Old BB] Revision
- [Old BB] Access privilege

The attributes are maintained in the **header component-specific (SR3146.9.2.15)** property window (page 128).

The **Target status** attribute defines if the related new recipe or workflow versions are automatically moved to the **Valid** status (page 43).

The **Comparison baseline handling** attribute allows to

- keep the comparison baseline of the old component to also be the baseline of the new component;
- update (default setting) the comparison baseline of the new component according to the most recent version/revision from which it was created;
- delete the comparison baseline, so that for the new component no comparison baseline is defined even if a comparison baseline was defined for the old component.

Workbench Object Operations (SR3146.12.6)

The following operations are supported for change requests:

- Create new change request
- Open existing (see **Open (SR3146.9.2.20+)** dialog (page 122))
- View locked change request (page 85)
- Save (current or all)
- Copy (Save as...)
- Rename
- Delete
- Close (current or all)

Open Dialog for Protected Objects (SR3146.12.6.2)

- For recent changes, see revision history (page 149).

A change request can hold objects that are protected by an access privilege, like the new or old building block or any object that is in the scope of the change request according the action list (master recipes, master workflows, and custom building blocks). In case a protected object is involved: The system shall only allow to open change requests whose

access privileges of any of the protected objects match the access privilege of the logged-in user.

Workbench Object Operations - View Locked Change Request (SR3146.12.6.1)

If a change request is locked by another user, the system allows to view the change request without the possibility to change it.

The change request can be copied (save as). Operations that do not change a change request are also available (e.g. status history).

Status Management of Change Requests (SR3146.12.2+)

Master Recipe-specific Change Requests (SR3146.12.2.2)

The **version control (SR3146.4.5)** function (page 141) supports the following default transitions for master recipe-specific change requests:

Transition (ID - From » To)	Signature (Access privilege, according to FSM)	Additional information
1.1 - Initial » Draft	---	Automatic transition.
1.2 - In review » Draft	Status_Trans_ChangeRequest_InReview-Draft	Default user groups: Recipe Author, Qualified Person, Supervisor (Shop Floor)
2.1 - Draft » In review	Status_Trans_ChangeRequest_Draft-InReview	Default user groups: Recipe Author, Qualified Person, Supervisor (Shop Floor)
3.1 - In review » Approved	Status_Trans_ChangeRequest_InReview-Approved	Default user group: Qualified Person
4.1 - Approved » In process	Status_Trans_ChangeRequest_Approved-InProcess	Automatic transition when the change request is being executed. Default user groups: Recipe Author, Qualified Person, Supervisor (Shop Floor)
5.1 - In process » Finished	---	Automatic transition when the change request execution is completed.
5.2 - Aborted » Finished	Status_Trans_ChangeRequest_Canceled-Finished	Default user groups: Recipe Author, Qualified Person, Supervisor (Shop Floor)
5.3 - In process » Finished	Status_Trans_ChangeRequest_InProcess-Finished	Manual transition, only available if the execution has been interrupted (e.g. resume during execution). Default user groups: Recipe Author,

Transition (ID - From » To)	Signature (Access privilege, according to FSM)	Additional information
		Qualified Person, Supervisor (Shop Floor)
6.1 - Finished » Archived	---	---
6.2 - Obsolete » Archived	---	---
7.1 - Draft » Obsolete	---	---
7.2 - In review » Obsolete	Status_Trans_ChangeRequest_InReview-Obsolete	Default user groups: Recipe Author, Qualified Person, Supervisor (Shop Floor)
7.3 - Approved » Obsolete	Status_Trans_ChangeRequest_Approved-Obsolete	Default user groups: Recipe Author, Qualified Person, Supervisor (Shop Floor)
7.4 - Aborted » Obsolete	Status_Trans_ChangeRequest_Canceled-Obsolete	Default user groups: Recipe Author, Qualified Person, Supervisor (Shop Floor)
7.5 - In process » Obsolete	Status_Trans_ChangeRequest_InProcess-Obsolete	Manual transition, only available if the execution has been interrupted (e.g. resume during execution). Default user groups: Recipe Author, Qualified Person, Supervisor (Shop Floor)
8.1 - In process » Aborted	Status_Trans_ChangeRequest_InProcess-Canceled	Automatic transition when the change request is canceled with the ESC key.

The diagram illustrates the **MESChangeRequestStatusGraph** flexible state model.

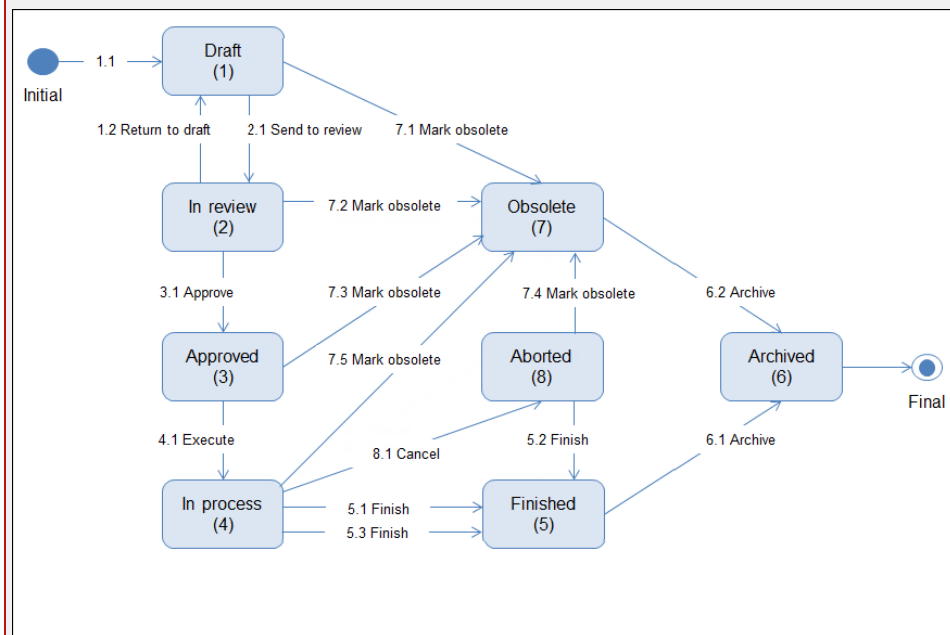


Figure 6: Master recipe-specific change request - status graph

Master Workflow-specific Change Requests (SR3146.12.2.3)

The **version control (SR3146.4.5)** function (page 141) supports the following default transitions for master workflow-specific change requests:

Transition (ID - From » To)	Signature (Access privilege, according to FSM)	Additional information
1.1 - Initial » Draft	---	Automatic transition.
1.2 - In review » Draft	Status_Trans_ChangeRequest_WF_InReview-Draft	Default user groups: Workflow Author, Qualified Person, Supervisor (Shop Floor)
2.1 - Draft » In review	Status_Trans_ChangeRequest_WF_Draft-InReview	Default user groups: Workflow Author, Qualified Person, Supervisor (Shop Floor)
3.1 - In review » Approved	Status_Trans_ChangeRequest_WF_InReview-Approved	Default user groups: Workflow Author, Qualified Person
4.1 - Approved » In process	Status_Trans_ChangeRequest_WF_Approved-InProcess	Automatic transition when the change request is being executed. Default user groups: Workflow Author, Qualified Person, Supervisor (Shop Floor)

Transition (ID - From » To)	Signature (Access privilege, according to FSM)	Additional information
5.1 - In process » Finished	---	Automatic transition when the change request execution is completed.
5.2 - Aborted » Finished	Status_Trans_ChangeRequest_WF_Canceled-Finished	Default user groups: Workflow Author, Qualified Person, Supervisor (Shop Floor)
5.3 - In process » Finished	Status_Trans_ChangeRequest_WF_InProcess-Finished	Manual transition, only available if the execution has been interrupted (e.g. resume during execution). Default user groups: Workflow Author, Qualified Person, Supervisor (Shop Floor)
6.1 - Finished » Archived	---	---
6.2 - Obsolete » Archived	---	---
7.1 - Draft » Obsolete	---	---
7.2 - In review » Obsolete	Status_Trans_ChangeRequest_WF_InReview-Obsolete	Default user groups: Workflow Author, Qualified Person, Supervisor (Shop Floor)
7.3 - Approved » Obsolete	Status_Trans_ChangeRequest_WF_Approved-Obsolete	Default user groups: Workflow Author, Qualified Person, Supervisor (Shop Floor)
7.4 - Aborted » Obsolete	Status_Trans_ChangeRequest_WF_Canceled-Obsolete	Default user groups: Workflow Author, Qualified Person, Supervisor (Shop Floor)
7.5 - In process » Obsolete	Status_Trans_ChangeRequest_WF_InProcess-Obsolete	Manual transition, only available if the execution has been interrupted (e.g. resume during execution). Default user groups: Workflow Author, Qualified Person, Supervisor (Shop Floor)
8.1 - In process » Aborted	Status_Trans_ChangeRequest_WF_InProcess-Canceled	Automatic transition when the change request is canceled with the ESC key.

The diagram illustrates the **MESChangeRequestMasterWorkflowStatusGraph** flexible state model.

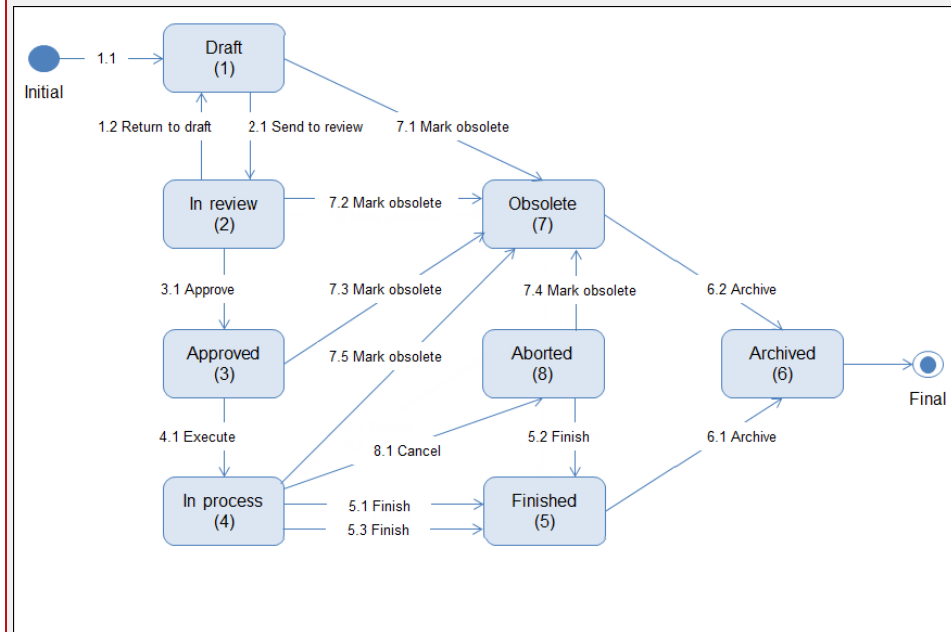


Figure 7: Master workflow-specific change request - status graph

Scope Validity Check (SR3146.12.2.1)

During the status transitions of a change request to **In review**, **Approved**, and **In process**, the system checks if the defined scope of the change request is still valid. This applies to the status of the new building block and to the status/version of any affected object (master recipe, master workflow, or custom building block).

Scope Definition (SR3146.12.3+)

The scope definition of a change request is composed of the **definition of old/new buildings blocks (SR3146.12.1)** task (page 89) and the **selection from the usage list (SR3146.12.3.2)** task (page 90).

Define Old/New Building Blocks (SR3146.12.3.1)

- For recent changes, see revision history (page 149).

The system allows to define an old building block and a new building block for a change request. During the execution of the change request, the old building block is replaced by the new building block.

The following rules apply:

- All building block levels of the procedural structure (page 13) are supported (procedure, unit procedure, operation, phase).
- An old building block must be in the **Approved** or **Archived** statuses (page 75).
- A new building block must be in the **Approved** status (page 75).
- The access privilege for the protection of a building block from unauthorized access has no impact on the replacement. It is always performed.

Select from Usage List (SR3146.12.3.2)

- For recent changes, see revision history (page 149).

Based on the definition of the old building block, the system allows to compile a **Usage list (SR3146.10.3)** overview (page 76) and to select a subset of artifacts to be changed within the scope of the change request.

The following rules apply:

- The **Compile usage list** function is only available for change requests in the **Draft** status (page 85).
- The usage list only includes master recipes and master workflows in the **Valid** or **Scheduled** statuses (page 43) and building blocks in the **Approved** or **Archived** statuses (page 75).
Only those usages are displayed for which the old building block was in the **Approved** status (page 75) when it was copied (i.e. **Draft** usages are filtered out).
- If the status of a master recipe, master workflow, or building block has changed since the last compilation of the usage list and the change request is still in the **Draft** status, the system notifies the user of the status change and updates the usage list. Example: Master recipes whose status has been changed to **Archived** are removed from the list.
- If an object (e.g. master recipe) selected for a change includes indirect usages of the old building block, the system selects in the usage list the affected custom building blocks to be changed within the scope of the change request.
- For objects that are protected by an access privilege: The system shall only provide objects whose access privilege matches the access privilege of the logged-in user.

Action List (SR3146.12.4)

The action list represents all master recipes, master workflows, and custom building blocks that are in the scope of the change request, i.e., that are selected in the usage list and hence require a new version/revision due to the defined building block replacement.

The action list provides the following information for each artifact:

- Object type
- Identifier
- Old object [version/revision and status before and after executing the change request]
- New object [version/revision and status]
- Progress
- Result
- Error

Execute Change Request (SR3146.12.5)

- For recent changes, see revision history (page 149).

During the automatic execution of a change request, new versions/revisions of the affected master recipes/master workflows/custom building blocks are created and the defined building blocks are replaced.

Related to the change request, the following rules apply:

- Only change requests in the **Approved** status (page 85) can be executed.
- The system sets the status of a change request automatically to **In process** when the change request is being executed.
- The system sets the status of a change request automatically to **Finished** when the change request execution is completed.

Related to the artifacts displayed in the **Action list (SR3146.12.4)** tab (page 90), the following rules apply:

- The replacement of building blocks follows the rules of the **Smart Replace (SR3146.9.10.10)** operation (page 32).
- The system automatically transitions new revisions of building blocks to the same status (page 75) the old revisions had.
- Master recipes and master workflows in the **Valid** or **Scheduled** status (page 43) are all considered to represent valid recipes and workflows. New versions will be transitioned according to their originally defined period of validity.
- The system automatically transitions new versions of master recipes and master workflows to the target status that is defined in the **header attributes (SR3146.12.1)** data (page 83) of the change request.
- The following descriptions reflect that the **Target status** attribute has been set to **All valid statuses**.

(If it has been set to **Edit** or **Verification**, new versions of master recipes and master workflows in the **Edit** or **Verification** statuses are created. Old versions of affected master recipes or master workflows are not changed, i.e. their statuses remain unchanged.)

- If the effective date of the change request is not defined (or today), an affected master recipe or master workflow in the **Valid** status is archived and a new version of the master recipe or master workflow in the **Valid** status is created.
 - If the effective date of the change request is in the future and is due while an affected master recipe or master workflow is in the **Valid** status, the master recipe's or master workflow's period of validity is shortened and a new version of the master recipe or master workflow in the **Scheduled** status is created. For this new version in the **Scheduled** status, the period of validity starts on the effective date of the change request and lasts until the end of the original period of validity.
 - If the effective date of the change request is due before the period of validity of an affected master recipe or master workflow in the **Scheduled** status starts, a new version of the master recipe or master workflow is created in the **Scheduled** status based on the original period of validity of the previous version.
 - If the effective date of the change request is in the future and is due while an affected master recipe or master workflow is in the **Scheduled** status, the master recipe's or master workflow's future period of validity is shortened and another version of the master recipe or master workflow in the **Scheduled** status is created. For this new version, the period of validity starts on the effective date of the change request and lasts until the end of the original period of validity.
 - For all status transitions related to master recipes and master workflows that require an electronic signature, a silent signature is recorded whose comment is filled with the data of the signing operator(s) that approved the change request.
 - In case an error occurs during execution of a change request, the result for the current artifact is recorded in the **Action list (SR3146.12.4)** tab (page 90) and the system proceeds with the next artifact to be changed.
Examples for error conditions: Status has changed since the compilation of the usage list, procedural structure is broken, etc.
- As only master recipes and master workflows in the **Valid** or **Scheduled** status (page 43) can be part of the action list, master recipes and master workflows can never be transitioned by a change request to a status that is higher than the original master recipe's or master workflow's status prior to the execution of the change request.

- For all status transitions of master recipes and master workflows that occur in the context of a change request execution, the master recipe- and master workflows-specific configurations related to signature privileges and e-mail notifications are ignored by the system.

Execution Status (SR3146.12.5.2)

The change request supports the following interim execution statuses:

- Not started
- Processing
- Canceling (triggered by the ESC key)

The change request supports the following final execution statuses:

- Canceled (subsequent to the **Canceling** interim execution status)
- Canceled (unexpected, after a resume of a running change request)
- Completed

Cancel Change Request (SR3146.12.5.1)

During the automatic execution of a change request, the execution of the change request can be canceled.

The following rules apply:

- The system completes the replacement-related steps of the artifacts that are displayed in the **Action list (SR3146.12.4)** tab (page 90) and are currently in process.
- All artifacts displayed in the **Action list (SR3146.12.4)** tab (page 90) that have not been started yet are marked as **Canceled**.
- The system sets the status of the change request automatically to **Aborted** when the execution of the change request has been canceled.

-
-
- Rockwell Software PharmaSuite® 8.4 - Functional Requirement Specification Recipe and Workflow Management
-
-

Master Recipe Report - Batch (SR3146.11+)

The master recipe report - batch represents the S88-compliant structure of a master recipe and its components in a printable way.

Print MR Report (SR3146.11.1)

The system allows the user to print a master recipe report that represents the S88-compliant structure of a master recipe. If desired, the comparison result can be excluded from the report.

MR Report - SFC Pagination (SR3146.11.1.1)

Displays the current graph structure (SFC flow) in a preview window for printing including pagination preview.

MR Report Header (SR3146.11.2)

The header of the MR report includes the following data:

- Placeholder for a customer logo
- Material ID and Short description
- MR identifier
- MR version

MR Report Footer (SR3146.11.3)

The footer of the MR report includes the following data:

- Rockwell Automation logo
- Printed by: <user name and login>
- Printed from: <station/work center>
- Database: <database name>
- Printed on: <local time and server time>
- Page x of y

MR Report Cover Page (SR3146.11.5)

The cover page of the MR report includes the following data:

- Report title
- Identifier
- Description
- Comparison baseline
- Material identifier and Short description
- Planned quantity
- Version and Status
- Valid from and Valid until
- Method and Method description

MR Report Sections (SR3146.11.4+)

The MR report includes the following sections:

MR report - Table of contents (SR3146.11.4.15)

The table of contents of the MR report.

MR report - Basic data (SR3146.11.4.1)

- For recent changes, see revision history (page [149](#)).

Section for the basic data of the master recipe, which includes:

- Identifier
- Description
- Comparison baseline
- Material identifier and Short description
- Planned quantity
- Minimum quantity and Maximum quantity
- Version and Status
- Valid from and Valid until
- Registration number
- Method and Method description
- Usage type

- Review mode
- Access privilege
- Packaging levels

MR report - Comments (SR3146.11.4.2)

Section for the comments related to the master recipe.

MR report - Approval record (SR3146.11.4.3)

Section for the approval record of the current version of the master recipe, which includes:

- For the creation of the recipe:
Action, role, logged-in user, reason for creation
- For the approval of the recipe:
Action, role, signatures with timestamp

MR report - Bill of materials (SR3146.11.4.4)

Formula section contains the bill of materials of the master recipe, which includes the following header attributes and columns:

- Header attributes:
 - Material [identifier and short description]
 - ERP BOM alternative
 - ERP BOM base quantity
 - ERP BOM released
- Columns:
 - Position
 - Material [identifier and short description]
 - Material type
 - Planned quantity
 - Fixed quantity

MR report - Bill of equipment (SR3146.11.4.16)

Formula section contains the bill of equipment of the master recipe. The header currently serves as a placeholder for a future table.

MR report - Work center assignments (SR3146.11.4.5)

Formula section contains the list of work centers of the master recipe, which includes the following columns:

- Type
- Identifier
- Description
- Used in recipe element

MR report - Material flow control (SR3146.11.4.6)

Formula section contains the material flow control data of the master recipe, which includes the following columns:

- MFC item
- Position
- Material [identifier and short description]
- Planned quantity (input)
- Planned quantity (output)
- Output unit procedure
- Input unit procedure
- Target MFC item

MR report - Critical process parameters (SR3146.11.4.7)

Formula section contains the critical process parameters of the master recipe, which includes the following columns:

- Identifier
- Content
- Path

MR report - Critical quality attributes (SR3146.11.4.8)

Formula section contains the critical quality attributes of the master recipe, which includes the following columns:

- Identifier
- Content
- Path

MR report - Key process parameters (SR3146.11.4.9)

Formula section contains the key process parameters of the master recipe, which includes the following columns:

- Identifier
- Content
- Path

MR report - Key quality attributes (SR3146.11.4.10)

Formula section contains the key quality attributes of the master recipe, which includes the following columns:

- Identifier
- Content
- Path

MR report - Procedure (SR3146.11.4.11)

➤ For recent changes, see revision history (page [149](#)).

The procedure-related data of the MR report includes the following details:

- Basic information
 - Element data
 - Custom source data
- List of steps
- List of transitions
- Graphical representation of the SFC flow

MR report - Unit procedure (SR3146.11.4.12)

For each unit procedure, the MR report includes the following details:

- Basic information
 - Element data
 - Custom source data
- Work center information
- Capability information (e.g. Detachable capability, Pause-enabled capability)
- Privilege information
- List of steps

- List of transitions
- Graphical representation of the SFC flow

MR report - Operation (SR3146.11.4.13)

For each operation, the MR report includes the following details:

- Basic information
 - Element data
 - Custom source data
- Station information
- Capability information (e.g. Detachable capability, Event-triggered capability)
- Privilege information
- List of steps
- List of transitions
- Graphical representation of the SFC flow

MR report - Phase (SR3146.11.4.14)

For each phase, the MR report includes the following details:

- Basic information
 - Element data
 - Custom source data
 - System source data
- Input material information
- Dispensing-specific material information (only for Dispense operations)
- Output material information
- Equipment requirement information
- Privilege information
- Process parameter information

MR report - Comparison with baseline (SR3146.11.4.17)

For each comparison-specific type, the MR report includes the following details:

- Mode
- Path
- Type

- Identifier
- Predecessor
- Successor
- Attribute
- New value
- Old value

If no differences have been detected between the master recipe and its comparison baseline, no comparison baseline is defined, or the comparison baseline is not available, the report includes an appropriate message instead of the table.

-
-
- Rockwell Software PharmaSuite® 8.4 - Functional Requirement Specification Recipe and Workflow Management
-
-

Master Recipe Report - Device (SR3146.14+)

The electronic device history record represents the S88-compliant structure of a master recipe and its components in a printable way.

Print MR Report (SR3146.14.1)

The system allows the user to print a master recipe report that represents the S88-compliant structure of a master recipe. If desired, the comparison result can be excluded from the report.

MR Report - SFC Pagination (SR3146.14.1.1)

Displays the current graph structure (SFC flow) in a preview window for printing including pagination preview.

MR Report Header (SR3146.14.2)

The header of the MR report includes the following data:

- Placeholder for a customer logo
- Material ID and Short description
- MR identifier
- MR version

MR Report Footer (SR3146.14.3)

The footer of the MR report includes the following data:

- Rockwell Automation logo
- Printed by: <user name and login>
- Printed from: <station/work center>
- Database: <database name>
- Printed on: <local time and server time>
- Page x of y

MR Report Cover Page (SR3146.14.5)

The cover page of the MR report includes the following data:

- Report title
- Identifier
- Description
- Comparison baseline
- Material identifier and Short description
- Planned quantity
- Version and Status
- Valid from and Valid until
- Method and Method description

MR Report Sections (SR3146.14.4+)

The MR report includes the following sections:

MR report - Table of contents (SR3146.14.4.15)

The table of contents of the MR report.

MR report - Basic data (SR3146.14.4.1)

Section for the basic data of the master recipe, which includes:

- Identifier
- Description
- Comparison baseline
- Material identifier and Short description
- Planned quantity
- Minimum quantity and Maximum quantity
- Version and Status
- Valid from and Valid until
- Registration number
- Method and Method description
- Review mode
- Access privilege

- Packaging levels

MR report - Comments (SR3146.14.4.2)

Section for the comments related to the master recipe.

MR report - Approval record (SR3146.14.4.3)

Section for the approval record of the current version of the master recipe, which includes:

- For the creation of the recipe:
Action, role, logged-in user, reason for creation
- For the approval of the recipe:
Action, role, signatures with timestamp

MR report - Bill of materials (SR3146.14.4.4)

Formula section contains the bill of materials of the master recipe, which includes the following header attributes and columns:

- Header attributes:
 - Material [identifier and short description]
 - ERP BOM alternative
 - ERP BOM base quantity
 - ERP BOM released
- Columns:
 - Position
 - Material [identifier and short description]
 - Material type
 - Planned quantity
 - Fixed quantity

MR report - Bill of equipment (SR3146.14.4.16)

Formula section contains the bill of equipment of the master recipe. The header currently serves as a placeholder for a future table.

MR report - Work center assignments (SR3146.14.4.5)

Formula section contains the list of work centers of the master recipe, which includes the following columns:

- Type
- Identifier
- Description
- Used in recipe element

MR report - Material flow control (SR3146.14.4.6)

Formula section contains the material flow control data of the master recipe, which includes the following columns:

- MFC item
- Position
- Material [identifier and short description]
- Planned quantity (input)
- Planned quantity (output)
- Output unit procedure
- Input unit procedure

MR report - Critical process parameters (SR3146.14.4.7)

Formula section contains the critical process parameters of the master recipe, which includes the following columns:

- Identifier
- Content
- Path

MR report - Critical quality attributes (SR3146.14.4.8)

Formula section contains the critical quality attributes of the master recipe, which includes the following columns:

- Identifier
- Content
- Path

MR report - Key process parameters (SR3146.14.4.9)

Formula section contains the key process parameters of the master recipe, which includes the following columns:

- Identifier
- Content
- Path

MR report - Key quality attributes (SR3146.14.4.10)

Formula section contains the key quality attributes of the master recipe, which includes the following columns:

- Identifier
- Content
- Path

MR report - Procedure (SR3146.14.4.11)

➤ For recent changes, see revision history (page [149](#)).

The procedure-related data of the MR report includes the following details:

- Basic information
 - Element data
 - Custom source data
- List of steps
- List of transitions
- Graphical representation of the SFC flow

MR report - Unit procedure (SR3146.14.4.12)

For each unit procedure, the MR report includes the following details:

- Basic information
 - Element data
 - Custom source data
- Work center information
- Capability information (e.g. Detachable capability, Pause-enabled capability)
- Privilege information
- List of steps

- List of transitions
- Graphical representation of the SFC flow

MR report - Operation (SR3146.14.4.13)

For each operation, the MR report includes the following details:

- Basic information
 - Element data
 - Custom source data
- Station information
- Capability information (e.g. Detachable capability)
- Privilege information
- List of steps
- List of transitions
- Graphical representation of the SFC flow

MR report - Phase (SR3146.14.4.14)

For each phase, the MR report includes the following details:

- Basic information
 - Element data
 - Custom source data
 - System source data
- Input material information
- Output material information
- Equipment requirement information
- Privilege information
- Process parameter information

MR report - Comparison with baseline (SR3146.14.4.17)

For each comparison-specific type, the MR report includes the following details:

- Mode
- Path
- Type
- Identifier

- Predecessor
- Successor
- Attribute
- New value
- Old value

If no differences have been detected between the master recipe and its comparison baseline, no comparison baseline is defined, or the comparison baseline is not available, the report includes an appropriate message instead of the table.

-
-
- Rockwell Software PharmaSuite® 8.4 - Functional Requirement Specification Recipe and Workflow Management
-
-

Master Workflow Report (SR3146.13+)

The master workflow report represents the S88-compliant structure of a master workflow and its components in a printable way.

Print MWF Report (SR3146.13.1)

The system allows the user to print a master workflow report that represents the S88-compliant structure of a master workflow. If desired, the comparison result can be excluded from the report.

MWF Report - SFC Pagination (SR3146.13.1.1)

Displays the current graph structure (SFC flow) in a preview window for printing including pagination preview.

MWF Report Header (SR3146.13.2)

The header of the MWF report includes the following data:

- Placeholder for a customer logo
- Processing name
- Master workflow identifier and Short description
- Master workflow version

MWF Report Footer (SR3146.13.3)

The footer of the MWF report includes the following data:

- Rockwell Automation logo
- Printed by: <user name and login>
- Printed from: <station/work center>
- Database: <database name>
- Printed on: <local time and server time>
- Page x of y

MWF Report Cover Page (SR3146.13.5)

The cover page of the MWF report includes the following data:

- Report title
- Identifier
- Processing name
- Short description
- Version and Status
- Valid from and Valid until

MWF Report Sections (SR3146.13.4+)

The MWF report includes the following sections:

MWF report - Table of contents (SR3146.13.4.13)

The table of contents of the MWF report.

MWF report - Basic data (SR3146.13.4.1)

- For recent changes, see revision history (page [149](#)).

Section for the basic data of the master workflow, which includes:

- Identifier
- Processing name
- Short description
- Description
- Comparison baseline
- Version and Status
- Valid from and Valid until
- Review mode
- Access privilege
- [Workflow] Type
- One-click startable
- Execution prefix
- Production-relevant
- Appendable during processing

MWF report - Comments (SR3146.13.4.2)

Section for the comments related to the master workflow.

MWF report - Approval record (SR3146.13.4.3)

Section for the approval record of the current version of the master workflow, which includes:

- For the creation of the workflow:
Action, role, logged-in user, reason for creation
- For the approval of the workflow:
Action, role, signatures with timestamp

MWF report - List of materials (SR3146.13.4.15)

Formula section contains the list of materials of the master workflow, which includes the following columns:

- Position
- Material [identifier and short description]
- Material type
- Planned quantity
- Fixed quantity

MWF report - Bill of equipment (SR3146.13.4.14)

Formula section contains the bill of equipment of the master workflow. The header currently serves as a placeholder for a future table.

MWF report - Work center assignments (SR3146.13.4.4)

Formula section contains the list of work centers of the master workflow, which includes the following columns:

- Type
- Identifier
- Description
- Used in workflow element

MWF report - Critical process parameters (SR3146.13.4.5)

Formula section contains the critical process parameters of the master workflow, which includes the following columns:

- Identifier
- Content
- Path

MWF report - Critical quality attributes (SR3146.13.4.6)

Formula section contains the critical quality attributes of the master workflow, which includes the following columns:

- Identifier
- Content
- Path

MWF report - Key process parameters (SR3146.13.4.7)

Formula section contains the key process parameters of the master workflow, which includes the following columns:

- Identifier
- Content
- Path

MWF report - Key quality attributes (SR3146.13.4.8)

Formula section contains the key quality attributes of the master workflow, which includes the following columns:

- Identifier
- Content
- Path

MWF report - Procedure (SR3146.13.4.9)

The procedure-related data of the MWF report includes the following details:

- Basic information
 - Element data
 - Custom source data
- List of steps

- Graphical representation of the SFC flow

MWF report - Unit procedure (SR3146.13.4.10)

For each unit procedure, the MWF report includes the following details:

- Basic information
 - Element data
 - Custom source data
- Work center information
- Capability information (e.g. Detachable capability, Pause-enabled capability)
- Privilege information
- List of steps
- List of transitions
- Graphical representation of the SFC flow

MWF report - Operation (SR3146.13.4.11)

For each operation, the MWF report includes the following details:

- Basic information
 - Element data
 - Custom source data
- Station information
- Capability information (e.g. Detachable capability)
- Privilege information
- List of steps
- List of transitions
- Graphical representation of the SFC flow

MWF report - Phase (SR3146.13.4.12)

For each phase, the MWF report includes the following details:

- Basic information
 - Element data
 - Custom source data
 - System source data

- Input material information
- Output material information
- Equipment requirement information
- Privilege information
- Process parameter information

MWF report - Comparison with baseline (SR3146.13.4.17)

For each comparison-specific type, the MWF report includes the following details:

- Mode
- Path
- Type
- Identifier
- Predecessor
- Successor
- Attribute
- New value
- Old value

If no differences have been detected between the master workflow and its comparison baseline, no comparison baseline is defined, or the comparison baseline is not available, the report includes an appropriate message instead of the table.

Designer Workbench (SR3146.9+)

- As requirements have been restructured compared to the original DRS document, this section also contains requirements of other requirement groups than SR3146.9+.

For recipe or workflow authors, the graphical Recipe and Workflow Designer workbench provides the means to create and maintain master recipes, master workflows, and custom building blocks (BB).

This section provides the requirements related to Designer workbench itself:

- **Workbench usability (SR3146.9.1+)** requirements (page 117)
- **Workbench UI elements (SR3146.9.2+)** requirements (page 121)
- **Version control (SR3146.4.5)** requirements (page 141)
- **Library concept (SR3146.9.11+)** requirements (page 143)

Workbench Usability (SR3146.9.1+)

The graphical Designer workbench supports usability aspects related to its graphical user interface and general handling. Both are targeted to minimize time to result significantly.

Building Block Image

The building block image provides access to the following data:

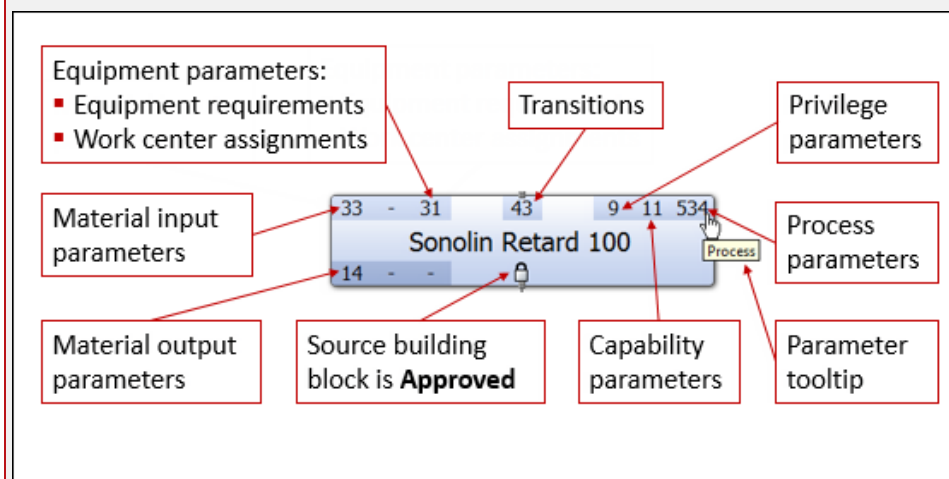


Figure 8: Building block image

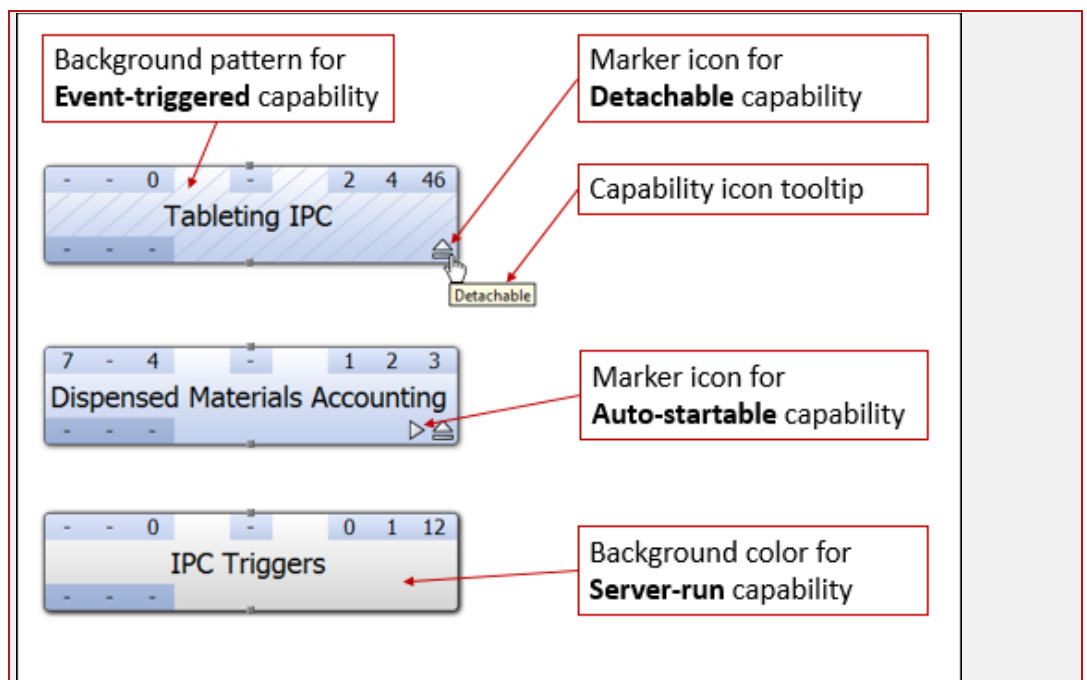


Figure 9: Operation building block image with capabilities

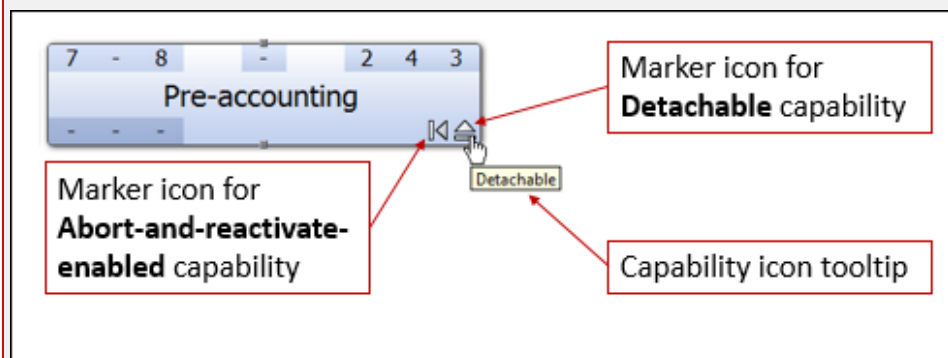


Figure 10: Unit procedure building block image with capabilities

Upper left section:

- Materials as process input (SR3146.9.5.1) data (page 16).
- Equipment parameter:
 Equipment requirement (class) (SR3146.9.5.5) data (page 19)
 Work center (SR3146.9.5.4) data (page 18), for unit procedures only
 Station (SR3146.9.5.8) data (page 18), for operations only.

Lower left section:

- Materials as process output (SR3146.9.6.1) data (page 24).

Upper right section:

- Signature privilege (SR3146.9.5.3) data (page 17).

- Capabilities
 - Detachable - Operation (SR3146.9.5.9.1) capability (page 20)
 - Event-triggered (SR3146.9.5.9.2) capability (page 21)
 - Escalation-enabled (SR3146.9.5.9.3) capability (page 21)
 - Trigger-enabled (SR3146.9.5.9.4) capability (page 22)
 - Server-run (SR3146.9.5.9.5) capability (page 22)
 - Auto-startable (SR3146.9.5.9.6) capability (page 22)
 - Pause-enabled (SR3146.9.5.10.1) capability (page 23)
 - Detachable - Unit procedure (SR3146.9.5.10.2) capability (page 23).
- Capabilities (Does not apply to Workflow Designer.)
 - Abort-and-reactivate-enabled (SR3146.9.5.10.3) capability (page 23).
- Process parameter (SR3146.9.7+) data (page 25).

Upper center section:

- Transition (SR3146.9.5.11) data (page 20).

Center section:

- Name of the building block.

Lower center section:

- Frozen (SR3146.9.13.2) marker (page 27) (if the used building block is in the Approved status).

Process View (SR3146.9.1.1)

A process-related view is available to support the creation and maintenance of procedural elements of the **procedural structure (SR3146.9.4)** hierarchy (page 13).

WORKBENCH - START AND END STEPS (SR3146.9.1.1.1)

On each level of the procedural structure, Start and End steps are available. Thus a recipe or workflow author can model branches and loops that include the first or final building blocks of a graph.

Examples: A process flow starts with parallel steps. The last step of a process flow can have a loop.

Navigation (SR3146.9.1.2)

The following navigation and correction tools are available:

- Zoom
- Map (bird eye)
- Explorer
- Pan current graph

- Hierarchical navigation
- Horizontal navigation (within one master recipe, master workflow, or custom BB)
- Undo last graph-related user actions (independent of the current Graph Window).
- Redo actions revoked by the **Undo** function.

- Zoom - Support of zoom in, zoom out, zoom to fit, and 100% of current graph.
- Map - Displays the current building block and its superordinate building blocks in three views and enables the navigation through the master recipe or master workflow. The view of the current Graph Window is highlighted. The Map applies to procedure level, unit procedure level, and operation level.
- Explorer - Displays a structural overview of the graph and enables the navigation through the master recipe, master workflow, or building block. The frozen marker indicates that a graph element is based on a source building block in the **Approved (SR3146.10.1)** status (page 75).
- Hierarchical navigation - Go up to superordinate procedural level and go back to previous procedural level.

Context menus (SR3146.9.1.3)

Context menus enhance the usability during the creation and maintenance of procedural elements of the **procedural structure (SR3146.9.4)** hierarchy (page 13).

Context menus are available for tab titles of the **Graph Window (SR3146.9.2.1)** UI element (page 121), the work area of the Graph Window, the building block images (page 117), and the **Material Flow Control (SR3146.9.8)** tab (page 37).

Docking (SR3146.9.1.4)

The docking configuration is available as a default layout and a user-specific layout.

In the default docking configuration the windows are arranged as follows:

- To the left of the Graph Window: **Map (SR3146.9.1.2)** tool (page 119) and **Explorer (SR3146.9.1.2)** tool (page 119).
- To the right of the Graph Window: **Setlist (SR3146.9.2.9)** UI element (page 124) and **property windows** (SR3146.9.2.12 (page 128), SR3146.9.2.13 (page 128), SR3146.9.2.15 (page 128)).
- At the bottom of the Graph Window: **Phase Preview (SR3146.9.2.11)** UI element (page 131), **Messages (SR3146.9.2.16+)** window (page 128), **Comparison (SR3146.9.2.19+)** window (page 130).

The following functions are available to maintain the docking configuration:

- Undo last layout changes related to the dockable windows.
- Redo last layout changes related to the dockable windows.
- Save the window layout (user-specific).
- Load the last saved window layout (user-specific).
- Reset the window layout to the default layout.

Tooltips (SR3146.9.1.5)

Tooltips are available for the following objects:

- transitions (page 14),
- parameter buttons in the building block images (page 117), and
- **MFC graph (SR3146.9.8.1)** components (page 37).

Workbench UI Elements (SR3146.9.2+)

The graphical Designer workbench provides several UI elements (e.g. Library Universe, Messages window, Parameter Panel).

Work Area (SR3146.9.2.1)

The work area panel consists of the Graph Window.

- The Graph Window can display more than one building block of a master recipe, master workflow, building block, or change request, each in a separate tab (subordinate tab).
- Subordinate tabs can be closed separately.
Does not apply to change requests.
- The tab title of subordinate tabs displays the building block identifiers.
Does not apply to change requests.
- An asterisk in the upper tab title is the indicator for unsaved master recipes, master workflows, building blocks, or change requests.
- An R in the upper tab title is the indicator for read-only master recipes, master workflows, building blocks, or changes requests.
- All tabs can be closed one at a time except for the current master recipe, master workflow, building block, or change request.
- The graph and its components can be arranged automatically in order to center sequences and branches with respect to their start components and to optimize the vertical spacing between individual components.

- The graph can be compacted automatically for optimal use of space.
- Grid lines can be shown or hidden.

Open Dialog (SR3146.9.2.20+)

- For recent changes, see revision history (page 149).

The **Open** dialog is related to **Master recipe (SR3146.9.10.1)** objects (page 28), **Master workflow (SR3146.9.10.9)** objects (page 29), **Custom building block (SR3146.9.10.7)** objects (page 74), and **Change request (SR3146.12+)** objects (page 83). It allows to select data records for opening in a resizable dialog. It provides an intelligent filter and search function.

- A quick search function applies to the search result and reduces the result according to the filter criteria.
- A filter provides various hierarchy levels. The selection of an item on a hierarchy level reduces the subsequent search results.
- The search result is available via a separate table.

OPEN DIALOG FOR PROTECTED OBJECTS (SR3146.9.2.20.1)

- For recent changes, see revision history (page 149).

For objects that are protected by an access privilege: The system shall only allow provide objects for selection whose access privilege matches the access privilege of the logged-in user.

Library Universe (SR3146.9.2.4+)

The Universe works as a library and allows to select data records for further processing in a resizable dialog.

The same library style is used for the following actions:

- Select a **material (SR3146.9.2.4.5)** data record (page 123) during the creation of a new master recipe.
- Select a **material (ERP BOM) (SR3146.9.2.4.2)** data record (page 123) during the creation of a new master recipe.
- Open a master recipe, master workflow, or custom building block by means of the **Open (SR3146.9.2.20+)** dialog (page 122).

LIBRARY FILTER & SEARCH (SR3146.9.2.4.3)

The Universe provides an intelligent filter and search function.

- A quick search function applies to the search result and reduces the result according to the filter criteria.
- A filter provides various hierarchy levels. The selection of an item on a hierarchy level reduces the subsequent search results.
- The search result is available via a separate table.
- Each item of the search result can be flagged as Setlist item in order to be available in the **Setlist (SR3146.9.2.9)** UI element (page [124](#)).

LIBRARY FOR BUILDING BLOCKS (SR3146.9.2.4.4)

- For recent changes, see revision history (page [149](#)).

The Universe provides building block-specific libraries to select procedures, unit procedures, operations, and phases.

For building block objects that are protected by an access privilege: The system shall only provide building block objects whose access privilege matches the access privilege of the logged-in user.

LIBRARY FOR MATERIALS (SR3146.9.2.4.5)

The Universe provides a material-specific library to select corresponding data records. The list of materials can be filtered per material type.

LIBRARY FOR MATERIALS (ERP BOM) (SR3146.9.2.4.2)

- For recent changes, see revision history (page [149](#)).
- Does not apply to Workflow Designer.

The workbench provides a material (ERP BOM)-specific library to select a product material for which an ERP BOM exists.

The library is only available during recipe creation in order to assign the product material to the recipe.

- For ERP BOMs that are protected by an access privilege: The system shall only allow to select product materials whose access privilege matches the access privilege of the logged-in user.
The access privilege of the ERP BOM is passed on to the master recipe. In the default configuration, the access privilege cannot be changed.
- For ERP BOMs that are not protected by an access privilege: The system shall allow to select an access privilege for the protection of the master recipe from unauthorized access for users with access privileges that allow to maintain protected recipes. The list of available access privileges shall be restricted to the access privileges of the logged-in user.

The library is not available in the **Universe (SR3146.9.2.4+)** UI element (page [122](#)).

LIBRARY FOR WORK CENTERS (SR3146.9.2.4.1)

The Universe provides a work center-specific library to select corresponding data records.

LIBRARY FOR STATIONS (SR3146.9.2.4.9)

The Universe provides a station-specific library to select corresponding data records.

LIBRARY FOR EQUIPMENT CLASSES (SR3146.9.2.4.7)

The Universe provides an equipment class-specific library to select corresponding data records. The list of equipment classes can be filtered per equipment class level.

LIBRARY FOR PROPERTY TYPES (SR3146.9.2.4.8)

The Universe provides a property type-specific library to select corresponding data records. The list of property types can be filtered per data type.

LIBRARY FOR SIGNATURE PRIVILEGES (SR3146.9.2.4.6)

The Universe provides a signature privilege-specific library to select corresponding data records.

Setlist (SR3146.9.2.9)

A Setlist provides pre-selected data records for building master recipes, master workflows, custom building blocks, and change requests.

- The data records (e.g. material, building blocks) are pre-selected in the **Universe (SR3146.9.2.4+)** UI element (page [122](#)).
- Capabilities are provided by the system.
- The Setlist is easily accessible during recipe, workflow, custom BB, and change requests building.
- Setlist items are added to the respective section of the recipe, workflow, and custom BB with a mouse click.
- Custom BBs are added to a change request with a mouse click.
- ERP BOM items are only available within Recipe Designer.

The Setlist provides functions for

- building sequences, selection branches, simultaneous branches, and joining elements according to the **procedural structure (SR3146.9.4)** hierarchy (page [13](#)),
- inverting the default inserting direction,

- (smart) replacing phases, operations, unit procedure, procedures,
- assigning materials and ERP BOM items as **process input (SR3146.9.5.1)** parameters (page 16) and **process output (SR3146.9.6.1)** parameters (page 24) to phases,
- assigning **work center (SR3146.9.5.4)** parameters (page 18) to unit procedures,
- assigning **station (SR3146.9.5.8)** parameters (page 18) to operations,
- assigning **equipment requirement (class) (SR3146.9.5.5)** parameters (page 19) as process input to phases,
- assigning **signature privilege (SR3146.9.5.3)** parameters (page 17) as process inputs to phases, and
- assigning **capability (SR3146.9.5.9.)** parameters (page 20) as process inputs to operations or unit procedures.

ERP BOM items are only available for master recipes based on an ERP BOM; in the **ERP BOM Item (SR3146.9.2.9.1)** section (page 125) of the Setlist.

SETLIST - ERP BOM ITEM (SR3146.9.2.9.1)

- Does not apply to Workflow Designer.

The **ERP BOM Item** section is filled automatically with BOM items (and their common attributes) of the ERP BOM on which a master recipe is based.

When an ERP BOM item is assigned to a recipe/building block element, the material data and common BOM item attributes are populated into the corresponding material input parameter.

Common attributes are position, material identifier and short description, quantity with UoM, and fixed quantity.

SETLIST - WORKFLOW-RELATED CONSISTENCY CHECKS (SR3146.9.2.9.2)

In Workflow Designer, the Setlist can hold custom building blocks that are not conform to the required master workflow structure.

- If the conflict cannot be solved within the master workflow, the system does not allow to add the custom building block to a master workflow (e.g.: a custom procedure holds more than one unit procedure).
- If the conflict can be solved within the master workflow, the system allows to add the custom building block to the master workflow and displays a corresponding error message in the **Messages (SR3146.9.2.16+)** window (page 128) (e.g. a custom operation holds a phase that is not marked for usage within a workflow).

SETLIST - WORKFLOW-RELATED DATA ADJUSTMENT (SR3146.9.2.9.3)

In Workflow Designer, the Setlist can hold custom building blocks with recipe-related attribute definitions. This applies to:

- **Dispense** attribute for operations/unit procedures (SR3146.9.4.4 operation type (page 15)).

If a custom building block is added to a master workflow, the recipe-related attribute definitions are ignored and removed. (They bear no meaning for workflows and cannot be edited in the context of a workflow.)

Parameter Panel (SR3146.9.2.7)

- For recent changes, see revision history (page 149).

A Parameter Panel allows to maintain existing parameters, transitions (with an identifier), and their attributes.

- If applicable, parameters can be deleted.
- Parameters and transitions can be locked and unlocked.
 - Locking of parameters (SR3146.9.13+) feature (page 27)
 - Locking of transitions (SR3146.9.14+) feature (page 27)
- The **Lock** attribute provides the following information:
 - Frozen: source building block is approved.
 - Locked: parameter/transition is locked.
 - Unlocked: parameter/transition is unlocked.
- Allows to filter per parameter type (e.g. process parameter, material parameter, transition).
- The Parameter Panel is resizable.
- Applies to material parameters:
 - In its vertical pane, the Parameter Panel displays each attribute of a parameter
 - In its horizontal pane, the Parameter Panel displays only the **Basic** attributes of a parameter.

Default attributes for each parameter are:

- Lock
- Hierarchy level
- Path

Parameters and transitions (with identifiers) are accessible via buttons on the building block images (page 117).

Cell editors facilitate easy editing of values: Measured Value editor (page 139), Multi-line text editor (page 139), Expression editor (SR3146.9.9.4) UI element (page 131), Date/Time Picker (SR3146.9.2.17.1) editor (page 139), Duration (SR3146.9.2.17.2) editor (page 139), Packaging Level Data (SR3146.9.2.17.3) editor (page 140), and further editors (page 140).

Materials, signature privileges, work center parameters, equipment requirement parameters, property types, capabilities, and dynamically created process parameters can be deleted as long as they are not locked by means of the Lock (SR3146.9.13+) attribute (page 27).

Applies to material parameters:
The vertical pane of the Parameter Panel is collapsible.

PARAMETER PANEL - VIEW ALL (SR3146.9.2.7.1)

The Parameter Panel can be accessed from all levels of the S88 building blocks. The visible parameters include all of the local parameters of the building block below the current building block level (view all per building block).

PARAMETER PANEL - SPECIFIC PROPERTY TYPES (SR3146.9.2.7.2)

Process input parameters can be refined by adding one or more parameter-specific property types and, if applicable, by defining a specific content for each property type.

This capability applies to equipment requirement parameters.

PARAMETER PANEL FILTER & SEARCH (SR3146.9.2.8)

The Parameter Panel provides an intelligent filter and search function.

- A quick search function applies to the search result and reduces the result according to the filter criteria.
- A filter provides various hierarchy levels. The selection of an item on a hierarchy level reduces the subsequent search results.
- The search result is available via a separate table.

Property Windows

- For recent changes, see revision history (page 149).

The workbench provides property windows specific to master recipe, master workflow, building block, and change request components.

Cell editors facilitate easy editing of values: Measured Value editor (page 139), Multi-line text editor (page 139), and Packaging Level Data (SR3146.9.2.17.3) editor (page 140).

Property window - Header (SR3146.9.2.15)

Provides the properties of the header component (master recipe, master workflow, custom building block, change request).

Property window - Element (SR3146.9.2.12)

Provides the properties of the element component within a master recipe, master workflow, or custom building block.

Property window - Source BB (SR3146.9.2.13)

Provides the properties of the source system/custom building block that was used as a template.

This applies to the usage of a BB as a master recipe, master workflow, or building block element, "Save as" of a BB, and "Create BB from selected".

Messages (SR3146.9.2.16+)

The **Messages** window displays messages of various categories related to the current master recipe, master workflow, or building block. They result from checks performed by the Designer workbench when master recipes, master workflows, or building blocks are being edited.

USER INTERFACE OF THE MESSAGES WINDOW

The tool panel displays the number of messages per message category and a breadcrumb trail that indicates the path to the element that is currently selected in the tree panel.

The tree panel displays a structural overview with all nodes of the master recipe, master workflow, or building block for which the checks have found messages. Navigating the tree highlights messages in the list panel and updates the breadcrumb trail.

The list panel displays a list of all messages that have been found for the master recipe, master workflow, or building block. Messages are prepended by unique message/message category identifiers.

Navigating the list highlights nodes in the tree panel, updates the breadcrumb trail, can select affected elements in the Graph Window, can select affected parameters in the Parameter Panel, can select the Material Flow Control tab, can select affected properties in the property windows, and can select expressions in the Expression editor of affected transitions, process parameter attributes, equipment requirement rules, if applicable.

The message category is indicated by an icon in the tool and the list panel.

The tab title also contains the number of messages per message category.

MESSAGES - VALIDATOR-RELATED (SR3146.9.2.16.1)

The messages are related to an ongoing validation (consistency checks) of the procedural structure and its parameters. The following message categories are supported:

- **Error**
Errors need to be resolved for approved custom building blocks and valid master recipes or master workflows.
- **Warning**
Warnings need to be resolved for valid master recipes or master workflows.
- **Information**
Validator result that does not require to be resolved.

For consistency check-related messages, see "Validator Consistency Checks (SR3146.9.4.3+)" (page 48).

MESSAGES - CATEGORY FILTER (SR3146.9.2.16.2)

The display of messages can be filtered by message category.

MESSAGES - NAVIGATION (SR3146.9.2.16.3)

The structural overview in the tree panel allows the navigation based on the tree structure of the master recipe, master workflow, or building block.

MESSAGES - RESOLVE ISSUE (SR3146.9.2.16.4)

- For recent changes, see revision history (page 149).

For messages caused by inconsistencies that can be resolved automatically, the system provides the **Resolve issue** function in the context menu of the related message.

The system allows to resolve a single issue or all issues at once.

The **Resolve issue** function is available for the following checks:

- **Validator - Phase parameter - ERP BOM attribute discrepancy (SR3146.9.4.3.56) check (page 55)**
- **Validator - Recipe - Planned quantity does not fit ERP BOM (SR3146.9.4.3.57) check (page 54)**
- **Validator - Recipe - Packaging Level Does Not Fit ERP BOM (SR3146.9.4.3.58) check (page 54)**
- **Validator - Material Parameter - Packaging Level Does Not Fit Material's Master Data (SR3146.9.4.3.109) check (page 60)**

Comparison (SR3146.9.2.19+)

The **Comparison** window displays the result of a comparison of the current master recipe, master workflow, or custom building block with the comparison baseline. See also "Comparing Master Recipes, Master Workflows, and Custom Building Blocks (SR3146.15+)" (page 79).

USER INTERFACE OF THE COMPARISON WINDOW

The tool panel displays the number of differences per difference category, an option list to define the objects whose differences are displayed, and a breadcrumb trail that indicates the path to the element that is currently selected in the tree panel.

The tree panel displays a structural overview with all nodes of the master recipe, master workflow, or custom building block for which the comparison has found differences. Navigating the tree highlights entries in the list panel and updates the breadcrumb trail.

The list panel displays a list of all differences that have been found for the master recipe, master workflow, or custom building block.

Navigating the list highlights nodes in the tree panel, updates the breadcrumb trail, can select affected elements in the Graph Window, can select affected parameters in the Parameter Panel, can select the Material Flow Control tab, can select affected properties in the property windows, and can select expressions in the Expression editor of affected transitions, process parameter attributes, equipment requirement rules, if applicable.

The difference category is indicated by an icon in the tool and the list panel.

The tab title also contains the number of differences per difference category.

Removed objects are displayed last under their parent nodes.

COMPARISON - DIFFERENCE CATEGORIES (SR3146.9.2.19.1)

The following difference categories are supported:

- Added
A parameter, parameter's child node (e.g. action of a capability), or procedural element has been added.
- Changed
An attribute, parameter, parameter's child node, procedural element, MFC-related setting, or the layout has been changed.
- Removed
A parameter, parameter's child node, or procedural element has been removed.

For the difference-specific handling, see "Comparison Algorithm (SR3146.15.1+)" (page 79).

COMPARISON - CATEGORY FILTER (SR3146.9.2.19.2)

The display of differences can be filtered by difference category.

COMPARISON - TYPE FILTER (SR3146.9.2.19.3)

The display of messages can be filtered by type.

COMPARISON - NAVIGATION (SR3146.9.2.19.4)

The structural overview in the tree panel allows the navigation based on the tree structure of the master recipe, master workflow, or custom building block.

COMPARISON - DIFFERENCE-SPECIFIC INFORMATION (SR3146.9.2.19.5)

The following difference-specific information is available:

- Type
- Identifier
- Predecessor
- Successor
- Attribute
- New value
- Old value

Preview Panel (SR3146.9.2.11)

The Preview panel consists of the Phase Preview which displays how the selected phase will be rendered when executed.

If the required information or input is not available yet, the preview shows a corresponding system message.

DEFAULT PREVIEW BUILDTIME (SR3146.9.9.4.1)

If a phase returns a preview-related error, the preview of the phase in the **Preview Panel (SR3146.9.2.11)** UI element (page [131](#)) is replaced by a default preview.

Example: Input parameters are not available.

Expression Editor (SR3146.9.9.4)

The Expression editor supports:

- References to output variables (e.g. used for information flow between phases)
- References to attributes of equipment properties (e.g. used for definition of equipment-related rules)

- Operators (see SR3146.9.9.4.4 (page 134))
- Functions (see SR3146.9.9.4.9 (page 138), SR3146.9.9.4.7 (page 135), SR3146.9.9.4.11 (page 135), and SR3146.9.9.6 (page 135))

The editor is available for **process parameter (SR3146.9.7+)** data (page 25), **rule (SR3146.9.5.7)** definitions (page 20), and **transition condition (SR3146.9.4.2)** definitions (page 15).

Intelligent auto-completion feature supports the selection of expression components.

The list of available expression elements displays indicators for specific elements with details as tooltip (e.g. parallel elements, event-triggered elements, successor elements).

Transition identifier of an expression is displayed within the transition in the Graph Window and as part of the respective tooltip.

Rule identifier of an expression is displayed within the process input requirement identifier.

Textual description of an expression is displayed as value in the parameter panel and as part of the tooltip for the transition in the Graph Window, respectively.

Concurrent syntax and semantics validation supports the definition of expressions and references.

Output references consist at least of a phase identifier and an output identifier. For references across operations and unit procedures, their identifiers are also components of the output reference.

Equipment property attribute references consist at least of an equipment reference, the property identifier, and the attribute identifier.

Blanks, tabs, and line breaks can be used for better readability.

Round brackets can be used to override or clarify an operators' precedence.

Single-line and multi-line comments can be added.

Indicators for errors (red underline of expression, error icon at affected line, global error status marker, error message as tooltip).

Code folding is supported to enable a clear structuring of large expressions (except for **Pnuts (SR3146.9.9.4.6)** scripting (page 138)).

EXPRESSIONS FOR TRANSITIONS (SR3146.9.9.4.2)

Conditions within a transition can be defined with the Expression editor.

In case a condition or a description is defined, the definition of a transition identifier is mandatory. Transition identifiers are generated by PharmaSuite and can be edited by the recipe or workflow author. They must be unique within their procedural element of the **procedural structure (SR3146.9.4)** hierarchy (page 13).

For transitions that only hold the default condition, identifiers are optional.

User-defined transition conditions apply in addition to the default condition of the **transition (SR3146.9.4.2)** condition (page 15).

The expression result must be of the **Boolean** type.

Default attributes for each transition are:

- Identifier
- Lock (SR3146.9.14+) attribute (page 27)
- Description
- Expression

EXPRESSIONS FOR PROCESS PARAMETERS (SR3146.9.9.4.3)

Process parameters can be defined with the Expression editor.

An editor-specific icon indicates the editor's availability.

The data type of the expression result must match the attribute's data type.

EXPRESSIONS FOR FLEXIBLE RULES (SR3146.9.9.4.10)

Flexible rules within a process input requirement can be defined with the Expression editor.

In case a flexible rule is added to a process input requirement, the definition of a rule identifier is mandatory. Rule identifiers are generated by PharmaSuite and can be edited by the recipe or workflow author. They must be unique within their phase.

The expression result must be of the **Boolean** type.

EXPRESSIONS FOR CONDITIONAL RULES (SR3146.9.9.4.12)

Conditional rules within a process input requirement can be defined with the Expression editor.

They are group-enabled, i.e. per default they apply as an equipment requirement to the identified parent entity and each of its child entities.

The Expression editor supports the following keywords to define a conditional rule:

- **provided that** allows to define conditions that must be fulfilled and
- **require** allows to define the checks themselves.

If the condition of the **provided that** section is not fulfilled for at least one entity of the group, the conditional rule fails.

In case a conditional rule is added to a process input requirement, the definition of a rule identifier is mandatory. Rule identifiers are generated by PharmaSuite and can be edited by the recipe or workflow author. They must be unique within their phase.

The expression result must be of the **Boolean** type.

EXPRESSION EDITOR OPERATIONS (SR3146.9.9.4.4)

The Operators tree panel provides the following operators:

- Arithmetic operators: + , - , * , /
- Comparison operators: ==, !=, <, <=, >, >=
- Logical operators: AND, OR, NOT
- Unary arithmetic operators: + , -
- String concatenation: +

A special handling applies to the Date and Duration data types:

- Date - Date = Duration
- Date (+, -) Duration = Date
- Duration (*, /) Scalar = Duration

A special handling applies to the Measured value (MV) data type:

- MV (*, /) Scalar = MV
- MV (+, -) MV = MV
- MV (/) MV = Scalar.

EXPRESSION EDITOR LITERALS (SR3146.9.9.4.5)

The following literals are supported:

- Long - integral numbers: e.g. 12345
- Float - floating point numbers (fractional numbers): e.g. 123.45Float
- BigDecimal - floating point numbers (fractional and integral numbers) with greater precision than Float: e.g. 123.45 or 123BigDecimal or 12.3e42
- Boolean: true, false
- String: e.g. "This is a string."
- Date and time: e.g. 07/19/2012 9:20:00 AM CEST
- Duration: e.g. 23d1h44min3s
- Measured value: e.g. 23kg

Additional literals:

- IMESS88Equipment: equipment data object
- PhaseDataReference collected phase data

EXPRESSION EDITOR FUNCTIONS (SR3146.9.9.4.7)

The Functions tree panel provides the following functions:

- Current Date and Time function: now()
- Convert to Unitless Number function: convertTo(arg1, arg2)
- Convert to BigDecimal function: convertToBigDecimal(arg1)
- Convert to Long function: convertToLong(arg1, arg2)
- Convert to MeasuredValue function: convertToMeasuredValue(arg1, [arg2])
- Convert to String for Display function: convertToDisplayString(arg1)
- String Contains function: containsItem(arg1, arg2, [arg3])
- Average function: average(arg1, arg2, ...)
- Minimum function: min(arg1, arg2, ...)
- Maximum function: max(arg1, arg2, ...)

- Null-handling functions: see SR3146.9.9.4.9 (page [138](#))
- Phase-specific functions: see SR3146.9.9.4.11 (page [135](#))
- Runtime context data functions: see SR3146.9.9.6 (page [135](#))
- Average function is only supported by Numeric data types.

EXPRESSION EDITOR - PHASE-SPECIFIC FUNCTIONS (SR3146.9.9.4.11)

The Functions tree panel provides the following function only if it is supported in the context of the specific phase (e.g. **Identify equipment** phase (SR0300+)):

- Reference to the currently identified equipment: currentEquipmentEntity (equipment context free)

EXPRESSION EDITOR - RUNTIME CONTEXT DATA (SR3146.9.9.6)

The Functions tree panel provides the following functions related to runtime context data:

- Order context
 - Order ID function: orderId

- Does not apply to Workflow Designer.
 - Order Scaling Factor function: orderScalingFactor
 - Produced Material ID function: producedMaterialId
 - Produced Material Short Description function: producedMaterialShortDescription
 - Produced Material Type function: producedMaterialType
 - Batch ID function: batchId
 - Batch Status function: batchStatus
- MFC position context

The functions of the MFC position context are only supported for Dispense, Inline Weighing, and Output Weighing; they are not supported by the Material Tracking phases (**Identify material, Produce material, Account material**).

 - Material ID function: materialId
 - Material Short Description function: materialShortDescription
 - Material Type function: materialType
 - Status function: mfcPositionStatus
 - Planned Quantity - Original function: originalPlannedQuantity
 - Planned Quantity - Execution function: plannedQuantity
 - Actual Quantity - Last split function: actualSplitQuantity
Returns the actual recorded quantity of all the last split positions that constituted the last target of a material input position. For material outputs, the function always returns null.
 - Actual Quantity - Total function: totalActualQuantity
Returns the total recorded quantity of a material position minus replaced quantity or quantity that was declared as waste.
 - Sum of Actual Input Quantities - Total function: sumTotalActualInputQuantity
Returns the total recorded quantity of all material input positions of the unit procedure or the specified material input positions of the procedure.
 - Prorate Factor function: prorateFactor
Does not apply to Recipe Designer - Device and Workflow Designer.
 - Yield function: yield
Does not apply to Recipe Designer - Device and Workflow Designer.
- Equipment context
 - Current equipment entity (runtime-related) function: currentEquipmentEntity

- Equipment is member of class function: `equipmentIsMemberOfClass(arg1, arg2)`
With an equipment object as first argument and an equipment class identifier as second argument.
- Equipment property attribute (runtime-related) function:
`equipmentPropertyAttribute(arg1, arg2 [, arg3, arg4])`
With an equipment object as first argument, a property type identifier as second argument, a property type attribute identifier as optional third argument, and an index as optional fourth argument.
- Equipment has property attribute function:
`equipmentHasPropertyAttribute(arg1, arg2 [, arg3, arg4])`
With an equipment object as first argument, a property type identifier as second argument, a property type attribute identifier as optional third argument, and an index as optional fourth argument.
- Equipment graph status (runtime-related) function:
`equipmentGraphStatus(arg1, arg2 [, arg3])`
With an equipment object as first argument, a graph purpose as second argument, and an attribute identifier (KEY or DISPLAY) as optional third argument.
- Equipment has graph function: `equipmentHasGraph(arg1, arg2)`
With an equipment object as first argument and a graph purpose as second argument.
- Equipment graph expiry date (runtime-related) function:
`equipmentGraphExpiryDate(arg1, arg2)`
With an equipment object as first argument and a graph purpose as second argument.
- Equipment FSM status function: `equipmentFSMStatus(arg1, arg2 [, arg3])`
With an equipment object as first argument, the identifier of a FlexibleStateModel property type as second argument, and as optional third argument either the "KEY" string or the "DISPLAY" string.
- Date and time of last equipment status change (logbook-related) function:
`lastEqTransitionTimestamp(arg1, arg2, arg3[, arg4])`
With an equipment object as first argument, the identifier of a FlexibleStateModel property type as second argument, the non-localized state name or the localized status of the entity as third argument, and as optional fourth argument either the "KEY" string or the "DISPLAY" string.
- Last product ID (logbook-related) function: `lastProduct`
- Location context
- Work Center ID function: `workCenterId`

- Work Center Description function: workCenterDescription
- Room ID function: roomId
- Room Status function: roomStatus
- Storage Area ID function: storageAreaId
- Storage Area Description function: storageAreaDescription

EXPRESSION EDITOR - IMPLICIT NULL HANDLING (SR3146.9.9.4.8)

The Expression editor implicitly handles NULL markers according to the SQL standard. This leads to the following behavior for expressions including NULL:

- Mathematical expressions result in NULL
Example: $3 + \text{NULL} = \text{NULL}$
- Aggregation functions eliminate NULL
Example: $\text{MIN}(3, \text{NULL}, 5) = 3$
- Comparative expressions result in NULL
Example: $4 > \text{NULL}$ results in NULL
- Boolean expressions follow the Three-valued logic (3VL).

For transitions, NULL is mapped to FALSE.

EXPRESSION EDITOR - EXPLICIT NULL HANDLING (SR3146.9.9.4.9)

The Functions tree panel provides the following functions for the explicit handling of NULL markers:

- Is undefined: isNull(arg1)
- Undefined - one of two: nvl(arg1, arg2)
- Undefined - one of many: coalesce(arg1, arg2, ...)

EXPRESSION EDITOR - PNUITS (SR3146.9.9.4.6)

The Expression editor provides Pnuts scripting.

The Pnuts mode does not support extensive syntax and semantics validation, access to the Operators and Functions tree panels, and code folding.

Other Editors (SR3146.9.2.17+)

The workbench provides various editors that support data entry. This applies to the various types of input data.

MEASURED VALUE EDITOR

The system provides a Measured Value editor for entering a value with a unit of measure.

- An editor-specific icon indicates the editor's availability.
- The editor is available in the **Parameter Panel (SR3146.9.2.7)** UI element (page [126](#)) and the property windows (page [127](#)).
- See also **Measured Value Technical Property Type (SR3071.8.7.14)** in "Functional Requirement Specification Data Management" [A1] (page [147](#)).

MULTI-LINE TEXT EDITOR

The system provides a Multi-line text editor for writing text.

- The editor can interpret HTML syntax.
The system supports the following HTML tags: ,
, , <u>, <i>, <sup>, <sub>, , , and .
Text that is formatted with HTML is displayed as plain text including the HTML tags in master recipe reports or master workflow reports. In batch reports, DHR reports, and workflow reports, however, formatting with the supported tags is interpreted and applied.
- An editor-specific icon indicates the editor's availability.
- The editor is available in the **Parameter Panel (SR3146.9.2.7)** UI element (page [126](#)) and the property windows (page [127](#)).

DATE/TIME PICKER EDITOR (SR3146.9.2.17.1)

The system provides a Date/Time Picker editor for entering a date and time value.

- An editor-specific icon indicates the editor's availability.
- The editor is available in the **Parameter Panel (SR3146.9.2.7)** UI element (page [126](#)) and the Change status dialogs.

DURATION EDITOR (SR3146.9.2.17.2)

The system provides a Duration editor for entering a duration value.

- An editor-specific icon indicates the editor's availability.
- The editor is available in the **Parameter Panel (SR3146.9.2.7)** UI element (page [126](#)) and the property windows.

- See also **Duration Technical Property Type (SR3071.8.7.16)** in "Functional Requirement Specification Data Management" [A1] (page [147](#)).

PACKAGING LEVEL DATA EDITOR (SR3146.9.2.17.3)

- For recent changes, see revision history (page [149](#)).

The system provides a Packaging Level Data editor for entering the packaging level data: Meaning, Contained number, Inventory level, and a Hide during execution option.

- An editor-specific icon indicates the editor's availability.
- The editor is available in the **Parameter Panel (SR3146.9.2.7)** UI element (page [126](#)) and the **Property Window - Header (SR3146.9.2.15)** UI element (page [128](#)).

FURTHER EDITORS

The editors listed in this section are specified in the "Functional Requirement Specification Data Management" [A1] (page [147](#)).

- **BigDecimal editor**
(see **BigDecimal Technical Property Type (SR3071.8.7.3)**)
- **Boolean editor**
(see **Boolean Technical Property Type (SR3071.8.7.4)**)
- **Expression editor**
(see **Expression Editor (Data Manager) (SR3071.8.1.4)**)
- **FlexibleAttributeDefinition editor**
(see **Flexible Attribute Definition Technical Property Type (SR3071.8.7.9)**)
- **FlexibleStateModel editor**
(see **Status Model Technical Property Type (SR3071.8.7.6)**)
- **FlexibleTagDefinition editor**
(see **Flexible Tag Definition Technical Property Type (SR3071.8.7.8)**)
- **Label Layout Selection editor**
(see **Label Layout Selection Editor (SR3071.8.8.1)**)
- **String editor**
(see **String Technical Property Type (SR3071.8.7.5)**)

Statistics (SR3146.9.2.18)

The system provides a statistical summary for recipes, workflows, and custom building blocks, which includes the number of used artifacts of a selected object.

STATISTICS VIEW (SR3146.9.2.18.1)

The statistics view provides details of the following artifacts:

- Unit procedures
- Operations
- Phases
- Process parameters
- Material parameters
- MFC items (with number of inputs, outputs, transfers)
- Privilege parameters
- Capability parameters
- Work center assignment parameters
- Equipment requirement parameters
- Additional property requirements
- Transitions with identifiers
- Process parameter input expressions

STATISTICS-BASED WARNING (SR3146.9.2.18.2)

When a recipe, workflow, or custom building block with more than 1,000 phases is saved, the system displays a warning.

A recipe, workflow, or custom building block of this magnitude exceeds the boundaries of PharmaSuite's design space and can lead to unexpected system responses.

Version Control (SR3146.4.5)

The following objects are under version control:

- master recipes,
- workflows, and
- change requests.

For the available statuses of master recipes and master workflows, see

- "Standard Statuses for Version Control (SR3146.6.4)" (page [43](#))

For the available statuses of change requests, see

- "Status Management of Change Requests (SR3146.12.2+)" (page [85](#))

Status Change (SR3146.4.5.1)

Recipe and Workflow Designer provides a status change function for master recipes, master workflows, and change requests. This includes:

- Definition of the desired transition (target status).
- Definition of the period of validity.
Does not apply to change requests.
- Overview of other versions of the same master recipe or master workflow, incl. their statuses and periods of validity.
Does not apply to change requests.

PERIOD OF VALIDITY (SR3146.6.4.2)

The period of validity of a valid object includes the **Valid from...** and **Valid until...** dates.

The start date of the validity period includes the entire day and thus begins at 12:00 AM while the end date of the validity period includes the entire day and thus lasts until 11:59 PM. This means that a validity period of exactly one day can be defined by selecting the same date for both start and end dates.

CONSISTENCY CHECKS (SR3146.6.4.3)

The status change checks verify the consistency of a desired status change and its related period of validity. In case of conflicts, corresponding messages are displayed.

Version History (SR3146.4.5.3)

Recipe and Workflow Designer provides a version history function for master recipes, master workflows, and change requests in order to track status transitions. For master recipes and master workflows, the version history function also tracks version transitions. The history includes transitions that took place in a different environment (e.g. prior to export or import of a master recipe or master workflow).

Configurable Version Control (SR3146.6)

The version control is configurable. The configuration includes:

- Definition of version graphs with statuses and transitions.
- Assignment of an electronic signature required to perform a status change.
Electronic signatures are configurable.
- Definition of a numbering scheme for versions.

Example:

- Initial version = **V1.0**
- Status change from **Verification** to **Edit**: new version = **V1.1**
- **V1.1** is copied: new version = **V2.0**

Flexible Version Graph Assignment (SR3146.6.1)

A status graph can be assigned per class.

Example: Class Master recipe - Standard graph

Flexible Definition of Approval Details (SR3146.6.7)

The details for status changes and approvals can be defined.

For master recipes and master workflows, this includes an e-mail notification mechanism on user level.

FLEXIBILITY OF USER GROUPS (SR3146.6.7.1)

The number of user groups available in the system defines the number of different user groups that can be assigned for approval.

For further details, see use cases and RQ of electronic signatures (Security module).

SERIAL APPROVAL (SR3146.6.7.3)

Serial approval of n user groups/users can be enforced.

Example 1: Group A approves, then group B can approve

For further details, see use cases and RQ of electronic signatures.

Library Concept (SR3146.9.11+)

The library concept for building blocks according to S88 includes Basic Building Blocks and Custom Building Blocks.

Basic Building Block Import (SR3146.9.11.1)

Basic building block libraries can be imported.

Custom Building Block (SR3146.9.11.3)

Building blocks can be saved as custom building blocks. This applies to

- phases,
- operations,
- unit procedures, and
- procedure-level building blocks.

-
-
- Rockwell Software PharmaSuite® 8.4 - Functional Requirement Specification Recipe and Workflow Management
-
-

Export/Import of Master Recipes, Master Workflows, and Building Blocks (SR1075.1+)

Developing a master recipe or master workflow is an ambitious task and should be carefully planned. Thorough testing is essential before they can be approved. It is a good practice to manage their development life-cycle across multiple PharmaSuite instances (e.g. development system, test system, productive system). For this purpose, master recipes, master workflows, and building blocks need to be transferred from a source system to a target system.

General Prerequisites and Usage of B2MML (SR1075.1.1)

- For recent changes, see revision history (page [149](#)).

Based on B2MML, the **PharmaSuite Export/Import Utility** tool allows to export and imports master recipes, master workflows, and building blocks maintained with Recipe and Workflow Designer.

For objects that are protected by an access privilege: The system shall only allow to export objects whose access privilege matches the access privilege of the logged-in user.

The following prerequisites apply:

- Both the source and target system must be of the same (maintenance) release.
- Available system building blocks must be identical in both the source and target system.
- Any master data object that is referenced in the source system by the master recipe, master workflow, or building block to be exported must be available in the target system with the same attributes (e.g. identifier, description, version).
- The versions of the underlying B2MML schemas must be identical.

Prerequisites For Import (SR1075.1.1.1)

Given the exported data, the **PharmaSuite Import Utility** tool shall check if the target system fulfills the prerequisites for the import. In case of potential conflicts, the import shall abort and the conflicts shall be recorded.

Checksum (SR1075.1.1.2)

The **PharmaSuite Import Utility** tool shall detect if the exported data file has been modified.

Handling of Status and Version (SR1075.1.1.3)

The **PharmaSuite Export/Import Utility** tool shall assure that in the target system the status and the version of the imported objects are identical compared to the source system.

Handling of Status Transition History (SR1075.1.1.4)

For a given imported master recipe, master workflow, or building block, the system shall provide information about the source system that holds the history data of the latest status transition prior to the export/import.

Reference Documents

The following documents are available from the Rockwell Automation Download Site.

No.	Document Title	Part Number
A1	PharmaSuite Functional Requirement Specification Data Management	PSFRSDM-RM004E-EN-E
A2	PharmaSuite Functional Requirement Specification Execution Framework	PSFRSEF-RM004E-EN-E
A3	PharmaSuite Functional Requirement Specification Runtime Data Management	PSFRSRT-RM004E-EN-E
A4	PharmaSuite Functional Requirement Specification Non-functional Requirements	PSFRSNF-RM001E-EN-E

TIP

To access the Rockwell Automation Download Site, you need to acquire a user account from Rockwell Automation Sales or Support.

-
-
- Rockwell Software PharmaSuite® 8.4 - Functional Requirement Specification Recipe and Workflow Management
-
-

Document Information

The document information covers various data related to the document.

Approval

This document has been approved electronically via the Rockwell Automation Document Management System (DMS). The required approvers of this document include the following:

Name	Role
Martin Dittmer	Product Manager
Steffen Landes	Development Manager
Martin Irmisch	Test Manager

In addition, the electronic document approval via DMS is confirmed by a handwritten signature of all approvers in the Quality Document when the release is completed. The Quality Document summarizes the quality-related planning activities and results of a PharmaSuite release.

Version Information

Object	Version
PharmaSuite	8.4
Functional Requirement Specification	1.1

Revision History

The following table describes the history of this document.

Changes related to the document:

Object	Description	Document
---	---	---

Changes related to "Recipe and Workflow Management" (page 3):

Object	Description	Document
Recipe and Workflow Management (SR3146+) (page 3)	Context information-related PharmaSuite supports the concept of confidential objects to protect the intellectual property of recipes, workflows, orders, and related data from unauthorized access.	1.0

Changes related to "Management of Master Recipes and Master Workflows" (page 7):

Object	Description	Document
Master Recipe - Header Attributes (SR3146.9.3.1) (page 9)	Update New attributes (L6 (Meaning, Contained number, Inventory level, Hide during execution), ..., L9 (Meaning, Contained number, Inventory level, Hide during execution)) to define packaging levels on recipe header level. New attribute (Access privilege): PharmaSuite supports the concept of confidential objects to protect the intellectual property of recipes, workflows, orders, and related data from unauthorized access.	1.0
Validator - XOR Branch (SR3146.9.4.3.2)	Deleted The validator consistency check is no longer applicable since selection branches (XOR) between unit procedures of master recipes are allowed.	1.0
Validator - Loop (SR3146.9.4.3.3) (page 48)	New The validator consistency check does not allow loops between unit procedures of a master recipe or (custom) building block.	1.0
Validator - Phase Parameter - Multiple ERP BOM Positions (SR3146.9.4.3.51) (page 55)	Update The validator consistency check also applies to procedure building blocks. In any case, the positions are not unique throughout the entire execution path.	1.0
Validator - Master Recipe and Material Parameter - Conflicting Inventory Levels (SR3146.9.4.3.52) (page 53)	New Validator consistency check to verify the definition of packaging levels.	1.0
Validator - Material Parameter - Conflicting Inventory Levels (SR3146.9.4.3.53) (page 59)	New Validator consistency check to verify the definition of packaging levels.	1.0
Validator - Recipe - Packaging Level Does Not Fit ERP BOM (SR3146.9.4.3.58) (page 54)	New Validator consistency check to verify the definition of packaging levels.	1.0

Object	Description	Document
Validator - Recipe - Packaging Level Does Not Fit ERP BOM - Resolve (SR3146.9.4.3.58.1) (page 54)	New Validator consistency check to verify the definition of packaging levels.	1.0
Validator - MFC - No Output Defined for Input (SR3146.9.4.3.63) (page 51)	Update The validator consistency check applies also to input switches.	1.0
Validator - MFC - No Input Defined for Output (SR3146.9.4.3.64) (page 51)	Update The validator consistency check applies also to output switches.	1.0
Validator - MFC - Uniqueness and Format of MFC Position (SR3146.9.4.3.68) (page 57)	Update The validator consistency check applies to an entire execution path.	1.0
Validator - Material - Incompatible UoM (SR3146.9.4.3.69) (page 57)	Update The validator consistency check also applies to relative tolerances. No change of code. The validator consistency check also applies to the target weight and its tolerances.	1.0
Validator - Material - One Tolerance Is Missing (SR3146.9.4.3.71) (page 57)	Update The validator consistency check also applies if a target weight has been defined and either its lower or upper tolerance is missing.	1.0
Validator - D_Material - Tolerances Not Defined (SR3146.9.4.3.74) (page 62)	Update The validator consistency check also applies if a target weight has been defined and no tolerances have been defined.	1.0
Validator - Material - Negative Values (SR3146.9.4.3.77) (page 58)	Update The validator consistency check also applies if a target weight has been defined and its quantity or tolerance is negative.	1.0
Validator - Packaging Level of Master Recipe (SR3146.9.4.3.107) (page 53)	New Validator consistency check to verify the definition of packaging levels.	1.0
Validator - Packaging Level of Material Parameter (SR3146.9.4.3.108) (page 59)	New Validator consistency check to verify the definition of packaging levels.	1.0
Validator - Material Parameter - Packaging Level Does Not Fit Material's Master Data (SR3146.9.4.3.109) (page 60)	New Validator consistency check to verify the definition of packaging levels.	1.0

Object	Description	Document
Validator - Material Parameter - Packaging Level Does Not Fit Material's Master Data - Resolve (SR3146.9.4.3.109.1) (page 60)	New Validator consistency check to verify the definition of packaging levels.	1.0
Validator - MFC - No Direct Output Defined (SR3146.9.4.3.117) (page 52)	Update The validator consistency check also applies if at least for one execution path no direct output is defined.	1.0
Validator - MFC - No Direct Output Defined - Warning (SR3146.9.4.3.117.1) (page 52)	Update The validator consistency check also applies if at least for one execution path no direct output is defined.	1.0
Validator - MFC - No Direct Output Defined - Error (SR3146.9.4.3.117.2) (page 52)	Update The validator consistency check also applies if at least for one execution path no direct output is defined.	1.0
Validator - MFC - Attributes of Material Parameters Do Not Match (SR3146.9.4.3.136) (page 52)	New Validator consistency check to verify the attributes of material parameters with the same MFC position.	1.0
Validator - BOM Position Not Used (SR3146.9.4.3.137) (page 58)	New Validator consistency check to verify the usage of BOM positions.	1.0
Validator - MFC - Outgoing or Incoming Connector Is Missing (SR3146.9.4.3.138) (page 52)	New Validator consistency check to verify the available outgoing and incoming connectors of switches.	1.0
Validator - MFC - More Than One Outgoing or Incoming Connector on One Execution Path (SR3146.9.4.3.139) (page 53)	New Validator consistency check to verify the available outgoing and incoming connectors of switches.	1.0
Process Inputs - Material Attributes (SR3146.9.5.1.1) (page 16)	Update New attributes to define a target weight and its tolerances (lower and upper, absolute and relative). New attributes (L0 (Meaning, Contained number, Inventory level, Hide during execution), ..., L5 (Meaning, Contained number, Inventory level, Hide during execution)) to define packaging levels on material parameter level.	1.0

Object	Description	Document
Process Outputs - Material Attributes (SR3146.9.6.1.1) (page 24)	<p>Update</p> <p>The attributes are arranged in the Basic data and Weighing groups to support its display in the horizontal and vertical panes of the Parameter Panel (SR3146.9.2.7) UI element. No change of code.</p> <p>New attributes to define a target weight and its tolerances (lower and upper, absolute and relative).</p> <p>New attributes (L0 (Meaning, Contained number, Inventory level, Hide during execution), ..., L5 (Meaning, Contained number, Inventory number, Hide during execution)) to define packaging levels on material parameter level.</p>	1.0
Graphical View (SR3146.9.8.1) (page 37)	<p>Context information-related</p> <p>Representation of switch nodes that are created when MFC inputs and outputs are merged.</p> <p>A double-click of an input switch, output switch, confluence, or connector shows the related material parameter(s) within the Parameter Panel (SR3146.9.2.7) UI element.</p> <p>Connectors represent the material flows between the switch nodes. Connectors between switch nodes are dashed lines if they represent alternative transfers.</p>	1.0
Tooltips (SR3146.9.8.2) (page 39)	<p>Update</p> <p>Tooltips are available for switches.</p>	1.0
Automatic Merge - Default MFC Definition (SR3146.9.8.3) (page 39)	<p>Editorial</p> <p>For clarification, Automatic Merge added to requirement name.</p>	1.0
Default MFC Definition between Unit Procedures (SR3146.9.8.5) (page 40)	<p>Update</p> <p>In case of selection branches (XOR), the merge results in alternative transfers.</p>	1.0
Manual Merge between Unit Procedures (SR3146.9.8.7) (page 41)	<p>Update</p> <p>In case of selection branches (XOR), one input can be merged into multiple alternative outputs and one output can be merged into multiple alternative inputs. In both cases alternative transfers are created.</p>	1.0
Split (SR3146.9.8.8) (page 42)	<p>Context information-related</p> <p>Split of input and output switches added.</p>	1.0
Workbench Object Operations - Master Recipe (SR3146.9.10.1) (page 28)	<p>Update</p> <p>Copy (Save as ...), Create new procedure: PharmaSuite supports the concept of confidential objects to protect the intellectual property of recipes, workflows, orders, and related data from unauthorized access.</p>	1.0

Object	Description	Document
Workbench Object Operations - New Master Recipe Dialog (SR3146.9.10.1.2) (page 29)	Update PharmaSuite supports the concept of confidential objects to protect the intellectual property of recipes, workflows, orders, and related data from unauthorized access.	1.0
Workbench Object Operations - Element (SR3146.9.10.6) (page 31)	Update PharmaSuite supports the concept of confidential objects to protect the intellectual property of recipes, workflows, orders, and related data from unauthorized access.	1.0
Workbench Object Operations - Master Workflow (SR3146.9.10.9) (page 29)	Update Copy (Save as ...), Create new procedure: PharmaSuite supports the concept of confidential objects to protect the intellectual property of recipes, workflows, orders, and related data from unauthorized access.	1.0
Workbench Object Operations - New Master Workflow Dialog (SR3146.9.10.9.1) (page 30)	Update PharmaSuite supports the concept of confidential objects to protect the intellectual property of recipes, workflows, orders, and related data from unauthorized access.	1.0
Workbench Object Operations - Replace Material (Product) (SR3146.9.10.11) (page 36)	Update Along with replacing the material (product), the system also updates the packaging levels of the materials and, in case the new material (product) is related to an ERP BOM, also the packaging levels of the master recipe. PharmaSuite supports the concept of confidential objects to protect the intellectual property of recipes, workflows, orders, and related data from unauthorized access.	1.0
Master Workflow - Header Attributes (SR3146.9.12.1) (page 11)	Update New attribute (Access privilege): PharmaSuite supports the concept of confidential objects to protect the intellectual property of recipes, workflows, orders, and related data from unauthorized access.	1.0
Workbench Object Operations - New Master Recipe Dialog (SR3146.9.10.1.2) (page 29)	Update The list of available materials shall be restricted by the access privileges of the logged-in user.	1.1
Workbench Object Operations - Replace Material (Product) (SR3146.9.10.11) (page 36)	Update The list of available materials shall be restricted by the access privileges of the logged-in user. If an access privilege is defined for the new ERP BOM or material that differs from the access privilege of the master recipe, the new access privilege is passed on to the master recipe.	1.1

Changes related to "Management of Building Blocks" (page 73):

Object	Description	Document
Usage List (SR3146.10.3) (page 76)	Update PharmaSuite supports the concept of confidential objects to protect the intellectual property of recipes, workflows, orders, and related data from unauthorized access.	1.0
Objects Information (SR3146.10.3.2) (page 77)	Update The usage list provides a list of changed objects whose objects information cannot be refreshed. No change of code.	1.0
Header Attributes (SR3146.10.4) (page 73)	Update New attribute (Access privilege): PharmaSuite supports the concept of confidential objects to protect the intellectual property of recipes, workflows, orders, and related data from unauthorized access.	1.0
Workbench Object Operations (SR3146.9.10.7) (page 74)	Update Save as new custom building block: PharmaSuite supports the concept of confidential objects to protect the intellectual property of recipes, workflows, orders, and related data from unauthorized access.	1.0
Workbench Object Operations - New Building Block Dialog (SR3146.9.10.7.2) (page 75)	New PharmaSuite supports the concept of confidential objects to protect the intellectual property of recipes, workflows, orders, and related data from unauthorized access. For the remaining content of the requirement applies: No change of code.	1.0

Changes related to "Comparing Master Recipes, Master Workflows, and Custom Building Blocks" (page 79):

Object	Description	Document
---	---	---

Changes related to "Management of Change Requests" (page 83):

Object	Description	Document
Header Attributes (SR3146.12.1) (page 83)	Update New attributes (Access privilege): PharmaSuite supports the concept of confidential objects to protect the intellectual property of recipes, workflows, orders, and related data from unauthorized access.	1.0

Object	Description	Document
Define Old/New Building Blocks (SR3146.12.3.1) (page 89)	Update PharmaSuite supports the concept of confidential objects to protect the intellectual property of recipes, workflows, orders, and related data from unauthorized access.	1.0
Select from Usage List (SR3146.12.3.2) (page 90)	Update PharmaSuite supports the concept of confidential objects to protect the intellectual property of recipes, workflows, orders, and related data from unauthorized access.	1.0
Execute Change Request (SR3146.12.5) (page 91)	Editorial For clarification, "Old versions of affected master recipes or master workflows are not changed, i.e. the status remains unchanged." added to the "The following descriptions reflect that the Target status attribute has been set to All valid statuses ." rule.	1.0
Open Dialog for Protected Objects (SR3146.12.6.2) (page 84)	New PharmaSuite supports the concept of confidential objects to protect the intellectual property of recipes, workflows, orders, and related data from unauthorized access.	1.0

Changes related to "Master Recipe Report - Batch" (page 95):

Object	Description	Document
MR Report - Basic Data (SR3146.11.4.1) (page 96)	Update PharmaSuite supports the concept of confidential objects to protect the intellectual property of recipes, workflows, orders, and related data from unauthorized access. PharmaSuite supports packaging levels.	1.0
MR Report - Procedure (SR3146.11.4.11) (page 99)	Update The report contains the list of transition between unit procedures that are modeled as exclusive alternative sequences (XOR).	1.0

Changes related to "Master Recipe Report - Device" (page 103):

Object	Description	Document
MR Report - Basic Data (SR3146.14.4.1) (page 104)	Update PharmaSuite supports the concept of confidential objects to protect the intellectual property of recipes, workflows, orders, and related data from unauthorized access. PharmaSuite supports packaging levels.	1.0

Object	Description	Document
MR Report - Procedure (SR3146.14.4.11) (page 107)	Update The report contains the list of transition between unit procedures that are modeled as exclusive alternative sequences (XOR).	1.0

Changes related to "Master Workflow Report" (page 111):

Object	Description	Document
MWF Report - Basic Data (SR3146.13.4.1) (page 112)	Update PharmaSuite supports the concept of confidential objects to protect the intellectual property of recipes, workflows, orders, and related data from unauthorized access.	1.0

Changes related to "Designer Workbench" (page 117):

Object	Description	Document
Library for Materials (ERP BOM) (SR3146.9.2.4.2) (page 123)	Update PharmaSuite supports the concept of confidential objects to protect the intellectual property of recipes, workflows, orders, and related data from unauthorized access.	1.0
Library for Building Blocks (SR3146.9.2.4.4) (page 123)	Update PharmaSuite supports the concept of confidential objects to protect the intellectual property of recipes, workflows, orders, and related data from unauthorized access.	1.0
Parameter Panel (SR3146.9.2.7) (page 126)	Context information-related PharmaSuite provides an editor specific to packaging levels.	1.0
Property Windows (page 127)	Context information-related PharmaSuite provides an editor specific to packaging levels.	1.0
Messages - Resolve Issue (SR3146.9.2.16.4) (page 129)	Context information-related The Resolve issue function is also available for the Validator - Recipe - Packaging level does not fit ERP BOM (SR3146.9.4.3.58) and Validator - Material Parameter - Packaging Level Does Not Fit Material's Master Data (SR3146.9.4.3.109) checks.	1.0
Packaging Level Data Editor (SR3146.9.2.17.3) (page 140)	New PharmaSuite provides an editor specific to packaging levels.	1.0
Open Dialog (SR3146.9.2.20+) (page 122)	New New grouping requirement. No change of code.	1.0

Object	Description	Document
Open Dialog for Protected Objects (SR3146.9.2.20.1) (page 122)	New PharmaSuite supports the concept of confidential objects to protect the intellectual property of recipes, workflows, orders, and related data from unauthorized access.	1.0

Changes related to "Export/Import of Master Recipes, Master Workflows, and Building Blocks" (page [145](#)):

Object	Description	Document
General Prerequisites and Usage of B2MML (SR1075.1.1) (page 145)	Update PharmaSuite supports the concept of confidential objects to protect the intellectual property of recipes, workflows, orders, and related data from unauthorized access.	1.0

B

BigDecimal editor • 140
 Boolean editor • 140
 Building block image • 117

C

Checks that are disabled for master recipes of the Cost center usage type • 65
 Comparing master recipes, master workflows, and custom building blocks (SR3146.15+) • 79
 Capability parameter-specific attributes (SR3146.15.1.4.6) • 81
 Comparing element attributes (SR3146.15.1.2) • 80
 Comparing header attributes (SR3146.15.1.1+) • 79
 Comparing parameter attributes (SR3146.15.1.4+) • 80
 Comparing procedural elements (SR3146.15.1.5+) • 81
 Comparing source attributes (SR3146.15.1.3) • 80
 Comparison algorithm (SR3146.15.1+) • 79
 Custom building block-specific header attributes (SR3146.15.1.1.3) • 80
 Defining a comparison baseline (SR3146.15.2) • 82
 Equipment requirement parameter-specific attributes (SR3146.15.1.4.3) • 81
 Location-specific characteristics (SR3146.15.1.5.3) • 82
 Master recipe-specific header attributes (SR3146.15.1.1.1) • 79
 Master workflow-specific header attributes (SR3146.15.1.1.2) • 79
 Material Flow Control-specific attributes (SR3146.15.1.4.2) • 80
 Material parameter-specific attributes (SR3146.15.1.4.1) • 80
 Performing a comparison (SR3146.15.3) • 82
 Privilege parameter-specific attributes (SR3146.15.1.4.5) • 81
 Process parameter-specific attributes (SR3146.15.1.4.7) • 81
 Relation-specific characteristics (SR3146.15.1.5.1) • 81
 Transition-specific attributes (SR3146.15.1.5.2) • 81
 Work center assignment parameter-specific attributes (SR3146.15.1.4.4) • 81
 Comparison window (user interface) • 130
 Compliance-related
 SR1075.1.1 - General prerequisites and usage of B2MML • 145
 SR1075.1.1.1 - Prerequisites for import • 145
 SR1075.1.1.2 - Checksum • 145
 SR1075.1.1.3 - Handling of status and version • 146
 SR1075.1.1.4 - Handling of status transition history • 146
 SR3071.6.2 - Access control • 4
 SR3071.6+ - Recipe Designer - Batch • 4
 SR3071.7.2 - Access control • 5
 SR3071.7+ - Workflow Designer • 5
 SR3071.9.2 - Access control • 5
 SR3071.9+ - Recipe Designer - Device • 5
 SR3146.10.1 - Status handling of building blocks • 75
 SR3146.10.4 - Header attributes • 73
 SR3146.11.4.3 - MR report - Approval record • 97
 SR3146.12.1 - Header attributes • 83
 SR3146.12.2.2 - Master recipe-specific change requests • 85
 SR3146.12.2.3 - Master workflow-specific change requests • 87
 SR3146.12.2+ - Status management of change requests • 85
 SR3146.13.4.3 - MWF report - Approval record • 113
 SR3146.14.4.3 - MR report - Approval record • 105
 SR3146.6 - Configurable version control • 142
 SR3146.6.4 - Standard statuses for version control • 43
 SR3146.6.4.1 - Default version graph for master recipes • 44
 SR3146.6.4.4 - Default version graph for master workflows • 46

SR3146.6.7 - Flexible definition of approval details • 143
SR3146.6.7.1 - Flexibility of user groups • 143
SR3146.6.7.3 - Serial approval • 143
SR3146.9.12.1 - Master workflow - Header attributes • 11
SR3146.9.12.2 - Master workflow - Review mode • 12
SR3146.9.2.4.6 - Library for signature privileges • 124
SR3146.9.3.1 - Master recipe - Header attributes • 9
SR3146.9.3.2 - Master recipe - Review mode • 11
SR3146.9.4.3+ - Validator consistency checks • 48
SR3146.9.5.3 - Process inputs - Signature privilege • 17
SR3146.9.5.3.1 - Process inputs - Signature privilege attributes • 17

Conventions (typographical) • 1

D

Designer workbench (SR3146.9+) • 117

E

Export/Import of Master Recipes, Master Workflows, and Building Blocks (SR1075.1+) • 145
Checksum (SR1075.1.1.2) • 145
General prerequisites and usage of B2MML (SR1075.1.1) • 145
Handling of status and version (SR1075.1.1.3) • 146
Handling of status transition history (SR1075.1.1.4) • 146
Prerequisites for import (SR1075.1.1.1) • 145

F

FlexibleAttributeDefinition editor • 140
FlexibleStateModel editor • 140
FlexibleTagDefinition editor • 140

I

Information flow (SR3146.9.9+) • 42
No Distributed Execution (SR3146.9.9.2.1) • 43
Output mechanism (SR3146.9.9.1) • 43
Referencing (SR3146.9.9.2) • 43
Validator (SR3146.9.9.3) • 43

L

Library concept (SR3146.9.11+) • 143
Basic building block import (SR3146.9.11.1) • 143
Custom building block (SR3146.9.11.3) • 143

M

Management of building blocks (SR3146.10+) • 73
Content (SR3146.10.3.1) • 77
Filter (SR3146.10.3.3) • 77
Header attributes (SR3146.10.4) • 73
Objects information (SR3146.10.3.2) • 77
Revisioning (SR3146.10.2) • 76
Status handling of building blocks (SR3146.10.1) • 75
Usage list (SR3146.10.3) • 76
Workbench object operations - New building block dialog (SR3146.9.10.7.2) • 75
Workbench object operations - View locked building block (SR3146.9.10.7.1) • 75
Workbench object operations (SR3146.9.10.7) • 74
Management of change requests (SR3146.12+) • 83
Action list (SR3146.12.4) • 90
Cancel Change Request (SR3146.12.5.1) • 93
Define old/new building blocks (SR3146.12.3.1) • 89
Execute change request (SR3146.12.5) • 91
Execution status (SR3146.12.5.2) • 93
Header attributes (SR3146.12.1) • 83
Master recipe-specific change requests (SR3146.12.2.2) • 85
Master workflow-specific change requests (SR3146.12.2.3) • 87
Open dialog for projected objects (SR3146.12.6.2) • 84
Scope definition (SR3146.12.3+) • 89
Scope validity check (SR3146.12.2.1) • 89
Select from usage list (SR3146.12.3.2) • 90
Status management of change requests (SR3146.12.2+) • 85
Workbench object operations - View locked change request (SR3146.12.6.1) • 85
Workbench object operations (SR3146.12.6) • 84
Management of master recipes and master workflows (SR3146.4+) • 7

- Master recipe model - S88 (SR3146.9.3) • 9
 - Master recipe - Header attributes (SR3146.9.3.1) • 9
 - Master recipe - Review mode (SR3146.9.3.2) • 11
- Master recipe report - Batch (SR3146.11+) • 95
 - MR report - Approval record (SR3146.11.4.3) • 97
 - MR report - Basic data (SR3146.11.4.1) • 96
 - MR report - Bill of equipment (SR3146.11.4.16) • 97
 - MR report - Bill of material (SR3146.11.4.4) • 97
 - MR report - Comments (SR3146.11.4.2) • 97
 - MR report - Comparison with baseline (SR3146.11.4.17) • 100
 - MR report - Critical process parameters (SR3146.11.4.7) • 98
 - MR report - Critical quality attributes (SR3146.11.4.8) • 98
 - MR report - Key process parameters (SR3146.11.4.9) • 99
 - MR report - Key quality attributes (SR3146.11.4.10) • 99
 - MR report - Material flow control (SR3146.11.4.6) • 98
 - MR report - Operation (SR3146.11.4.13) • 100
 - MR report - Phase (SR3146.11.4.14) • 100
 - MR report - Procedure (SR3146.11.4.11) • 99
 - MR report - SFC pagination (SR3146.11.1.1) • 95
 - MR report - Table of contents (SR3146.11.4.15) • 96
 - MR report - Unit procedure (SR3146.11.4.12) • 99
 - MR report - Work center assignments (SR3146.11.4.5) • 98
 - MR report cover page (SR3146.11.5) • 96
 - MR report footer (SR3146.11.3) • 95
 - MR report header (SR3146.11.2) • 95
 - MR report sections (SR3146.11.4+) • 96
 - Print MR report (SR3146.11.1) • 95
- Master recipe report - Device (SR3146.14+) • 103
 - MR report - Approval record (SR3146.14.4.3) • 105
 - MR report - Basic data (SR3146.14.4.1) • 104
 - MR report - Bill of equipment (SR3146.14.4.16) • 105
 - MR report - Bill of materials (SR3146.14.4.4) • 105
 - MR report - Comments (SR3146.14.4.2) • 105
 - MR report - Comparison with baseline (SR3146.14.4.17) • 108
 - MR report - Critical process parameters (SR3146.14.4.7) • 106
 - MR report - Critical quality attributes (SR3146.14.4.8) • 106
 - MR report - Key process parameters (SR3146.14.4.9) • 107
 - MR report - Key quality attributes (SR3146.14.4.10) • 107
 - MR report - Material flow control (SR3146.14.4.6) • 106
 - MR report - Operation (SR3146.14.4.13) • 108
 - MR report - Phase (SR3146.14.4.14) • 108
 - MR report - Procedure (SR3146.14.4.11) • 107
 - MR report - SFC pagination (SR3146.14.1.1) • 103
 - MR report - Table of contents (SR3146.14.4.15) • 104
 - MR report - Unit procedure (SR3146.14.4.12) • 107
 - MR report - Work center assignments (SR3146.14.4.5) • 106
 - MR report cover page (SR3146.14.5) • 104
 - MR report footer (SR3146.14.3) • 103
 - MR report header (SR3146.14.2) • 103
 - MR report sections (SR3146.14.4+) • 104
 - Print MR report (SR3146.14.1) • 103
- Master workflow model - S88 (SR3146.12.3) • 11
 - Master workflow - Header attributes (SR3146.9.12.1) • 11
 - Master workflow - Review mode (SR3146.9.12.2) • 12
- Master workflow report (SR3146.13+) • 111
 - MWF report - Approval record (SR3146.13.4.3) • 113
 - MWF report - Basic data (SR3146.13.4.1) • 112
 - MWF report - Bill of equipment (SR3146.13.4.14) • 113
 - MWF report - Comments (SR3146.13.4.2) • 113
 - MWF report - Comparison with baseline (SR3146.13.4.17) • 116
 - MWF report - Critical process parameters (SR3146.13.4.5) • 114
 - MWF report - Critical quality attributes (SR3146.13.4.6) • 114
 - MWF report - Key process parameters (SR3146.13.4.7) • 114

MWF report - Key quality attributes (SR3146.13.4.8) • 114

MWF report - List of materials (SR3146.13.4.15) • 113

MWF report - Operation (SR3146.13.4.11) • 115

MWF report - Phase (SR3146.13.4.12) • 115

MWF report - Procedure (SR3146.13.4.9) • 114

MWF report - SFC pagination (SR3146.13.1.1) • 111

MWF report - Table of contents (SR3146.13.4.13) • 112

MWF report - Unit procedure (SR3146.13.4.10) • 115

MWF report - Work center assignments (SR3146.13.4.4) • 113

MWF report cover page (SR3146.13.5) • 112

MWF report footer (SR3146.13.3) • 111

MWF report header (SR3146.13.2) • 111

MWF report sections (SR3146.13.4+) • 112

Print MWF report (SR3146.13.1) • 111

Material Flow Control (MFC) (SR3146.9.8) • 37

Automated update of MFC representation (SR3146.9.8.10) • 42

Automatic merge - Default MFC definition (SR3146.9.8.3) • 39

Default MFC definition between unit procedures (SR3146.9.8.5) • 40

Default MFC definition within a unit procedure (SR3146.9.8.4) • 40

Graphical view (SR3146.9.8.1) • 37

Manual merge between unit procedures (SR3146.9.8.7) • 41

Manual merge of MFC inputs with the final output. (SR3146.9.8.13) • 41

Manual merge within a unit procedure (SR3146.9.8.6) • 40

MFC definition linked to procedure level (SR3146.9.8.12) • 42

Reset (SR3146.9.8.9) • 42

Split (SR3146.9.8.8) • 42

Table view (SR3146.9.8.2) • 38

Tooltips (SR3146.9.8.2) • 39

Measured value Editor • 139

Messages window (user interface) • 128

Multi-line text editor • 139

P

Procedural structure - Building blocks • 13

Procedural structure - Links and branches • 14

Procedural structure - Transitions • 14

Procedural structure (SR3146.9.4) • 13

Bundle output variables (SR3146.9.7.4.2) • 26

Bundle process parameters (SR3146.9.7.4.1) • 26

Dynamic creation of process parameters (SR3146.9.7.4) • 26

Frozen parameters (SR3146.9.13.2) • 27

Frozen transitions (SR3146.9.14.2) • 27

Locked parameters (SR3146.9.13.1) • 27

Locked transitions (SR3146.9.14.1) • 27

Locking of parameters (SR3146.9.13+) • 27

Locking of transitions (SR3146.9.14+) • 27

Procedural structure - Building blocks • 13

Procedural structure - Conditions (SR3146.9.4.2) • 15

Procedural structure - Hidden phase (SR3146.9.4.5) • 15

Procedural structure - Links and branches • 14

Procedural structure - Operation type (SR3146.9.4.4) • 15

Procedural structure - Parallel execution (SR3146.9.4.1) • 14

Procedural structure - Transitions • 14

Process inputs - Abort-and-reactivate-enabled capability (SR3146.9.5.10.3) • 23

Process inputs - Auto-startable capability (SR3146.9.5.9.6) • 22

Process inputs - Capabilities on operation level (SR3146.9.5.9) • 20

Process inputs - Capabilities on unit procedure level (SR3146.9.5.10) • 23

Process inputs - Detachable capability - Operation (SR3146.9.5.9.1) • 20

Process inputs - Detachable capability - Unit procedure (SR3146.9.5.10.2) • 23

Process inputs - Equipment requirement (class) (SR3146.9.5.5) • 19

Process inputs - Equipment requirement (property type) (SR3146.9.5.6) • 19

Process inputs - Equipment requirement (rule)
 (SR3146.9.5.7) • 20

Process inputs - Escalation-enabled capability
 (SR3146.9.5.9.3) • 21

Process inputs - Event-triggered capability
 (SR3146.9.5.9.2) • 21

Process inputs - Material (SR3146.9.5.1) • 16

Process inputs - Material attributes (SR3146.9.5.1.1) • 16

Process inputs - Pause-enabled capability
 (SR3146.9.5.10.1) • 23

Process inputs - Requirement attributes (class)
 (SR3146.9.5.5.1) • 19

Process inputs - Requirement attributes (property type)
 (SR3146.9.5.6.1) • 19

Process inputs - Requirement attributes (rule)
 (SR3146.9.5.7.1) • 20

Process inputs - Server-run capability (SR3146.9.5.9.5) • 22

Process inputs - Signature privilege (SR3146.9.5.3) • 17

Process inputs - Signature privilege attributes
 (SR3146.9.5.3.1) • 17

Process inputs - Station (SR3146.9.5.8) • 18

Process inputs - Station attributes (SR3146.9.5.8.1) • 19

Process inputs - Transition (SR3146.9.5.11) • 20

Process inputs - Transition attributes (SR3146.9.5.11.1) • 20

Process inputs - Trigger-enabled capability
 (SR3146.9.5.9.4) • 22

Process inputs - Work center (SR3146.9.5.4) • 18

Process inputs - Work center attributes
 (SR3146.9.5.4.1) • 18

Process inputs (SR3146.9.5) • 15

Process outputs - Material (SR3146.9.6.1) • 24

Process outputs - Material attributes (SR3146.9.6.1.1) • 24

Process outputs (SR3146.9.6+) • 24

Process parameters - Classification (SR3146.9.7.1) • 26

Process parameters (SR3146.9.7+) • 25

Property windows • 127

R

Recipe and workflow management (SR3146+) • 3

Access control (SR3071.6.2) • 4

Access control (SR3071.7.2) • 5

Access control (SR3071.9.2) • 5

Recipe Designer - Batch (SR3071.6+) • 4

Recipe Designer - Device (SR3071.9+) • 5

Workflow Designer (SR3071.7+) • 5

S

SR1075.1.1 - General prerequisites and usage of B2MML • 145

SR1075.1.1.1 - Prerequisites for import • 145

SR1075.1.1.2 - Checksum • 145

SR1075.1.1.3 - Handling of status and version • 146

SR1075.1.1.4 - Handling of status transition history • 146

SR1075.1+ - Export/Import of Master Recipes, Master Workflows, and Building Blocks • 145

SR3071.6.2 - Access control • 4

SR3071.6+ - Recipe Designer - Batch • 4

SR3071.7.2 - Access control • 5

SR3071.7+ - Workflow Designer • 5

SR3071.9.2 - Access control • 5

SR3071.9+ - Recipe Designer - Device • 5

SR3146.10.1 - Status handling of building blocks • 75

SR3146.10.2 - Revisioning • 76

SR3146.10.3 - Usage list • 76

SR3146.10.3.1 - Content • 77

SR3146.10.3.2 - Objects information • 77

SR3146.10.3.3 - Filter • 77

SR3146.10.4 - Header attributes • 73

SR3146.10+ - Management of building blocks • 73

SR3146.11.1 - Print MR report • 95

SR3146.11.1.1 - MR report - SFC pagination • 95

SR3146.11.2 - MR report header • 95

SR3146.11.3 - MR report footer • 95

SR3146.11.4.1 - MR report - Basic data • 96

SR3146.11.4.10 - MR report - Key quality attributes • 99

SR3146.11.4.11 - MR report - Procedure • 99

SR3146.11.4.12 - MR report - Unit procedure • 99

SR3146.11.4.13 - MR report - Operation • 100	SR3146.13.2 - MWF report header • 111
SR3146.11.4.14 - MR report - Phase • 100	SR3146.13.3 - MWF report footer • 111
SR3146.11.4.15 - MR report - Table of contents • 96	SR3146.13.4.1 - MWF report - Basic data • 112
SR3146.11.4.16 - MR report - Bill of equipment • 97	SR3146.13.4.10 - MWF report - Unit procedure • 115
SR3146.11.4.17 - MR report - Comparison with baseline • 100	SR3146.13.4.11 - MWF report - Operation • 115
SR3146.11.4.2 - MR report - Comments • 97	SR3146.13.4.12 - MWF report - Phase • 115
SR3146.11.4.3 - MR report - Approval record • 97	SR3146.13.4.13 - MWF report - Table of contents • 112
SR3146.11.4.4 - MR report - Bill of materials • 97	SR3146.13.4.14 - MWF report - Bill of equipment • 113
SR3146.11.4.5 - MR report - Work center assignments • 98	SR3146.13.4.15 - MWF report - List of materials • 113
SR3146.11.4.6 - MR report - Material flow control • 98	SR3146.13.4.17 - MWF report - Comparison with baseline • 116
SR3146.11.4.7 - MR report - Critical process parameters • 98	SR3146.13.4.2 - MWF report - Comments • 113
SR3146.11.4.8 - MR report - Critical quality attributes • 98	SR3146.13.4.3 - MWF report - Approval record • 113
SR3146.11.4.9 - MR report - Key process parameters • 99	SR3146.13.4.4 - MWF report - Work center assignments • 113
SR3146.11.4+ - MR report sections • 96	SR3146.13.4.5 - MWF report - Critical process parameters • 114
SR3146.11.5 - MR report cover page • 96	SR3146.13.4.6 - MWF report - Critical quality attributes • 114
SR3146.11+ - Master recipe report - Batch • 95	SR3146.13.4.7 - MWF report - Key process parameters • 114
SR3146.12.1 - Header attributes • 83	SR3146.13.4.8 - MWF report - Key quality attributes • 114
SR3146.12.2.1 - Scope validity check • 89	SR3146.13.4.9 - MWF report - Procedure • 114
SR3146.12.2.2 - Master recipe-specific change requests • 85	SR3146.13.4+ - MWF report sections • 112
SR3146.12.2.3 - Master workflow-specific change requests • 87	SR3146.13.5 - MWF report cover page • 112
SR3146.12.2+ - Status management of change requests • 85	SR3146.13+ - Master workflow report • 111
SR3146.12.3.1 - Define old/new building blocks • 89	SR3146.14.1 - Print MR report • 103
SR3146.12.3.2 - Select from usage list • 90	SR3146.14.1.1 - MR report - SFC pagination • 103
SR3146.12.3+ - Scope definition • 89	SR3146.14.2 - MR report header • 103
SR3146.12.4 - Action list • 90	SR3146.14.3 - MR report footer • 103
SR3146.12.5 - Execute change request • 91	SR3146.14.4.1 - MR report - Basic data • 104
SR3146.12.5.1 - Cancel change request • 93	SR3146.14.4.10 - MR report - Key quality attributes • 107
SR3146.12.5.2 - Execution status • 93	SR3146.14.4.11 - MR report - Procedure • 107
SR3146.12.6 - Workbench object operations (change request) • 84	SR3146.14.4.12 - MR report - Unit procedure • 107
SR3146.12.6.1 - Workbench object operations - View locked change request • 85	SR3146.14.4.13 - MR report - Operation • 108
SR3146.12.6.2 - Open dialog for protected objects • 84	SR3146.14.4.14 - MR report - Phase • 108
SR3146.12+ - Management of change requests • 83	SR3146.14.4.15 - MR report - Table of contents • 104
SR3146.13.1 - Print MWF report • 111	SR3146.14.4.16 - MR report - Bill of equipment • 105
SR3146.13.1.1 - MWF report - SFC pagination • 111	SR3146.14.4.17 - MR report - Comparison with baseline • 108
	SR3146.14.4.2 - MR report - Comments • 105

- SR3146.14.4.3 - MR report - Approval record • 105
- SR3146.14.4.4 - MR report - Bill of materials • 105
- SR3146.14.4.5 - MR report - Work center assignments • 106
- SR3146.14.4.6 - MR report - Material flow control • 106
- SR3146.14.4.7 - MR report - Critical process parameters • 106
- SR3146.14.4.8 - MR report - Critical quality attributes • 106
- SR3146.14.4.9 - MR report - Key process parameters • 107
- SR3146.14.4+ - MR report sections • 104
- SR3146.14.5 - MR report cover page • 104
- SR3146.14+ - Master recipe report - Device • 103
- SR3146.15.1.1.1 - Master recipe-specific header attributes • 79
- SR3146.15.1.1.2 - Master workflow-specific header attributes • 79
- SR3146.15.1.1.3 - Custom building block-specific header attributes • 80
- SR3146.15.1.1+ - Comparing header attributes • 79
- SR3146.15.1.2 - Comparing element attributes • 80
- SR3146.15.1.3 - Comparing source attributes • 80
- SR3146.15.1.4.1 - Material parameter-specific attributes • 80
- SR3146.15.1.4.2 - Material Flow Control-specific attributes • 80
- SR3146.15.1.4.3 - Equipment requirement parameter-specific attributes • 81
- SR3146.15.1.4.4 - Work center assignment parameter-specific attributes • 81
- SR3146.15.1.4.5 - Privilege parameter-specific attributes • 81
- SR3146.15.1.4.6 - Capability parameter-specific attributes • 81
- SR3146.15.1.4.7 - Process parameter-specific attributes • 81
- SR3146.15.1.4+ - Comparing parameter attributes • 80
- SR3146.15.1.5.1 - Relation-specific characteristics • 81
- SR3146.15.1.5.2 - Transition-specific attributes • 81
- SR3146.15.1.5.3 - Location-specific characteristics • 82
- SR3146.15.1.5+ - Comparing procedural elements • 81
- SR3146.15.1+ - Comparison algorithm • 79
- SR3146.15.2 - Defining a comparison baseline • 82
- SR3146.15.3 - Performing a comparison • 82
- SR3146.15+ - Comparing master recipes, master workflows, and custom building blocks • 79
- SR3146.4.5 - Version control • 141
- SR3146.4.5.1 - Status change • 142
- SR3146.4.5.3 - Version history • 142
- SR3146.4+ - Management of master recipes and master workflows • 7
- SR3146.6 - Configurable version control • 142
- SR3146.6.1 - Flexible version graph assignment • 143
- SR3146.6.4 - Standard statuses for version control • 43
- SR3146.6.4.1 - Default version graph for master recipes • 44
- SR3146.6.4.2 - Period of validity • 142
- SR3146.6.4.3 - Consistency checks • 142
- SR3146.6.4.4 - Default version graph for master workflows • 46
- SR3146.6.7 - Flexible definition of approval details • 143
- SR3146.6.7.1 - Flexibility of user groups • 143
- SR3146.6.7.3 - Serial approval • 143
- SR3146.9.1.1 - Process view • 119
- SR3146.9.1.1.1 - Workbench - Start and end steps • 119
- SR3146.9.1.2 - Navigation • 119
- SR3146.9.1.3 - Context menus • 120
- SR3146.9.1.4 - Docking • 120
- SR3146.9.1.5 - Tooltips • 121
- SR3146.9.1+ - Workbench usability • 117
- SR3146.9.10.1 - Workbench object operations - Master recipe • 28
- SR3146.9.10.1.2 - Workbench object operations - New master recipe dialog • 29
- SR3146.9.10.1.3 - Workbench object operations - Prepare (deep save) • 29
- SR3146.9.10.1.4 - Workbench object operations - View locked master recipe • 29
- SR3146.9.10.10 - Workbench object operations - Smart replace • 32
- SR3146.9.10.11 - Workbench object operations - Replace material (product) • 36

- SR3146.9.10.6 - Workbench object operations - Element • 31
- SR3146.9.10.7 - Workbench object operations (custom building block) • 74
- SR3146.9.10.7.1 - Workbench object operations - View locked building block • 75
- SR3146.9.10.7.2 - Workbench object operations - New building block dialog • 75
- SR3146.9.10.9 - Workbench object operations - Master workflow • 29
- SR3146.9.10.9.1 - Workbench object operations - New master workflow dialog • 30
- SR3146.9.10.9.2 - Workbench object operations - Prepare (deep save) • 30
- SR3146.9.10.9.3 - Workbench object operations - View locked master workflow • 31
- SR3146.9.10+ - Workbench object operations • 28
- SR3146.9.11.1 - Basic building block import • 143
- SR3146.9.11.3 - Custom building block • 143
- SR3146.9.11+ - Library concept • 143
- SR3146.9.12 - Master workflow model - S88 • 11
- SR3146.9.12.1 - Master workflow - Header attributes • 11
- SR3146.9.12.2 - Master workflow - Review mode • 12
- SR3146.9.13.1 - Locked parameters • 27
- SR3146.9.13.2 - Frozen parameters • 27
- SR3146.9.13+ - Locking of parameters • 27
- SR3146.9.14.1 - Locked transitions • 27
- SR3146.9.14.2 - Frozen transitions • 27
- SR3146.9.14+ - Locking of transitions • 27
- SR3146.9.2.1 - Work area • 121
- SR3146.9.2.11 - Preview panel • 131
- SR3146.9.2.12 - Property window - Element • 128
- SR3146.9.2.13 - Property window - Source BB • 128
- SR3146.9.2.15 - Property window - Header • 128
- SR3146.9.2.16 - Messages • 128
- SR3146.9.2.16.1 - Messages - Validator-related • 129
- SR3146.9.2.16.2 - Messages - Category filter • 129
- SR3146.9.2.16.3 - Messages - Navigation • 129
- SR3146.9.2.16.4 - Messages - Resolve issue • 129
- SR3146.9.2.17.1 - Date/Time Picker editor • 139
- SR3146.9.2.17.2 - Duration editor • 139
- SR3146.9.2.17.3 - Packaging Level Data editor • 140
- SR3146.9.2.17+ - Other editors • 138
- SR3146.9.2.18 - Statistics • 140
- SR3146.9.2.18.1 - Statistics view • 141
- SR3146.9.2.18.2 - Statistics-based warning • 141
- SR3146.9.2.19.1 - Comparison - Difference categories • 130
- SR3146.9.2.19.2 - Comparison - Category filter • 131
- SR3146.9.2.19.3 - Comparison - Type filter • 131
- SR3146.9.2.19.4 - Comparison - Navigation • 131
- SR3146.9.2.19.5 - Comparison - Difference-specific information • 131
- SR3146.9.2.19+ - Comparison • 130
- SR3146.9.2.20.1 - Open dialog for protected objects • 122
- SR3146.9.2.20+ - Open dialog • 122
- SR3146.9.2.4.1 - Library for work centers • 124
- SR3146.9.2.4.2 - Library for materials (ERP BOM) • 123
- SR3146.9.2.4.3 - Library filter & search • 122
- SR3146.9.2.4.4 - Library for building blocks • 123
- SR3146.9.2.4.5 - Library for materials • 123
- SR3146.9.2.4.6 - Library for signature privileges • 124
- SR3146.9.2.4.7 - Library for equipment classes • 124
- SR3146.9.2.4.8 - Library for property types • 124
- SR3146.9.2.4.9 - Library for stations • 124
- SR3146.9.2.4+ - Library Universe • 122
- SR3146.9.2.7 - Parameter Panel • 126
- SR3146.9.2.7.1 - Parameter Panel - View all • 127
- SR3146.9.2.7.2 - Parameter Panel - Specific property type • 127
- SR3146.9.2.8 - Parameter Panel filter & search • 127
- SR3146.9.2.9 - Setlist • 124
- SR3146.9.2.9.1 - Setlist - ERP BOM item • 125
- SR3146.9.2.9.2 - Setlist - Workflow-related consistency checks • 125
- SR3146.9.2.9.3 - Setlist - Workflow-related data adjustment • 126
- SR3146.9.2+ - Workbench UI elements • 121
- SR3146.9.3 - Master recipe model - S88 • 9
- SR3146.9.3.1 - Master recipe - Header attributes • 9
- SR3146.9.3.2 - Master recipe - Review mode • 11
- SR3146.9.4 - Procedural structure • 13

- SR3146.9.4.1 - Procedural structure - Parallel execution • 14
- SR3146.9.4.2 - Procedural structure - Conditions • 15
- SR3146.9.4.3.1 - Validator - Dummy elements • 48
- SR3146.9.4.3.10 - Validator - Minimum required material parameters • 56
- SR3146.9.4.3.10.1 - Validator - Minimum required material parameters - Warning • 56
- SR3146.9.4.3.10.2 - Validator - Minimum required material parameters - Error • 56
- SR3146.9.4.3.100 - Validator - Non-production MR - More than one unit procedure or operation • 65
- SR3146.9.4.3.101 - Trigger phase without event-triggered operation • 71
- SR3146.9.4.3.102 - Server-run validator - Conflicting station assignment • 71
- SR3146.9.4.3.103 - Server-run validator - Phase cannot be run on a server • 71
- SR3146.9.4.3.104 - Server-run validator - Phase needs to be run on a server • 71
- SR3146.9.4.3.105 - Incompatible capability assignment between different levels • 69
- SR3146.9.4.3.106 - Usage of capability not allowed for simultaneous branches • 69
- SR3146.9.4.3.107 - Validator - - Packaging level of master recipe • 53
- SR3146.9.4.3.108 - Validator - - Packaging level of material parameter • 59
- SR3146.9.4.3.109 - Validator - Material parameter - Packaging level does not fit material's master data • 60
- SR3146.9.4.3.109.1 - Validator - Material parameter - Packaging level does not fit material's master data - Resolve • 60
- SR3146.9.4.3.110 - Validator - MFC - Output defined for device-specific recipe or workflow • 52
- SR3146.9.4.3.111 - Validator - Recipe - Min/Max quantities without UoM • 53
- SR3146.9.4.3.112 - Validator - Workflow - Only appendable if production-relevant • 65
- SR3146.9.4.3.113 - Validator - MFC - Unallowed planned quantity mode • 58
- SR3146.9.4.3.113.1 - Validator - MFC - Unallowed planned quantity mode - Warning • 58
- SR3146.9.4.3.113.2 - Validator - MFC - Unallowed planned quantity mode - Error • 58
- SR3146.9.4.3.114 - Validator - MFC - Fixed quantity not allowed • 59
- SR3146.9.4.3.115 - Validator - MFC - None planned quantity mode for Dispense unit procedure • 61
- SR3146.9.4.3.116 - Validator - MFC - Planned quantities of input and output parameters do not match • 59
- SR3146.9.4.3.117 - Validator - MFC - No direct output defined • 52
- SR3146.9.4.3.117.1 - Validator - MFC - No direct output defined - Warning • 52
- SR3146.9.4.3.117.2 - Validator - MFC - No direct output defined - Error • 52
- SR3146.9.4.3.118 - Validator - MFC - Final output for Dispense unit procedures • 61
- SR3146.9.4.3.118.1 - Validator - MFC - Final output for Dispense unit procedures - Warning • 61
- SR3146.9.4.3.118.2 - Validator - MFC - Final output for Dispense unit procedures - Error • 61
- SR3146.9.4.3.119 - Validator - MFC - Planned quantity mode is As produced - Tolerances are missing • 59
- SR3146.9.4.3.12 - Validator - Outgoing link • 48
- SR3146.9.4.3.120 - Validator - MFC - Unallowed weighing method • 62
- SR3146.9.4.3.121 - Validator - MFC - Planned quantity mode is As defined - Tolerances are missing • 62
- SR3146.9.4.3.122 - Validator - MFC - Planned quantity mode is As defined - Planned quantity is missing • 63
- SR3146.9.4.3.123 - Validator - MFC - Weighing material type not defined • 63
- SR3146.9.4.3.124 - Validator - MFC - MFC position is missing • 63
- SR3146.9.4.3.125 - Validator - MFC - Ingoing MFC transfer item – Unallowed weighing material type (Active) • 63
- SR3146.9.4.3.127 - Validator - D_Material - Planned potency not defined • 64

SR3146.9.4.3.128 - Validator - MFC - Fixed quantity not allowed for Active or Compensator • 59	SR3146.9.4.3.25.1 - Validator - Phase parameter - Broken link - Warning • 50
SR3146.9.4.3.129 - Validator - MFC - Planned quantity mode is As defined - Tolerances are missing • 64	SR3146.9.4.3.25.2 - Validator - Phase parameter - Broken link - Error • 50
SR3146.9.4.3.13 - Validator - Start and end step • 48	SR3146.9.4.3.27 - Validator - Transition condition is missing • 49
SR3146.9.4.3.130 - Validator - MFC - Planned quantity mode is As defined - Planned quantity is missing • 64	SR3146.9.4.3.3 - Validator - Loop • 48
SR3146.9.4.3.131 - Validator - MFC - Planned quantity mode is As defined - Allowed weighing method is missing • 64	SR3146.9.4.3.31 - Validator - Transition condition - Syntax • 49
SR3146.9.4.3.132 - Validator - D_Material - Default weighing method not allowed • 61	SR3146.9.4.3.35 - Validator - Dispense operation • 60
SR3146.9.4.3.133 - Validator - MFC (final output) - Planned quantity does not fit master recipe • 56	SR3146.9.4.3.36 - Validator - Phase parameter - Invalid reference • 50
SR3146.9.4.3.134 - Validator - MFC – Ingoing MFC transfer item – Unallowed weighing material type (Compensator, Filler) • 63	SR3146.9.4.3.37 - Validator - Transition condition - Broken link • 49
SR3146.9.4.3.135 - Validator - MFC – Ingoing MFC transfer item – Unallowed planned quantity mode • 64	SR3146.9.4.3.37.1 - Validator - Transition condition - Broken link - Warning • 49
SR3146.9.4.3.136 - Validator - MFC - Attributes of material parameters do not match • 52	SR3146.9.4.3.37.2 - Validator - Transition condition - Broken link - Error • 49
SR3146.9.4.3.137 - Validator - BOM position not used • 58	SR3146.9.4.3.38 - Validator - Transition condition - Invalid reference • 49
SR3146.9.4.3.138 - Validator - MFC - Outgoing or incoming connector is missing • 52	SR3146.9.4.3.40 - Station assignment does not match a given work center assignment • 68
SR3146.9.4.3.139 - Validator - MFC - More than one outgoing or incoming connector on one execution path • 53	SR3146.9.4.3.42 - Station assignment is missing for a given work center assignment • 68
SR3146.9.4.3.14 - Validator - Simultaneous branch • 48	SR3146.9.4.3.43 - Unusable station • 68
SR3146.9.4.3.15 - Validator - Selection branch • 48	SR3146.9.4.3.44 - Unusable work center • 68
SR3146.9.4.3.17 - Validator - Simultaneous branch - Unsuitable components • 48	SR3146.9.4.3.45 - Usage of capability not allowed for Dispense • 68
SR3146.9.4.3.18 - Validator - Errors in preceding components • 49	SR3146.9.4.3.46 - Incompatible capability assignment • 69
SR3146.9.4.3.19 - Validator - Start step path • 49	SR3146.9.4.3.47 - Capability assignment is missing • 70
SR3146.9.4.3.20 - Validator - End step path • 49	SR3146.9.4.3.5 - Validator - Double signature for sequential phase completion • 66
SR3146.9.4.3.21 - Validator - Phase parameter - Syntax • 50	SR3146.9.4.3.50 - Validator - Phase parameter - Missing ERP BOM position • 55
SR3146.9.4.3.25 - Validator - Phase parameter - Broken link • 50	SR3146.9.4.3.51 - Validator - Phase parameter - Multiple ERP BOM positions • 55
	SR3146.9.4.3.52 - Validator - Master recipe and material parameter - Conflicting inventory level • 53
	SR3146.9.4.3.53 - Validator - Material parameter - Conflicting inventory level • 59

- SR3146.9.4.3.55 - Validator - Phase parameter - Non-ERP
BOM position • 55
- SR3146.9.4.3.56 - Validator - Phase parameter - ERP
BOM attribute discrepancy • 55
- SR3146.9.4.3.56.1 - Validator - Phase parameter - ERP
BOM attribute discrepancy - Resolve • 55
- SR3146.9.4.3.57 - Validator - Recipe - Planned quantity
does not fit ERP BOM • 54
- SR3146.9.4.3.57.1 - Validator - Recipe - Planned quantity
does not fit ERP BOM - Resolve • 54
- SR3146.9.4.3.58 - Validator - Recipe - Packaging level
does not fit ERP BOM • 54
- SR3146.9.4.3.58.1 - Validator - Recipe - Packaging level
does not fit ERP BOM - Resolve • 54
- SR3146.9.4.3.6 - Validator - (Sequential) phase completion
signatures • 66
- SR3146.9.4.3.61 - Validator - MFC - Final output does not
match material (product) • 50
- SR3146.9.4.3.62 - Validator - MFC - Multiple outputs for
non-Dispense unit procedures • 50
- SR3146.9.4.3.62.1 - Validator - MFC - Multiple outputs
for non-Dispense unit procedures - Warning • 50
- SR3146.9.4.3.62.2 - Validator - MFC - Multiple outputs
for non-Dispense unit procedures - Error • 50
- SR3146.9.4.3.63 - Validator - MFC - No output defined for
input • 51
- SR3146.9.4.3.63.1 - Validator - MFC - No output defined
for input - Warning • 51
- SR3146.9.4.3.63.2 - Validator - MFC - No output defined
for input - Error • 51
- SR3146.9.4.3.64 - Validator - MFC - No input defined
MFC output • 51
- SR3146.9.4.3.64.1 - Validator - MFC - No input defined
for output - Warning • 51
- SR3146.9.4.3.64.2 - Validator - MFC - No input defined
for output - Error • 51
- SR3146.9.4.3.65 - Validator - MFC - Transfer input for
Dispense unit procedures • 60
- SR3146.9.4.3.65.1 - Validator - MFC - Transfer input for
Dispense unit procedures - Warning • 60
- SR3146.9.4.3.65.2 - Validator - MFC - Transfer input for
Dispense unit procedures - Error • 60
- SR3146.9.4.3.66 - Validator - MFC - Input without
planned quantity • 56
- SR3146.9.4.3.66.1 - Validator - MFC - Input without
planned quantity - Warning • 56
- SR3146.9.4.3.66.2 - Validator - MFC - Input without
planned quantity - Error • 56
- SR3146.9.4.3.67 - Validator - Recipe - Reference quantity
not defined • 53
- SR3146.9.4.3.68 - Validator - MFC - Uniqueness and
format of MFC position • 57
- SR3146.9.4.3.69 - Validator - Material - Incompatible
UoM • 57
- SR3146.9.4.3.7 - Server-run validator - No phase
completion signature or mandatory comment allowed •
72
- SR3146.9.4.3.70 - Validator - Material - UoM of planned
potency • 57
- SR3146.9.4.3.71 - Validator - Material - One tolerance is
missing • 57
- SR3146.9.4.3.72 - Validator - D_Material - Default
weighing method not defined • 61
- SR3146.9.4.3.73 - Validator - D_Material - Allowed
weighing methods not defined • 61
- SR3146.9.4.3.74 - Validator - D_Material - Tolerances not
defined • 62
- SR3146.9.4.3.75 - Validator - D_Material - Weighing
material type not defined • 62
- SR3146.9.4.3.76 - Validator - D_Material - Planned
potency not defined • 62
- SR3146.9.4.3.77 - Validator - Material - Negative values •
58
- SR3146.9.4.3.78 - Validator - MFC - MFC position is
missing • 58
- SR3146.9.4.3.78.1 - Validator - MFC - MFC position is
missing - Warning • 58
- SR3146.9.4.3.78.2 - Validator - MFC - MFC position is
missing - Error • 58
- SR3146.9.4.3.79 - Usage of capability not allowed for
workflows • 68

SR3146.9.4.3.8 - Validator - Usage type • 65	SR3146.9.4.3.94 - Escalation validator - Duration is missing for enabled exceptions • 70
SR3146.9.4.3.80 - Validator - Workflow - Phase is incompatible • 65	SR3146.9.4.3.95 - Escalation validator - Durations are conflicting • 70
SR3146.9.4.3.81 - Validator - Workflow - Workflow type not defined • 65	SR3146.9.4.3.96 - Escalation validator - No duration defined • 70
SR3146.9.4.3.82 - Validator - Property definition is inconsistent • 66	SR3146.9.4.3.96.1 - Escalation validator - No duration defined - Warning • 70
SR3146.9.4.3.83 - Validator - Minimum required equipment parameters • 66	SR3146.9.4.3.96.2 - Escalation validator - No duration defined - Error • 70
SR3146.9.4.3.83.1 - Validator - Minimum required equipment parameters - Warning • 66	SR3146.9.4.3.97 - Trigger validator - No valid trigger reference defined • 70
SR3146.9.4.3.83.2 - Validator - Minimum required equipment parameters - Error • 66	SR3146.9.4.3.97.1 - Trigger validator - No valid trigger reference defined - Warning • 70
SR3146.9.4.3.85 - Validator - Equipment class requirements are not unique • 66	SR3146.9.4.3.97.2 - Trigger validator - No valid trigger reference defined - Error • 70
SR3146.9.4.3.86 - Validator - Equipment requirement identifier is not unique • 67	SR3146.9.4.3.98 - Trigger validator - Trigger identifiers are not unique • 71
SR3146.9.4.3.87 - Validator - Equipment requirement rule - Identifier is missing • 67	SR3146.9.4.3.99 - Trigger validator - Triggers are missing • 71
SR3146.9.4.3.88 - Validator - Equipment requirement - BOE position is not unique • 66	SR3146.9.4.3+ - Validator consistency checks • 48
SR3146.9.4.3.89 - Validator - Equipment requirement rule - Syntax error • 67	SR3146.9.4.4 - Procedural structure - Operation type • 15
SR3146.9.4.3.9 - Validator - Exception signature supported • 65	SR3146.9.4.5 - Procedural structure - Hidden phase • 15
SR3146.9.4.3.90 - Validator - Equipment requirement rule - Broken link • 67	SR3146.9.5 - Process inputs • 15
SR3146.9.4.3.90.1 - Validator - Equipment requirement rule - Broken link - Warning • 66	SR3146.9.5.1 - Process inputs - Material • 16
SR3146.9.4.3.90.2 - Validator - Equipment requirement rule - Broken link - Error • 67	SR3146.9.5.1.1 - Process inputs - Material attributes • 16
SR3146.9.4.3.91 - Validator - Equipment requirement rule - Invalid reference • 67	SR3146.9.5.10 - Process inputs - Capabilities on unit procedure level • 23
SR3146.9.4.3.92 - Validator - Equipment requirement rule - Rule expression is missing • 67	SR3146.9.5.10.1 - Process inputs - Pause-enabled capability • 23
SR3146.9.4.3.93 - Capability assignment does not match usage type of privilege parameter • 69	SR3146.9.5.10.2 - Process inputs - Detachable capability - Unit procedure • 23
SR3146.9.4.3.93.1 - Capability assignment does not match usage type of privilege parameter - Warning • 69	SR3146.9.5.10.3 - Process inputs - Abort-and-reactivate-enabled capability • 23
SR3146.9.4.3.93.2 - Capability assignment does not match usage type of privilege parameter - Error • 69	SR3146.9.5.11 - Process inputs - Transition • 20
	SR3146.9.5.11.1 - Process inputs - Transition attributes • 20
	SR3146.9.5.3 - Process inputs - Signature privilege • 17
	SR3146.9.5.3.1 - Process inputs - Signature privilege attributes • 17
	SR3146.9.5.4 - Process inputs - Work center • 18

- SR3146.9.5.4.1 - Process inputs - Work center attributes • 18
- SR3146.9.5.5 - Process inputs - Equipment requirement (class) • 19
- SR3146.9.5.5.1 - Process inputs - Requirement attributes (class) • 19
- SR3146.9.5.6 - Process inputs - Equipment requirement (property type) • 19
- SR3146.9.5.6.1 - Process inputs - Requirement attributes (property type) • 19
- SR3146.9.5.7 - Process inputs - Equipment requirement (rule) • 20
- SR3146.9.5.7.1 - Process inputs - Requirement attributes (rule) • 20
- SR3146.9.5.8 - Process inputs - Station • 18
- SR3146.9.5.8.1 - Process inputs - Station attributes • 19
- SR3146.9.5.9 - Process inputs - Capabilities on operation level • 20
- SR3146.9.5.9.1 - Process inputs - Detachable capability - Operation • 20
- SR3146.9.5.9.2 - Process inputs - Event-triggered capability • 21
- SR3146.9.5.9.3 - Process inputs - Escalation-enabled capability • 21
- SR3146.9.5.9.4 - Process inputs - Trigger-enabled capability • 22
- SR3146.9.5.9.5 - Process inputs - Server-run capability • 22
- SR3146.9.5.9.6 - Process inputs - Auto-startable capability • 22
- SR3146.9.6.1 - Process outputs - Material • 24
- SR3146.9.6.1.1 - Process outputs - Material attributes • 24
- SR3146.9.6+ - Process outputs • 24
- SR3146.9.7.1 - Process parameters - Classification • 26
- SR3146.9.7.4 - Dynamic creation of process parameters • 26
- SR3146.9.7.4.1 - Bundle process parameters • 26
- SR3146.9.7.4.2 - Bundle output variables • 26
- SR3146.9.7+ - Process parameters • 25
- SR3146.9.8 - Material Flow Control (MFC) • 37
- SR3146.9.8.1 - Graphical view • 37
- SR3146.9.8.1.1 - Table view • 38
- SR3146.9.8.10 - Automated update of MFC representation • 42
- SR3146.9.8.12 - MFC definition linked to procedure level • 42
- SR3146.9.8.13 - Manual merge of MFC inputs with the final output • 41
- SR3146.9.8.2 - Tooltips • 39
- SR3146.9.8.3 - Automatic merge - Default MFC definition • 39
- SR3146.9.8.4 - Default MFC definition within a unit procedure • 40
- SR3146.9.8.5 - Default MFC definition between unit procedures • 40
- SR3146.9.8.6 - Manual merge within a unit procedure • 40
- SR3146.9.8.7 - Manual merge between unit procedures • 41
- SR3146.9.8.8 - Split • 42
- SR3146.9.8.9 - Reset • 42
- SR3146.9.9.1 - Output mechanism • 43
- SR3146.9.9.2 - Referencing • 43
- SR3146.9.9.2.1 - No Distributed Execution • 43
- SR3146.9.9.3 - Validator • 43
- SR3146.9.9.4 - Expression editor • 131
- SR3146.9.9.4.1 - Default preview buildtime • 131
- SR3146.9.9.4.10 - Expressions for flexible rules • 133
- SR3146.9.9.4.11 - Expression editor - Phase-specific functions • 135
- SR3146.9.9.4.12 - Expressions for conditional rules • 133
- SR3146.9.9.4.2 - Expressions for transitions • 132
- SR3146.9.9.4.3 - Expressions for process parameters • 133
- SR3146.9.9.4.4 - Expression editor operations • 134
- SR3146.9.9.4.5 - Expression editor literals • 134
- SR3146.9.9.4.6 - Expression editor - Pnuts • 138
- SR3146.9.9.4.7 - Expression editor functions • 135
- SR3146.9.9.4.8 - Expression editor - Implicit null handling • 138
- SR3146.9.9.4.9 - Expression editor - Explicit null handling • 138
- SR3146.9.9.6 - Expression editor - Runtime context data • 135

SR3146.9.9+ - Information flow • 42
SR3146.9+ - Designer workbench • 117
SR3146+ - Recipe and workflow management • 3
Standard statuses for version control (SR3146.6.4) • 43
 Default version graph for master recipes (SR3146.6.4.1)
 • 44
 Default version graph for master workflows
 (SR3146.6.4.4) • 46
String editor • 140

V

Validator consistency checks (SR3146.9.4.3+) • 48
 Capability assignment does not match usage type of
 privilege parameter - Error (SR3146.9.4.3.93.2) • 69
 Capability assignment does not match usage type of
 privilege parameter - Warning (SR3146.9.4.3.93.1) •
 69
 Capability assignment does not match usage type of
 privilege parameter (SR3146.9.4.3.93) • 69
 Capability assignment is missing (SR3146.9.4.3.47) •
 70
 Checks that are disabled for master recipes of the Cost
 center usage type • 65
 Escalation validator - Duration is missing for enabled
 exceptions (SR3146.9.4.3.94) • 70
 Escalation validator - Durations are conflicting
 (SR3146.9.4.3.95) • 70
 Escalation validator - No duration defined - Error
 (SR3146.9.4.3.96.2) • 70
 Escalation validator - No duration defined - Warning
 (SR3146.9.4.3.96.1) • 70
 Escalation validator - No duration defined
 (SR3146.9.4.3.96) • 70
 Incompatible capability assignment (SR3146.9.4.3.46) •
 69
 Incompatible capability assignment between different
 levels (SR3146.9.4.3.105) • 69
 Server-run validator - Conflicting station assignment
 (SR3146.9.4.3.102) • 71
 Server-run validator - No phase completion signature or
 mandatory comment allowed (SR3146.9.4.3.7) • 72

Server-run validator - Phase cannot be run on a server
 (SR3146.9.4.3.103) • 71
Station assignment does not match a given work center
 assignment (SR3146.9.4.3.40) • 68
Station assignment is missing for a given work center
 assignment (SR3146.9.4.3.42) • 68
Trigger phase without event-triggered operation
 (SR3146.9.4.3.101) • 71
Trigger validator - No valid trigger reference defined -
 Error (SR3146.9.4.3.97.2) • 70
Trigger validator - No valid trigger reference defined -
 Warning (SR3146.9.4.3.97.1) • 70
Trigger validator - No valid trigger reference defined
 (SR3146.9.4.3.97) • 70
Trigger validator - Trigger identifiers are not unique
 (SR3146.9.4.3.98) • 71
Trigger validator - Triggers are missing
 (SR3146.9.4.3.99) • 71
Unusable station (SR3146.9.4.3.43) • 68
Unusable work center (SR3146.9.4.3.44) • 68
Usage of capability not allowed for Dispense
 (SR3146.9.4.3.45) • 68
Usage of capability not allowed for simultaneous
 branches (SR3146.9.4.3.106) • 69
Usage of capability not allowed for workflows
 (SR3146.9.4.3.79) • 68
Validator - (Sequential) phase completion signatures
 (SR3146.9.4.3.6) • 66
Validator - BOM position not used (SR3146.9.4.3.137)
 • 58
Validator - D_Material - Allowed weighing methods
 not defined (SR3146.9.4.3.73) • 61
Validator - D_Material - Default weighing method not
 allowed (SR3146.9.4.3.132) • 61
Validator - D_Material - Default weighing method not
 defined (SR3146.9.4.3.72) • 61
Validator - D_Material - Planned potency not defined
 (SR3146.9.4.3.127) • 64
Validator - D_Material - Planned potency not defined
 (SR3146.9.4.3.76) • 62

- Validator - D_Material - Tolerances not defined
(SR3146.9.4.3.74) • 62
- Validator - D_Material - Weighing material type not
defined (SR3146.9.4.3.75) • 62
- Validator - Dispense operation (SR3146.9.4.3.35) • 60
- Validator - Double signature for sequential phase
completion (SR3146.9.4.3.5) • 66
- Validator - Dummy elements (SR3146.9.4.3.1) • 48
- Validator - End step path (SR3146.9.4.3.20) • 49
- Validator - Equipment class requirements are not
unique (SR3146.9.4.3.85) • 66
- Validator - Equipment requirement - BOE position is
not unique (SR3146.9.4.3.88) • 66
- Validator - Equipment requirement identifier is not
unique (SR3146.9.4.3.86) • 67
- Validator - Equipment requirement rule - Broken link -
Error (SR3146.9.4.3.90.2) • 67
- Validator - Equipment requirement rule - Broken link -
Warning (SR3146.9.4.3.90.1) • 67
- Validator - Equipment requirement rule - Broken link
(SR3146.9.4.3.90) • 67
- Validator - Equipment requirement rule - Identifier is
missing (SR3146.9.4.3.87) • 67
- Validator - Equipment requirement rule - Invalid
reference (SR3146.9.4.3.91) • 67
- Validator - Equipment requirement rule - Rule
expression is missing (SR3146.9.4.3.92) • 67
- Validator - Equipment requirement rule - Syntax error
(SR3146.9.4.3.89) • 67
- Validator - Errors in preceding components
(SR3146.9.4.3.18) • 49
- Validator - Exception signature supported
(SR3146.9.4.3.9) • 65
- Validator - Loop (SR3146.9.4.3.3) • 48
- Validator - Master recipe and material parameter -
Conflicting inventory level (SR3146.9.4.3.52) • 53
- Validator - Material - Incompatible UoM
(SR3146.9.4.3.69) • 57
- Validator - Material - Negative values
(SR3146.9.4.3.77) • 58
- Validator - Material - One tolerance is missing
(SR3146.9.4.3.71) • 57
- Validator - Material - UoM of planned potency
(SR3146.9.4.3.70) • 57
- Validator - Material parameter - Conflicting inventory
level (SR3146.9.4.3.53) • 59
- Validator - Material parameter - Packaging level does
not fit material's master data - Resolve
(SR3146.9.4.3.109.1) • 60
- Validator - Material parameter - Packaging level does
not fit material's master data (SR3146.9.4.3.109) •
60
- Validator - MFC - Attributes of material parameters do
not match (SR3146.9.4.3.136) • 52
- Validator - MFC - Final output does not match material
(product) (SR3146.9.4.3.61) • 50
- Validator - MFC - Final output for Dispense unit
procedures - Error (SR3146.9.4.3.118.2) • 61
- Validator - MFC - Final output for Dispense unit
procedures - Warning (SR3146.9.4.3.118.1) • 61
- Validator - MFC - Final output for Dispense unit
procedures (SR3146.9.4.3.118) • 61
- Validator - MFC - Fixed quantity not allowed
(SR3146.9.4.3.114) • 59
- Validator - MFC - Fixed quantity not allowed for
Active or Compensator (SR3146.9.4.3.128) • 59
- Validator - MFC – Ingoing MFC transfer item –
Unallowed planned quantity mode
(SR3146.9.4.3.135) • 64
- Validator - MFC – Ingoing MFC transfer item –
Unallowed weighing material type (Active)
(SR3146.9.4.3.125) • 63
- Validator - MFC – Ingoing MFC transfer item –
Unallowed weighing material type (Compensator,
Filler) (SR3146.9.4.3.134) • 63
- Validator - MFC - Input without planned quantity -
Error (SR3146.9.4.3.66.2) • 56
- Validator - MFC - Input without planned quantity -
Warning (SR3146.9.4.3.66.1) • 56
- Validator - MFC - Input without planned quantity
(SR3146.9.4.3.66) • 56

Validator - MFC - MFC position is missing - Error (SR3146.9.4.3.78.2) • 58	Validator - MFC - Planned quantities of input and output parameters do not match (SR3146.9.4.3.116) • 59
Validator - MFC - MFC position is missing - Warning (SR3146.9.4.3.78.1) • 58	Validator - MFC - Planned quantity mode is As defined - Allowed weighing method is missing (SR3146.9.4.3.131) • 64
Validator - MFC - MFC position is missing (SR3146.9.4.3.124) • 63	Validator - MFC - Planned quantity mode is As defined - Planned quantity is missing (SR3146.9.4.3.122) • 63
Validator - MFC - MFC position is missing (SR3146.9.4.3.78) • 58	Validator - MFC - Planned quantity mode is As defined - Planned quantity is missing (SR3146.9.4.3.130) • 64
Validator - MFC - More than one outgoing or incoming connector on one execution path (SR3146.9.4.3.139) • 53	Validator - MFC - Planned quantity mode is As defined - Tolerances are missing (SR3146.9.4.3.121) • 62
Validator - MFC - Multiple outputs for non-Dispense unit procedures - Error (SR3146.9.4.3.62.2) • 50	Validator - MFC - Planned quantity mode is As defined - Tolerances are missing (SR3146.9.4.3.129) • 64
Validator - MFC - Multiple outputs for non-Dispense unit procedures - Warning (SR3146.9.4.3.62.1) • 50	Validator - MFC - Planned quantity mode is As produced - Tolerances are missing (SR3146.9.4.3.119) • 59
Validator - MFC - Multiple outputs for non-Dispense unit procedures (SR3146.9.4.3.62) • 50	Validator - MFC - Transfer input for Dispense unit procedures - Error (SR3146.9.4.3.65.2) • 60
Validator - MFC - No direct output defined - Error (SR3146.9.4.3.117.2) • 52	Validator - MFC - Transfer input for Dispense unit procedures - Warning (SR3146.9.4.3.65.1) • 60
Validator - MFC - No direct output defined - Warning (SR3146.9.4.3.117.1) • 52	Validator - MFC - Transfer input for Dispense unit procedures (SR3146.9.4.3.65) • 60
Validator - MFC - No direct output defined (SR3146.9.4.3.117) • 52	Validator - MFC - Unallowed planned quantity mode - Error (SR3146.9.4.3.113.2) • 58
Validator - MFC - No input defined for output - Error (SR3146.9.4.3.64.2) • 51	Validator - MFC - Unallowed planned quantity mode - Warning (SR3146.9.4.3.113.1) • 58
Validator - MFC - No input defined for output - Warning (SR3146.9.4.3.64.1) • 51	Validator - MFC - Unallowed planned quantity mode (SR3146.9.4.3.113) • 58
Validator - MFC - No input defined for output (SR3146.9.4.3.64) • 51	Validator - MFC - Unallowed weighing method (SR3146.9.4.3.120) • 62
Validator - MFC - No output defined for input - Error (SR3146.9.4.3.63.2) • 51	Validator - MFC - Uniqueness and format of MFC position (SR3146.9.4.3.68) • 57
Validator - MFC - No output defined for input - Warning (SR3146.9.4.3.63.1) • 51	Validator - MFC - Weighing material type not defined (SR3146.9.4.3.123) • 63
Validator - MFC - No output defined for input (SR3146.9.4.3.63) • 51	Validator - MFC (final output) - Planned quantity does not fit master recipe (SR3146.9.4.3.133) • 56
Validator - MFC - None planned quantity mode for Dispense unit procedure (SR3146.9.4.3.115) • 61	Validator - Minimum required equipment parameters - Error (SR3146.9.4.3.83.2) • 66
Validator - MFC - Outgoing or incoming connector is missing (SR3146.9.4.3.138) • 52	
Validator - MFC - Output defined for device-specific recipe or workflow (SR3146.9.4.3.110) • 52	

- Validator - Minimum required equipment parameters - Warning (SR3146.9.4.3.83.1) • 66
- Validator - Minimum required equipment parameters (SR3146.9.4.3.83) • 66
- Validator - Minimum required material parameters - Error (SR3146.9.4.3.10.2) • 56
- Validator - Minimum required material parameters - Warning (SR3146.9.4.3.10.1) • 56
- Validator - Minimum required material parameters (SR3146.9.4.3.10) • 56
- Validator - Non-production MR - More than one unit procedure or operation (SR3146.9.4.3.100) • 65
- Validator - Outgoing link (SR3146.9.4.3.12) • 48
- Validator - Phase parameter - Broken link - Error (SR3146.9.4.3.25.2) • 50
- Validator - Phase parameter - Broken link - Warning (SR3146.9.4.3.25.1) • 50
- Validator - Phase parameter - Broken link (SR3146.9.4.3.25) • 50
- Validator - Phase parameter - ERP BOM attribute discrepancy - Resolve (SR3146.9.4.3.56.1) • 55
- Validator - Phase parameter - ERP BOM attribute discrepancy (SR3146.9.4.3.56) • 55
- Validator - Phase parameter - Invalid reference (SR3146.9.4.3.36) • 50
- Validator - Phase parameter - Missing ERP BOM position (SR3146.9.4.3.50) • 55
- Validator - Phase parameter - Multiple ERP BOM positions (SR3146.9.4.3.51) • 55
- Validator - Phase parameter - Non-ERP BOM position (SR3146.9.4.3.55) • 55
- Validator - Phase parameter - Syntax (SR3146.9.4.3.21) • 50
- Validator - Property definition is inconsistent (SR3146.9.4.3.82) • 66
- Validator - Recipe - Min/Max quantities without UoM (SR3146.9.4.3.111) • 53
- Validator - Recipe - Packaging level does not fit ERP BOM - Resolve (SR3146.9.4.3.58.1) • 54
- Validator - Recipe - Planned quantity does not fit ERP BOM - Resolve (SR3146.9.4.3.57.1) • 54
- Validator - Recipe - Planned quantity does not fit ERP BOM (SR3146.9.4.3.57) • 54
- Validator - Recipe - Reference quantity not defined (SR3146.9.4.3.67) • 53
- Validator - Selection branch (SR3146.9.4.3.15) • 48
- Validator - Simultaneous branch - Unsuitable components (SR3146.9.4.3.17) • 48
- Validator - Simultaneous branch (SR3146.9.4.3.14) • 48
- Validator - Start and end step (SR3146.9.4.3.13) • 48
- Validator - Start step path (SR3146.9.4.3.19) • 49
- Validator - Transition condition - Broken link - Error (SR3146.9.4.3.37.2) • 49
- Validator - Transition condition - Broken link - Warning (SR3146.9.4.3.37.1) • 49
- Validator - Transition condition - Broken link (SR3146.9.4.3.37) • 49
- Validator - Transition condition - Invalid reference (SR3146.9.4.3.38) • 49
- Validator - Transition condition - Syntax (SR3146.9.4.3.31) • 49
- Validator - Transition condition is missing (SR3146.9.4.3.27) • 49
- Validator - Usage type (SR3146.9.4.3.8) • 65
- Validator - Workflow - Only appendable if production-relevant (SR3146.9.4.3.112) • 65
- Validator - Workflow - Phase is incompatible (SR3146.9.4.3.80) • 65
- Validator - Workflow - Workflow type not defined (SR3146.9.4.3.81) • 65
- Validator consistency checks (SR3146.9.4.3+)_Server-run validator - Phase needs to be run on a server (SR3146.9.4.3.104) • 71
- Validator -Packaging level of master recipe (SR3146.9.4.3.107) • 53
- Validator -Packaging level of material parameter (SR3146.9.4.3.108) • 59
- Validator -Recipe - Packaging level does not fit ERP BOM (SR3146.9.4.3.58) • 54
- Version control (SR3146.4.5) • 141
- Configurable version control (SR3146.6) • 142
- Consistency checks (SR3146.6.4.3) • 142

Flexibility of user groups (SR3146.6.7.1) • 143
Flexible definition of approval details (SR3146.6.7) • 143
Flexible version graph assignment (SR3146.6.1) • 143
Period of validity (SR3146.6.4.2) • 142
Serial approval (SR3146.6.7.3) • 143
Status change (SR3146.4.5.1) • 142
Version history (SR3146.4.5.3) • 142

W

Workbench object operations - Link and transition • 32
Workbench object operations (SR3146.9.10+) • 28
Workbench object operations - Element (SR3146.9.10.6) • 31
Workbench object operations - Link and transition • 32
Workbench object operations - Master recipe (SR3146.9.10.1) • 28
Workbench object operations - Master workflow (SR3146.9.10.9) • 29
Workbench object operations - New master recipe dialog (SR3146.9.10.1.2) • 29
Workbench object operations - New master workflow dialog (SR3146.9.10.9.1) • 30
Workbench object operations - Prepare (deep save) (SR3146.9.10.1.3) • 29
Workbench object operations - Prepare (deep save) (SR3146.9.10.9.2) • 30
Workbench object operations - Replace material (product) (SR3146.9.10.11) • 36
Workbench object operations - Smart replace (SR3146.9.10.10) • 32
Workbench object operations - View locked master recipe (SR3146.9.10.1.4) • 29
Workbench object operations - View locked master workflow (SR3146.9.10.9.3) • 31
Workbench UI elements (SR3146.9.2+) • 121
BigDecimal editor • 140
Boolean editor • 140
Comparison - Category filter (SR3146.9.2.19.2) • 131
Comparison - Difference categories (SR3146.9.2.19.1) • 130

Comparison - Difference-specific information (SR3146.9.2.19.5) • 131
Comparison - Navigation (SR3146.9.2.19.4) • 131
Comparison - Type filter (SR3146.9.2.19.3) • 131
Comparison (SR3146.9.2.19+) • 130
Date/Time Picker editor (SR3146.9.2.17.1) • 139
Default preview buildtime (SR3146.9.9.4.1) • 131
Duration editor (SR3146.9.2.17.2) • 139
Expression editor - Explicit null handling (SR3146.9.9.4.9) • 138
Expression editor - Implicit null handling (SR3146.9.9.4.8) • 138
Expression editor - Phase-specific functions (SR3146.9.9.4.11) • 135
Expression editor - Pnuts (SR3146.9.9.4.6) • 138
Expression editor - Runtime context data (SR3146.9.9.6) • 135
Expression editor (SR3146.9.9.4) • 131
Expression editor functions (SR3146.9.9.4.7) • 135
Expression editor literals (SR3146.9.9.4.5) • 134
Expression editor operations (SR3146.9.9.4.4) • 134
Expressions for conditional rules (SR3146.9.9.4.12) • 133
Expressions for flexible rules (SR3146.9.9.4.10) • 133
Expressions for process parameters (SR3146.9.9.4.3) • 133
Expressions for transitions (SR3146.9.9.4.2) • 132
FlexibleAttributeDefinition editor • 140
FlexibleStateModel editor • 140
FlexibleTagDefinition editor • 140
Library filter & search (SR3146.9.2.4.3) • 122
Library for building blocks (SR3146.9.2.4.4) • 123
Library for equipment classes (SR3146.9.2.4.7) • 124
Library for materials (ERP BOM) (SR3146.9.2.4.2) • 123
Library for materials (SR3146.9.2.4.5) • 123
Library for property types (SR3146.9.2.4.8) • 124
Library for signature privileges (SR3146.9.2.4.6) • 124
Library for stations (SR3146.9.2.4.9) • 124
Library for work centers (SR3146.9.2.4.1) • 124
Library Universe (SR3146.9.2.4+) • 122

- Measured value Editor • 139
- Messages - Category filter (SR3146.9.2.16.2) • 129
- Messages - Navigation (SR3146.9.2.16.3) • 129
- Messages - Resolve issue (SR3146.9.2.16.4) • 129
- Messages - Validator-related (SR3146.9.2.16.1) • 129
- Messages (SR3146.9.2.16+) • 128
- Multi-line text editor • 139
- Open dialog (SR3146.9.2.20+) • 122
- Open dialog for proptected objects (SR3146.9.2.20.1) • 122
- Other editors (SR3146.9.2.17+) • 138
- Pacakging Level Data editor (SR3146.9.2.17.3) • 140
- Parameter Panel - Specific property types (SR3146.9.2.7.2) • 127
- Parameter Panel - View all (SR3146.9.2.7.1) • 127
- Parameter Panel (SR3146.9.2.7) • 126
- Parameter Panel filter & search (SR3146.9.2.8) • 127
- Preview panel (SR3146.9.2.11) • 131
- Property window - Element (SR3146.9.2.12) • 128
- Property window - Header (SR3146.9.2.15) • 128
- Property window - Source BB (SR3146.9.2.13) • 128
- Property windows • 127
- Setlist - ERP BOM item (SR3146.9.2.9.1) • 125
- Setlist - Workflow-related consistency checks (SR3146.9.2.9.2) • 125
- Setlist - Workflow-related data adjustment (SR3146.9.2.9.3) • 126
- Setlist (SR3146.9.2.9) • 124
- Statistics (SR3146.9.2.18) • 140
- Statistics view (SR3146.9.2.18.1) • 141
- Statistics-based warning (SR3146.9.2.18.2) • 141
- String editor • 140
- User interface of the Comparison window • 130
- User interface of the Messages window • 128
- Work area (SR3146.9.2.1) • 121
- Workbench usability (SR3146.9.1+) • 117
 - Building block image • 117
 - Context menus (SR3146.9.1.3) • 120
 - Docking (SR3146.9.1.4) • 120
 - Navigation (SR3146.9.1.2) • 119
 - Process view (SR3146.9.1.1) • 119
 - Tooltips (SR3146.9.1.5) • 121
 - Workbench - Start and end steps (SR3146.9.1.1.1) • 119