

Quality Document



PharmaSuite®



Allen-Bradley • Rockwell Software

**Rockwell
Automation**

This document contains information which is confidential and proprietary.

Copyright Notice Copyright © 2017 Rockwell Automation, Inc. All rights reserved.
This document and any accompanying Rockwell Software products are copyrighted by Rockwell Automation Technologies, Inc. Any reproduction and/or distribution without prior written consent from Rockwell Automation Technologies, Inc. is strictly prohibited. Please refer to the license agreement for details.

Trademark Notices FactoryTalk, PharmaSuite, Rockwell Automation, Rockwell Software, and the Rockwell Software logo are registered trademarks of Rockwell Automation, Inc.

The following logos and products are trademarks of Rockwell Automation, Inc.:

FactoryTalk Shop Operations Server, FactoryTalk ProductionCentre, FactoryTalk Administration Console, FactoryTalk Automation Platform, and FactoryTalk Security.

Operational Data Store, ODS, Plant Operations, Process Designer, Shop Operations, Rockwell Software CPGSuite, and Rockwell Software AutoSuite.

Disclaimer ROCKWELL PROVIDES INFORMATION ON AN "AS IS" BASIS. ROCKWELL will not be liable to RECIPIENT for any damages arising out of RECIPIENT's use or reliance on the information contained in this document. Information in this document containing technology strategies and business plans is for planning purposes only. ROCKWELL may change or cancel its plans at any time. Therefore, use of such ROCKWELL information is at RECIPIENT's own risk. Further, this document does not represent any warranty or promise by Rockwell Automation, or any of its business units, to provide any service or support to the use of this document. This document is intended to be executed by the customer.

Contents

Introduction.....	1
Life-Cycle Methodology.....	2
General Development Model.....	2
Requirements Engineering	2
Roles and Responsibilities.....	2
Development: Design and Coding.....	2
Test Strategy.....	3
Validation Activities.....	3
User and Technical Documentation	3
Risk Management	3
Configuration Management	3
Traceability	3
Quality Management	3
Project Monitoring and Control	3
Project Closure	3
Further Development of the System	4
Quality Risks.....	4
Specification Artifacts	5
Test Strategy.....	6
Principles of Test Execution	6
Structure of Test Documentation	6
Document Management System	7
Evaluation of Test Results	7
Repetition of Tests	7
Testing Phases	7
Release Criteria.....	8
Test Environment and Tools	8

Test Plan	8
Test Reports	8
Validation Activities	9
Traceability	10
Quality Management	11
Product and Process Quality	11
Anomaly Management	11
Test Planning	12
Test Approach	12
Testing Scope	14
Requirements Traceability	15
Test Team and Responsibilities	15
Test Environment	16
Test Summary	17
Feature Tests	17
Reliability Tests	17
Regression Testing	17
Final Acceptance Testing	17
Glossary & Abbreviations	18
References	19
Approvals	21
Revision History	22
Appendix A: System Boundary Definition	23
Current Situation	23
Controls for GMP Data	23
PharmaSuite's Boundary	24
Recommendation	25
Appendix B: Test Equipment	26
Server Hardware Components	26
Client Hardware Components	27
Server Software Components	28

Client Software Components	28
Other Hardware.....	28
Appendix C: List of Documents and Approvals	29
Handwritten Signatures for Release Documents.....	34
Handwritten Signatures for Executed Tests.....	35

Introduction

This document summarizes the quality-related planning activities and results for PharmaSuite 8.4 to support further customer-related qualification and validation activities. Other planning activities and results are not covered within this document.

Development of PharmaSuite® followed

- Good Engineering Practices
- Good Testing Practices
- Good Documentation Practices
- Supplier Good Practices

Regarding quality assurance and verification, PharmaSuite® has been developed following a scalable science-based verification process that is fully compliant with the Quality Risk Management as described in ICH Q9 [1] and the risk-based test approach as described in GAMP® 5 [2].

A 2-step approach was used to analyze the system impact and assess the risk of any system change.

The basic development and verification approach as well as the product life-cycle are mainly defined in two internal documents, the Project Plan [18] and the Testing Guideline [15] (with project-specific aspects described in the Test Plan [20]).

Note: Internal documents are confidential and will only be shown to external parties in case a Non-Disclosure Agreement is in place.

The underlying processes are defined in a mature Quality Management System [4] and complemented by various guidelines specific to PharmaSuite®.

Note: Specific terms and abbreviations are explained in chapter *Glossary & Abbreviations*.

The system has been qualified and released to market.

Life-Cycle Methodology

General Development Model

PharmaSuite 8.4 has been developed using an iterative development approach based on (i) iterations (typically lasting three weeks each), (ii) a prioritized product backlog, and (iii) review, retrospective, and planning meetings per iteration, see *Figure 1*. As a result, each iteration produces a potentially shippable product. The last four iterations were reserved for final verification and clean-up activities.



Figure 1: Iterative development approach

Requirements Engineering

The specification artifacts are described in detail in a separate chapter *Specification Artifacts* on page 5.

Roles and Responsibilities

A list of stakeholders and all personnel involved in designing, coding, testing, and documenting PharmaSuite 8.4 is provided in section 3.1 of the Project Plan [18].

Development: Design and Coding

Procedural aspects regarding design reviews, code reviews, the development environment, and tools are documented in the Architecture Specification [26]. Furthermore, the underlying processes for design and code review and their project-specific adaptations are described in sections 5.6.6 and 5.6.7, respectively, of the Development Guideline [14].

Test Strategy

The test strategy is described in detail in a separate chapter *Test Strategy* on page 6.

Validation Activities

These activities are described in detail in a separate chapter *Validation Activities* on page 9.

User and Technical Documentation

For a list of user and technical manuals, please refer to *Appendix C: List of Documents and Approvals*.

User and technical manuals are also available as WebHelp, for most user manuals with context-sensitive access from the PharmaSuite® application. All reviews have been conducted by subject matter experts. Reviews are informal, whereas all approvals are documented via DMS [5] and [6].

Risk Management

Risk management followed the risk management process as described in section 5.10 of the Development Guideline [14].

Configuration Management

Configuration management is described in section 5.11 of the Development Guideline [14] and in the Configuration Management Plan [19].

Traceability

Traceability is described in detail in a separate chapter *Traceability* on page 10.

Quality Management

Quality Management is described in detail in a separate chapter *Quality Management* on page 11.

Project Monitoring and Control

More details on Project Monitoring and Control, based upon the approach described in chapter *Life-Cycle Methodology*, can be found in section 5.13 of the Development Guideline [14].

Project Closure

The formal approval of the Project Closure Report [25] officially closed the project.

Further Development of the System

The remaining anomalies will be evaluated for resolution in future versions of PharmaSuite®.

Quality Risks

Known issues relevant for end users have been documented in the project's Release Notes [24].

Note: The qualification of PharmaSuite 8.4 has been completed following a risk-based test approach focusing on all changes made to the PharmaSuite® system since the previous release. Thus, all qualification activities are based on PharmaSuite 8.3 and rely on the quality status of this release. All changes have been analyzed and tested with regard to their risk and impact on the system.

Specification Artifacts

Requirements for PharmaSuite 8.4 have been developed in a 2-step approach.

In a first step, certain artifacts, using various tools (e.g. mind maps, graphics, drawings, text documents), have been created as supporting documents for items on the Product Backlog (PBL). This has been done in order to support the required effort estimation and knowledge dissemination.

As soon as a certain backlog item had been committed for implementation within an iteration, more specific details were provided. During iterations, a spreadsheet has been used to hold the entire set of requirements. Finally, all of the relevant specification details, including use cases and requirements, have been compiled within the Supported Platforms Guide [27], the Detailed Requirements Specification (DRS) documents ([28] – [32]) and the Functional Requirement Specification (FRS) documents ([33] – [49]).

Requirements are subject to version control and base-lined per approved document. Baseline requirements are modified using change control per relevant QMS procedures. Changes have been tracked through the change control mechanism of the respective documents.

As a result, the Supported Platforms Guide and all DRS and FRS documents of PharmaSuite 8.4 ([27] – [49]) include all requirements and list all changes that apply to this release (compared to previous releases).

All specification documents have been posted on the Document Management Systems (DMS; either [5] or [6]) for formal approval.

Test Strategy

The software test process is generally described in the Software Testing Procedure [9]. The document outlines the testing process as well as the entry and exit criteria of the different test types.

This chapter describes the framework requirements and methods of test planning. The detailed planning of test activities and their execution have been described in a Risk-Based Testing Guideline [15], with project-specific deviations defined in the Test Plan [20].

Principles of Test Execution

The objective of testing is the verification of the correspondence between the system and its specification. In order to achieve this objective, testing tried to reveal any deviation in the implemented functionality against requirements, use cases, and other specification artifacts. Anomalies found in the scope of testing have been recorded and analyzed. The majority has been resolved and verified according to the change control procedures and guidelines ([16] and [17]) until a high level of confidence had been assured that the system meets the specified requirements.

Structure of Test Documentation

The highest level document is the Project Plan [18]. It lays out the general framework of requirements for test execution and for the structure of test plans and individual testing phases as well as the general test scope.

The Risk-based Testing Guideline [15] in combination with the project-specific Test Plan [20] constitutes the second level. It describes all activities of the testing phase. The testing phase reflects important milestones in the product development. A Test Summary Report was generated at the end of the testing phase [21]. The Test Summary Report corresponds to the Test Plan and gives a summary of all test activities and their results for that test phase.

Third level items are the System Impact Analyses, Risk Assessments, as well as Test Cases. The System Impact Analysis along with the Risk Assessment identifies and assesses the risk of changes related to the system. Based on the results, the required test scope has been derived. Details are defined in the Risk-Based Testing Guideline [15].

Prior to test execution, all test documents were reviewed and approved.

Document Management System

All quality-relevant documents and project documentation are subject to the organization's Document Management Systems (either [5] or [6]). Each document is numbered and version-controlled.

All documentation of PharmaSuite 8.4 is generally stored following Rockwell Automation's Record Retention Procedure [13].

Evaluation of Test Results

A combination of the type and frequency of anomalies is applied as metrics for the software quality during the testing phase. This requires a suitable classification of anomalies.

The triage team is responsible for conducting an initial analysis on the anomaly submitted to determine the appropriate FMEA value following the FMEA guidelines as described in the Defect Management Guideline [16].

Repetition of Tests

The repetition of a testing phase after a correction cycle followed the principle of regression testing (for details see Software Testing Procedure [9]). The tests to be executed during the repeated execution of a testing phase usually represent a subset of the original test cases. The selection process was geared by the System Impact Analysis and Risk Assessment, as described in detail in the Risk-Based Testing Guideline [15].

Implementation of corrective actions and regression testing was repeated until all anomalies violating the acceptance criteria (see [16] for details) were corrected and verified.

Testing Phases

The testing of PharmaSuite 8.4 was conducted in phases. The typical course of testing was as follows:

- Creation of the Test Plan [20] (incl. review and approval)
- Creation or update of Test Cases [22] (incl. review and approval)
- Test execution and recording of results
- Classification of detected anomalies and definition of corrective actions
- Creation of a report summarizing the test results of the testing phase
- Review and approval of the Test Summary Report [21]

Prior to testing, an installation qualification for the software has been carried out with the objective to document the proper installation of software for all work stations and servers used for the tests. This was executed based on the guidelines defined in the Risk-Based Testing Guideline [15].

Release Criteria

The test phase was accepted when all defined test cases had been completed (incl. review and approval) and the results of the execution of all tests were accepted and thereby no CRs violating the project's release criteria (as defined per [16]) were left.

Note: PharmaSuite 8.4 provides add-on functionality to ProductionCentre and thus uses existing functionality of ProductionCentre.

During verification and qualification of PharmaSuite 8.4, testing has been focused on PharmaSuite 8.4 requirements. ProductionCentre was considered to be a qualified baseline for PharmaSuite 8.4.

Test Environment and Tools

The list of hardware and software used for testing has been captured in the Test Summary Report [21] and included in this document in *Appendix B: Test Equipment*.

Test Plan

A project-specific Test Plan [20] has been created, reviewed, and approved before start of test. See chapter *Test Planning* for more details.

Test Reports

After finalization of the system test, a report summarizing the test results of the testing phase has been created, reviewed, and approved [21]. See chapter *Test Summary* for more details.

In addition, all test cases have been listed in the Final Test Status Report [22].

Validation Activities

Note: The term validation refers to the definition of the Validation Procedure [10], based upon the CMMI definition [3]. It does neither refer to the definition of GAMP [2] (Computer System Validation: Achieving and maintaining compliance with applicable GxP regulations and fitness for intended use (...)) nor to regulations like 21 CFR 820.70(i), as a supplier can obviously not fulfill these requirements. However, PharmaSuite 8.4 was launched in a qualified status, thereby allowing customers to leverage both, Rockwell Automation's verification and validation activities.

All validation activities for PharmaSuite 8.4 followed the Validation Procedure [10]. In general, validation has been performed by extended peer reviews while creating new and revising existing requirements. The peer reviews have been performed

- by relevant stakeholders of all disciplines (Requirements Engineering, Product Management, Development, and Testing) involved in the project;
- on detailed requirements (or other comparable specification artifacts like Use Cases etc.) level;
- according to the Peer Review Procedure [11].

Furthermore, newly developed features have been presented to the engineering team, to the product owners, and to integrators after completion of each iteration. This review-per-demonstration process ensured to receive feedback on new features as soon as possible. Also, intermediate product builds including new features have been shown to prospective customers as well.

Traceability

Traceability is maintained partially within a traceability matrix document, and partially within the different types of documents:

- SRs reference to MRs.¹
- Work items reference to related SRs implemented and anomalies fixed.
- Code and design review protocols are documented at work items (and defects, as needed) in the defect management tool.
- Each work item references in which build it was integrated.²
- Work items and / or anomaly CRs list the changes to the system.
- Test cases reference to the SRs they are going to verify.

The traceability matrix document [23] is maintained per iteration as part of the “Done” state, and finally uploaded to the Document Management System [6] for approval (prior to project completion). This document shows the mapping of SRs to test cases, thereby helping to assure that all requirements have been verified.

¹ For releases 1.0 and 1.1 only, i.e. no new market requirements have been developed for PharmaSuite 8.4 and traceability from SRs to MRs has not been maintained anymore for PharmaSuite 8.4. Rationale: By using the iterative life-cycle methodology and the product backlog, a different method is applied to gear the prioritization of feature implementation. Furthermore, MRs have not been used for formal tracking and have not been part of the formal testing approach, but were used for documentation and project control purposes only.

² Build numbers are increasing over time. Beware: Work items may reference to a build that did not match acceptance criteria. However, the next number of a build meeting acceptance criteria also contains that work item.

Quality Management

Product and Process Quality

Product and process quality is assured through regular assessments.

Anomaly Management

Anomaly tracking has been done within Jira per the Defect Management Guideline [16]. The triage team reviewed and evaluated anomalies in regular meetings (frequency as needed) according to the Defect Management Guideline [16]. The triage team was led by the Test Manager. In this role, the Test Manager had the final decision in case of disagreement of the triage team. Further members of the triage team can be found in section 3.1 of the Project Plan [18].

Test Planning

This chapter discloses important information from the Risk-based Testing Guideline [15]. Parts that have not been transferred to this document may be reviewed during audits, as needed.

Test Approach

Testing was conducted to verify the stability of the product and the functionality based on the specification artifacts (e.g. requirements, use cases, etc.). The following sub-sections describe the types of testing that have been conducted on PharmaSuite 8.4. The definitions for each testing type are defined in the Software Testing Procedure [9] and the Software Process Document Glossary [8].

For the verification of PharmaSuite 8.4, a risk-based testing approach has been applied (see [15]). In general the approach is defined as follows:

- Each new feature and functionality had to be analyzed through a risk assessment to determine the risk priority.
- All system changes (i.e. the corresponding Change Requests) needed to be analyzed within a System Impact Analysis (SIA) and a risk assessment had to be conducted to determine the necessity and rigor of testing.
- The risk-based approach was not applied to installation tests, installation qualification, unit tests, performance tests, and the final acceptance tests, as a variation of the test depth (level of detail) was not foreseen for these types of test cases.

Based on the various risk assessments, test cases have been created or revised to define the individual test scope. Sometimes, several features or functionalities have been combined into one test scenario.

Risk Evaluation

The risk-based approach is considered to establish a scalable science-based verification process and to ensure that the work and testing effort spent focuses on risky areas. Therefore, the risk evaluation was divided in two interrelated processes:

1. The first step of the risk evaluation is to identify all system areas affected with potential risk within the System Impact Analysis (SIA).
2. The second step is to rate the potential risk to get the risk priority as an indicator for the resulting testing or mitigation strategy.

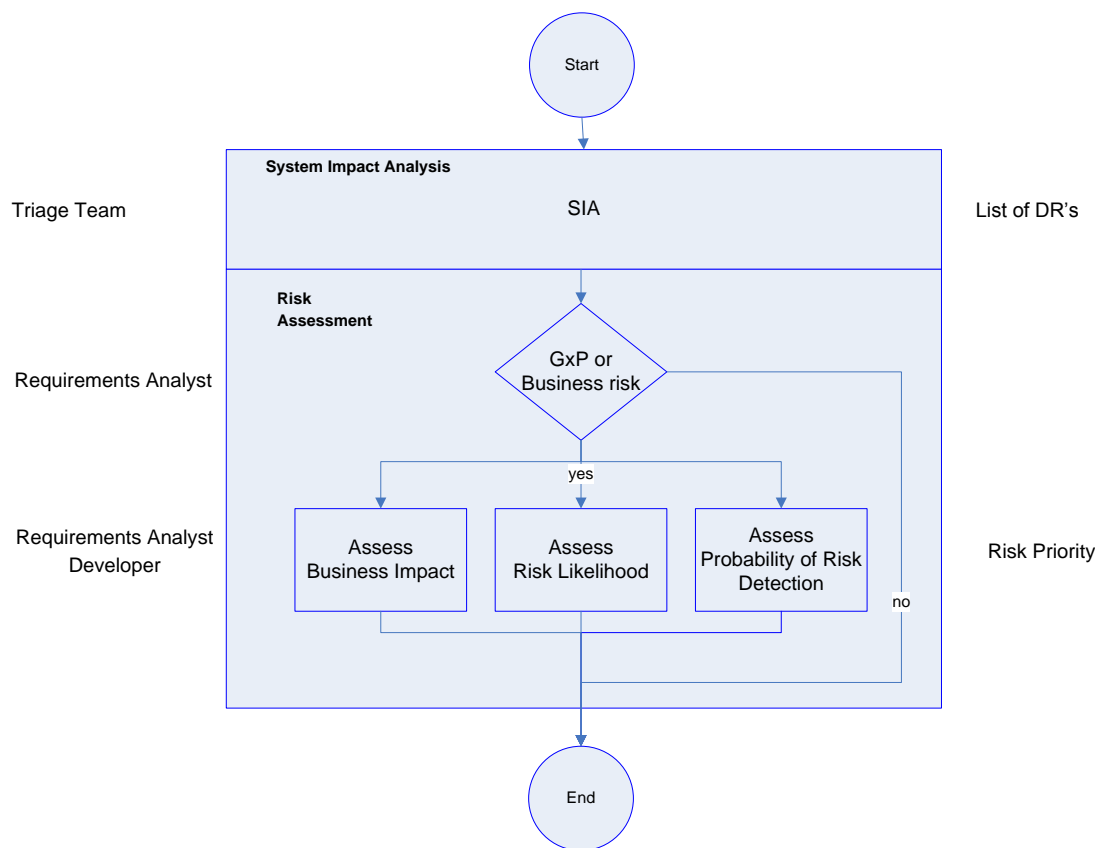


Figure 2: Risk Evaluation Process

SYSTEM IMPACT ANALYSIS (SIA)

All changes made to the system have to be listed and analyzed within the System Impact Analysis. This means every change (e.g. bug fixes, new functionalities) needs to be analyzed to identify all system areas with potential risk (e.g. by side effects) caused by the changes. The System Impact Analysis identifies affected system areas as a first indication of the appropriate scope of testing.

RISK ASSESSMENT (RA) PROCESS

The risk assessment aims to classify the risk to the different functionalities of PharmaSuite 8.4 using pre-defined assessment criteria.

Adapted from this, a Subject Matter Expert derived the appropriate risk priority for each function. The risk assessment was performed based on specification artifacts identified during the System Impact Analysis process.

The assessment of the GxP and Business Risk as well as the Business Impact are artifact-specific attributes and therefore independent of the change. These classifications have been documented per specification artifact in the System Impact Analysis (SIA) & Risk Assessment (RA).

Guidelines for risk assessment are described in detail in the Testing Guideline [15].

Testing Scope

Unit Tests

JUnit has been used for unit testing where applicable (up to and including the service layer). Therefore, the public interfaces of MES components as described in the Architecture Specification [26] have been tested. Unit testing followed the Unit Testing Guideline [12]. The tests were managed within Subversion like ordinary source code.

More details regarding Unit Tests can be found in the Risk-based Testing Guideline [15].

Installation Test

The installation test was aimed to check the correct implementation and description of the PharmaSuite 8.4 installation process including the configuration for all components.

The installation process was based on the description of PharmaSuite 8.4's Installation Manual [50].

Feature Tests

For the testing of PharmaSuite 8.4 new features and functionalities, a risk assessment has been created in which the potential risks of the features and the impact of it with regard to existing functions (SIA) has been evaluated ([51] and [52]).

The feature tests were based on the specification artifacts (mainly requirements). For results, refer to section *Feature Tests* in chapter *Test Summary*, page 17.

Note: Parts of the feature tests have been realized as automated tests. However, independent of the used technique, there is no difference to the claim of the tests. Both, automated as well as manual tests have been treated as equally qualified test methods.

Performance Tests

Performance tests (as well as stress tests, boundary tests, limit tests, and response tests) are implicitly covered within the feature tests. The performance test scope was covered within the System Impact Analysis [51] and [52].

Reliability Tests

The objective of reliability testing was to show the system running 12+ hours, performing operations like insertions and deletions of master data. Therefore, special scripts have been run to verify the reliability as described. Reliability tests have been performed in the qualified test environment.

Integration Tests

Integration tests for new features or functionalities as well as CRs have been covered within the feature or regression tests. The integration test scope was covered within the System Impact Analysis [51] and [52].

System Tests

System tests as described in the Software Testing Procedure [9] have been covered within the corresponding risk evaluation and the resulting regression tests. Therefore, no additional system tests were required.

Customer Acceptance Testing

Customer acceptance tests were not planned for PharmaSuite 8.4.

Regression Tests

Regression testing was performed by re-running feature tests that had been executed in previous releases against the current system. The objective was to ensure that the functionality of the system continues to work as expected after all features have been implemented. The scope has been derived by and described in the respective risk evaluation ([51] and [52]).

Note: Parts of the regression tests have been realized as automated tests. However, independent of the used technique, there is no difference to the claim of the tests. Both, automated as well as manual tests have been treated as equally qualified test methods.

CR verification

Change requests (CRs) integrated in the system were objects of this test, as long as they were not covered by other tests. The verification of CRs ensures that anomalies have been resolved satisfactorily.

Anomalies violating the defined system acceptance criteria had to be verified on a controlled test environment.

Final Acceptance Tests (FAT)

The final acceptance test represented a high-level, overarching test scenario in the course of PharmaSuite 8.4 qualification. It examined the functional capability of the system installed based on the final build³. The objective of the final acceptance test was to examine the system applicability in the operational area and to confirm the product satisfies the criteria for the customer.

Integrated Solutions Testing

Within this project, no integrated solution tests were planned.

Requirements Traceability

A traceability matrix [23] has been generated to trace the test cases to the appropriate detailed requirement(s) and other specification artifacts. The matrix has been developed throughout the testing activities as new test cases were created and executed or as specifications changed.

Test Team and Responsibilities

The test team members and responsibilities are called out in the Project Plan [18].

³ “Final” in the sense that this is the final build of the current release, the build that is going to be shipped.

Test Environment

For the execution of tests, dedicated test hardware (clients and servers) available in the test lab of Rockwell Automation campus Karlsruhe has been used. The test lab provides a segregated network managed by the test team. All test hardware used for official tests has been installed according to the definitions in the detailed requirement specifications.

For more details, refer to *Appendix B: Test Equipment*.

Installation Qualification

The objective of the installation qualification was to check the suitability of client and server hardware as far as dedicated to PharmaSuite 8.4, as well as software components related to PharmaSuite 8.4 installation. The installation qualification was performed on the test environment before start of formal testing.

For this specific test, no acceptance criteria were defined. The software installed on the clients and servers must correspond to the definitions in the Supported Platforms Guide [27]. Otherwise, the work station or server must not be used within the official test. The software installation qualification refers to all clients and servers mentioned in the Supported Platforms Guide [27].

Test Summary

This chapter discloses important information from the Test Summary Report [21]. Parts that have not been transferred to this document may be reviewed during audits, as needed.

Feature Tests

226 test cases have been completed for new feature testing in 6 test sessions. Of those 226 tests, 220 tests passed. The test coverage achieved for the specific requirements addressed in the release is 100% for all new functionalities, i.e. all requirements are covered by one or more tests (see [22]).

Reliability Tests

One test running 12+ hours has been successfully executed.

Regression Testing

49 test cases have been completed for regression testing in 4 test sessions. Of those 49 tests, 48 passed.

Final Acceptance Testing

26 test cases have been completed for final acceptance testing in 1 test session. Of those 26 tests, all tests passed.

Glossary & Abbreviations

Term	Definition
<i>CMMI</i>	Capability Maturity Model Integration A process improvement approach that provides organizations with the essential elements for effective process improvement.
<i>CR</i>	Change Request A record that describes a request to change the system and may be an anomaly, an enhancement, etc.
<i>Done state</i>	Final state of a product backlog item indicating that all related tasks (e.g. specification, implementation, test, documentation) have been completed.
<i>DRS</i>	Detailed Requirements Specification Collection of detailed software requirements (SRs)
<i>FRS</i>	Functional Requirements Specification Includes detailed software requirements (SRs)
<i>IDS, IDSs</i>	Integrated Design Specification(s)
<i>Jira</i>	Anomaly management system used by Rockwell Automation
<i>MES</i>	Manufacturing Execution System
<i>MR, MRs</i>	Market Requirement(s) (high level)
<i>PBL</i>	Product Backlog Prioritized list of items to be realized in PharmaSuite
<i>SR, SRs</i>	Software Requirement(s)
<i>Triage</i>	Process of deciding upon further processing of CRs (resolving, deferring, etc.).
<i>Triage team</i>	Team of subject matter experts, led by the Test Manager, performing the triage process.

References

- [1] ICH Q9 Quality Risk Management; International Conference On Harmonisation Of Technical Requirements For Registration Of Pharmaceuticals For Human Use; Step 4, version 09-Nov-2005
- [2] GAMP® 5: A Risk-Based Approach to Compliant GxP Computerized Systems; ISPE; Feb 2008
- [3] CMMI Guidelines for Process Integration and Product Improvement, Mary Beth Chrissis, Mike Konrad, Sandy Shrum, SEI
- [4] Rockwell Software Quality Management System (QMS),
<https://rockwellautomation.sharepoint.com/teams/AS/ISPB/PMO/MES/pharmasuite/LSQMS/>
- [5] Rockwell Software Document Management System (DMS),
<http://softwaredms.na.home.ra-int.com>
- [6] Rockwell Automation Document Management System (DMS), basic functionality included in central SAP
- [7] Rockwell Automation Testing Environment, <https://aptest.ra.rockwell.com/atm/login.mpl>

Copyright of all following documents by Rockwell Automation, Inc.:

- [8] Software Development Process Glossary, Doc No. 102295, Rev. 2.1
- [9] Software Testing Procedure, Doc No. 102378, Rev. 3.4
- [10] Validation Procedure, Doc No. 102634, Rev. 2.1
- [11] Peer Review Guideline, Doc No. 102556, Rev. 1.5
- [12] Unit Testing Guideline, Doc No. 101890, Rev. 1.2
- [13] Record Retention Procedure, Doc. No. 103837, Rev. 2.2
- [14] LS Release Development Guideline, Doc No. 105112, Rev. 3.0
- [15] LS Risk-Based Testing Guideline, Doc No. 104706, Rev. 2.0
- [16] LS Defect Management Guideline, Doc No. 105543, Rev. 2.0
- [17] LS Work Item Management Guideline, Doc No. 105547, Rev. 2.0
- [18] PharmaSuite 8.4 – Project Plan, DIR 10003202609/PRJ, Rev. 1.3
- [19] PharmaSuite 8.4 – Configuration Management Plan, DIR 10003202610/PLN, Rev. 1.1
- [20] PharmaSuite 8.4 – Test Plan, DIR 10003202611/PLN, Rev. 1.0
- [21] PharmaSuite 8.4 – Test Summary Report, DIR 10003547210/REV, Rev. 1.0
- [22] PharmaSuite 8.4 – Final Test Status Report , DIR 10003547212/REV, Rev. 1.0
- [23] PharmaSuite 8.4 – Traceability Matrix, DIR 10003202613/FRM, Rev. 1.0
- [24] PharmaSuite 8.4 – Release Notes, DIR 10003224098/PUB, Rev. 1.1

- [25] PharmaSuite 8.4 – Project Closure Report, DIR 10003547206/REV, Rev. 1.0
- [26] PharmaSuite 8.4 – Architecture Specification, DIR 10003202612/SPC, Rev. 1.0
- [27] PharmaSuite 8.4 – Supported Platforms Guide, DIR 10003202615/SPC, Rev. 1.1
- [28] PharmaSuite 8.4 – DRS Enhancements, DIR 10003202616/SPC, Rev. 1.0
- [29] PharmaSuite 8.0 – DRS Information Management, Doc No. 106810, Rev.1.0
- [30] PharmaSuite 8.0 – DRS Dispense, Doc No. 106811, Rev. 1.0
- [31] PharmaSuite 8.0 – DRS Inventory Management, Doc No. 106815, Rev. 1.0
- [32] PharmaSuite 8.0 – DRS Production Management, Doc No. 106817, Rev. 1.0
- [33] PharmaSuite 8.4 – FRS Execution Framework, DIR 10003236391/SPC, Rev. 1.1
- [34] PharmaSuite 8.4 – FRS Data Management, DIR 10003224063/SPC, Rev. 1.0
- [35] PharmaSuite 8.4 – FRS Production Execution, DIR 10003526985/SPC, Rev. 1.0
- [36] PharmaSuite 8.4 – FRS Runtime Data Management, DIR 10003253129/SPC, Rev. 1.0
- [37] PharmaSuite 8.4 – FRS Recipe and Workflow Management, DIR 10003236392/SPC, Rev. 1.1
- [38] PharmaSuite 8.4 – FRS Review and Approval, DIR 10003224065/SPC, Rev. 1.0
- [39] PharmaSuite 8.4 – FRS Non-functional Requirements, DIR 10003256030/SPC, Rev. 1.0
- [40] PharmaSuite 8.4 – FRS Workflow Phases, DIR 10003447532/SPC, Rev. 1.0
- [41] PharmaSuite 8.4 – FRS Material Tracking Phases, DIR 10003303218/SPC, Rev. 1.0
- [42] PharmaSuite 8.4 – FRS IPC Phases, DIR 10003303216/SPC, Rev. 1.1
- [43] PharmaSuite 8.4 – FRS Dispense and Inline Weighing, DIR 10003224073/SPC, Rev. 1.0
- [44] PharmaSuite 8.4 – FRS Output Weighing, DIR 10003224077/SPC, Rev. 1.0
- [45] PharmaSuite 8.4 – FRS Equipment Automation Phases, DIR 10003303213/SPC, Rev. 1.0
- [46] PharmaSuite 8.4 – FRS Equipment Tracking Phases, DIR 10003303214/SPC, Rev. 1.0
- [47] PharmaSuite 8.4 – FRS EBR Phases, DIR 10003303200/SPC, Rev. 1.0
- [48] PharmaSuite 8.4 – FRS DCS Phases, DIR 10003246501/SPC, Rev. 1.1
- [49] DCS 2.1 – FRS DCS Adapter, DIR 10003246500/SPC, Rev. 1.1
- [50] PharmaSuite 8.4 – TM Installation – Enterprise Edition, DIR 10003224081/PUB, Rev. 1.2
- [51] PharmaSuite 8.4 – Artifact Specific Risk Assessment, DIR 10003202617/REV, Rev. 1.2
- [52] PharmaSuite 8.4 – System Impact Analysis (SIA) & Risk Assessment (RA), DIR 10003202619/REV, Rev. 1.2
- [53] Implementation of 21 CFR Part 11 – PharmaSuite 8.4, DIR 10003305443/PUB, Rev. 1.0

Approvals

Approvals are captured electronically on the organization's Document Management System [6]. The required approvers of this document include the following:

Name	Role
Martin Dittmer	Product Manager
Steffen Landes	Development Manager
Martin Irmisch	Test Manager

Revision History

Revision History

The following table describes the history of this document. Each version has been approved per Document Management System [6].

Version	Author	Description
1.0	Martin Irmisch	Initial document creation

Appendix A: System Boundary Definition

Appendix A: System Boundary Definition

Current Situation

PharmaSuite provides audit trail as required by regulations in different ways, see “Implementation of 21 CFR Part 11 – PharmaSuite 8.4” [53].

In addition, raw data is recorded for creation, modification, and deletion of objects in *Process Designer*. However, this data is neither exposed nor usable as is for audit trail review, reporting etc., but requires additional steps to become meaningful (*not subject to standard PharmaSuite product offering*).

Audit trail data for objects in *Process Designer* are not required by any regulation⁴. The rationale is explained in further detail in subsequent sections.

Process Designer allows to implement a fine-granular authorization scheme on two levels: Access can be granted or denied on *object type* (e.g. "Lists") or *individual object* (e.g. list "CostCenters") level. Furthermore, different privileges can be assigned for different actions like creation, modification, and deletion.

Controls for GMP Data

The basic question is: “Do all changes to data and records that (may) have an impact on data integrity or product quality need to be subject to a computer-generated audit trail or not?” From a compliance perspective, the answer is clearly “No”.

Instead, a risk assessment should be applied. The result of this risk assessment identifies the critical data and records, describes the criticality and risk of uncontrolled changes and deletion, and defines risk management measures to eliminate, mitigate, or accept an identified risk.

Implementing (or activating) an audit trail is one amongst several possible measures, but clearly not the only one. Alternative measures must provide adequate controls and include, but are not limited to

1. the application of rigorous, documented change or version control (in combination with access control) or
2. security measures to prevent or avoid accidental and unauthorized changes (e.g. storage on read-only media).

⁴ With the exception of users, user groups, and access privileges, see section *Exceptions*.

PharmaSuite's Boundary

Drawing the Line

As any other system, PharmaSuite has a boundary. Within this boundary, PharmaSuite fully complies with GMP aspects, specifically the tracking of changes to GMP data that may have an impact on data integrity or product quality.

Some data that is managed in PharmaSuite – it should be distinguished from other data managed in other applications (e.g. *Process Designer*), thereby defining the boundary. This chapter should help with the disambiguation.

Inside

All master and runtime data that can be created, altered, or deleted *by means of the PharmaSuite UI* is considered *inside*. “Implementation of 21 CFR Part 11 – PharmaSuite 8.4” [53] describes the various manifestations of audit trail for this data – the ‘classical’ audit trail for material, the transaction history for inventory objects, the change history and logbook for equipment, the batch record for produced batches, etc. All this data resides *inside* the boundary of PharmaSuite.

Outside

All configuration and build-time data that *cannot* be created, altered, or deleted by means of the PharmaSuite UI (this data is typically setup once during system deployment) is considered *outside*. This data is subject to system validation and installation qualification processes. Changes in the operational phase may be recorded and made available via ‘audit trail’, but this is not a regulatory requirement. However, changes need to be controlled by means of good engineering practices.

This applies to many objects managed in *Process Designer* – incl. the application configuration and lists, but also forms, reports, flexible state models etc. As changes to these objects can have a severe impact on the PharmaSuite application itself, access to this data needs to be strictly limited to dedicated persons. This is no different than access control to other deployment-relevant components like database or application server.

Example

As an example and to highlight the different dimension and scope of a change, alterations to a master recipe can be made visible by the comparison capabilities in PharmaSuite and will be reviewed and approved before going into production. While an error (that would affect all future orders based on this master recipe) would already be serious, all other master recipes are not affected at all when changing a master recipe. In contrast to this, changes to the *MESMasterRecipeVersionGraph* flexible state model (FSM; the default version graph for master recipes) have a global impact on the review and approval workflow of all master recipes. This kind of change clearly requires different handling.

Exceptions

Users, user groups, and access privileges are an exception to the above definition: While they are managed by *Process Designer's UI*, they are included by definition in the scope of PharmaSuite. As a consequence, an audit trail-specific UI is available in PharmaSuite, and rigorous access control needs to be in place to limit an administrator's access to these objects.

Recommendation

When combining the two previous sections, it becomes obvious that (changes to) GMP data *inside* PharmaSuite is handled differently compared to (changes to) GMP data *outside* this boundary: Whereas a comprehensive audit trail is available for the first, data outside of PharmaSuite requires rigorous, documented controls.

The basic assumption is that changes to configuration and build-time data are not daily business, but rather exceptions, justifying a more rigorous and securer process.

Implementation of Controls

The various options to set up, administer, and configure PharmaSuite are described in several *Technical Manuals*. After the initial system setup (and validation), any changes to configuration and build-time data made outside of PharmaSuite require a defined change control process.

Separation of Duties

The change control process applied to data in *Process Designer* should be reinforced by making adequate use of *Process Designer's* access privileges, thereby implementing a scheme to support clear *separation of duties*. As an example, an administrator that manages users and user groups does not need privileges to change any other object. On *individual object* level, a person may be able to modify a *custom* list (e.g. add a new list entry), but must not be allowed to alter a *system* list.

Appendix B: Test Equipment

Appendix B: Test Equipment

The following sections describe the hardware and software components used for testing of PharmaSuite 8.4. This does not mean that other hardware or software is unsuitable to run PharmaSuite 8.4.

In general, PharmaSuite 8.4 can run on any standard hardware that is sufficiently equipped regarding CPU, RAM, etc., so the listed servers and clients should be taken as examples only.

For suitable software components, please refer to the Supported Platforms Guide [27].

If different hardware or software components are used, the impact and risk of the different components should be evaluated and testing conducted, as appropriate.

Server Hardware Components

All testing has been performed with logical servers running on Virtual Machines (VMs) using VMWare ESXi (nine logical server VMs for application (JBoss) and database (Oracle) servers and four server VMs for migration, running both). The five physical servers used for the system tests are listed in the following table:

Server Name	CPU ⁵	HDD	RAM
QAS003	2x Intel XEON 5110 / 1.6 GHz (4 cores)	400 GB	20 GB
QAS004	2x Intel XEON 5110 / 1.6 GHz (4 cores)	350 GB	20 GB
EUDEKARQAESXI1	2x Intel XEON X5660 / 2.8 GHz (24 cores)	1.9 TB	96 GB
EUDEKARQAESXI2	2x Intel XEON E5-2640 / 2.5 GHz (24 cores)	1.9 TB	262 GB
EUDEKARQAESXI4	2x Intel XEON E5-2620 / 2.4 GHz (24 cores)	5.46 TB	130 GB

The virtual machines for the twelve logical application (JBoss) and database (Oracle) servers running on the above-mentioned HW servers have been configured as follows:

VM Name	Server Name	Purpose	HDD	RAM
QAS001	EUDEKARQAESXI2	JBoss	40 GB	4 GB
QAS002	QAS003	JBoss	40 GB	4 GB
QAS003 ⁶	QAS004	JBoss	40 GB	4 GB
QAS004 ⁹	EUDEKARQAESXI2	Oracle	40 GB	4 GB
QAS005	QAS003	Oracle	40 GB	4 GB

⁵ The number of cores is based on hyper-threading, i.e. the physical number of cores is half of it.

⁶ The “identical” names do not create any conflict, as the name spaces of logical VMs and HW servers is isolated.

VM Name	Server Name	Purpose	HDD	RAM
QAS006	QAS004	Oracle	40 GB	4 GB
QAS008	EUDEKARQAESXI1	Oracle	40 GB	4 GB
QAS012	EUDEKARQAESXI2	JBoss	40 GB	4 GB
QAS013	EUDEKARQAESXI2	JBoss	40 GB	4 GB
QAS014	EUDEKARQAESXI4	Oracle	40 GB	4 GB
QAS015	EUDEKARQAESXI4	JBoss	40 GB	4 GB
QAS016	EUDEKARQAESXI4	JBoss	40 GB	4 GB

The virtual machines for the four logical migration servers running on the above-mentioned HW servers have been configured as follows:

VM Name	Server Name	Purpose	HDD	RAM
QAS007	EUDEKARQAESXI4	JBoss + Oracle	40 GB	8 GB
QAS009	EUDEKARQAESXI4	JBoss + Oracle	40 GB	8 GB
QAS010	EUDEKARQAESXI4	JBoss + Oracle	40 GB	8 GB
QAS011	EUDEKARQAESXI4	JBoss + Oracle	40 GB	8 GB

Client Hardware Components

The physical VMWare ESXi servers EUDEKARQAESXI1, EUDEKARQAESXI2, and EUDEKARQAESXI4 (see above) also host the 22 client VMs (EUDEKARQAESXI1 hosts 11 logical clients, EUDEKARQAESXI2 seven, and EUDEKARQAESXI4 4). All logical clients are configured with disk space (HDD) of 20 GB and memory (RAM) of 2 GB.

In addition, 2 physical clients have been used, e.g. for tests requiring access to physical hardware like barcode scanner or scale. Both physical clients are equipped with an Intel Core i5-2400 3.1 GHz CPU, disk space (HDD) of 250 GB, and memory (RAM) of 8 GB.

Server Software Components

Please note that the language configuration of the operating system for the test system must be English⁷.

Database Server

Type of Software	Software and Version
Operating System	Microsoft Windows 2012 Server R2 Standard
DBMS	Oracle 12c Server (Enterprise), 12.1.0.2

Application Server

Type of Software	Software and Version
Operating System	Microsoft Windows 2012 Server R2 Standard
Application Server	JBoss Enterprise Application Platform 7.0

Client Software Components

Type of Software	Software and Version
Operating System	Microsoft Windows 10 Professional
Browser	Internet Explorer 11

Other Hardware

Type of Hardware	Supplier	Type
Printer	Hewlett Packard	LaserJet 4350dtn
	PDFCreator	PDFCreator ⁸
Label Printer	Zebra	110 Xi III Plus
Scales	Sartorius	LA 2200
Scanner	Honeywell	4820
	Handheld	3800g

⁷ Locale = English (U.S.)

⁸ Printer emulator / soft printer to store PDF files.

Appendix C: List of Documents and Approvals

Appendix C: List of Documents and Approvals

Many documents of PharmaSuite 8.4 are listed in chapter *References*.

Table 1 below lists *all* documents that have been created or updated within the scope of PharmaSuite 8.4. The title has been abbreviated in this list, as the prefix “PharmaSuite 8.4 –” applies to all documents⁹.

Approval of all documents was performed electronically via Rockwell’s Document Management System [6], which also manages information regarding revision numbers as well as list of approvers.

#	Title	DIR	Revision ¹⁰	Approval ¹¹	S. Muench	M. Irmisch	S. Landes	M. Dittmer
1	Project Plan	10003202609/PRJ	1.3	18-Dec-2017	–	✓	✓	✓
2	Configuration Management Plan	10003202610/PLN	1.1	17-Dec-2017	–	✓	✓	–
3	Test Plan	10003202611/PLN	1.0	03-Jul-2017	✓	–	✓	✓
4	Glossary	10003224097/PUB	1.0	01-Dec-2017	–	✓	✓	✓
5	Architecture Specification	10003202612/SPC	1.0	21-Dec-2017	–	✓	✓	–

⁹ Except for the documents related to *DCS 2.1*, references 18, 35, and 57.

¹⁰ The revision shown is the latest document revision that has been approved within the course of PharmaSuite 8.4.

¹¹ The date shown is the date of the latest given approval for the indicated revision.

#	Title	DIR	Revision ¹⁰	Approval ¹¹	S. Muench	M. Irmisch	S. Landes	M. Dittmer
6	Supported Platforms Guide	10003202615/SPC	1.1	21-Dec-2017	–	✓	✓	✓
7	DRS Enhancements	10003202616/SPC	1.0	10-Nov-2017	–	✓	✓	✓
8	IDS Recipe, Order, and Workflow Order Management	10003610697/SPC	1.0	17-Dec-2017	–	✓	✓	–
9	IDS Exception Review Dashboard	10003484290/SPC	1.0	15-Dec-2017	–	✓	✓	–
10	IDS Data Manager	10003285813/SPC	1.0	15-Dec-2017	–	✓	✓	–
11	IDS EBR Execution	10003246565/SPC	1.0	18-Dec-2017	–	✓	✓	–
12	IDS Dispense Application	10003302778/SPC	1.0	15-Dec-2017	–	✓	✓	–
13	IDS Recipe and Workflow Designer	10003308714/SPC	1.0	18-Dec-2017	–	✓	✓	–
14	IDS EBR Data Model	10003288464/SPC	1.0	15-Dec-2017	–	✓	✓	–
15	IDS EBR Exception Recording	10003318832/SPC	1.0	15-Dec-2017	–	✓	✓	–
16	IDS GHS Symbols	10003285838/SPC	1.0	18-Dec-2017	–	✓	✓	–
17	IDS Protection of Intellectual Property	10003534415/SPC	1.0	17-Dec-2017	–	✓	✓	–
18	DCS 2.1 – IDS DCS Adapter	10003308671/SPC	1.0	15-Dec-2017	–	✓	✓	–
19	FRS Data Management	10003224063/SPC	1.0	10-Nov-2017	–	✓	✓	✓
20	FRS Execution Framework	10003236391/SPC	1.1	15-Dec-2017	–	✓	✓	✓
21	FRS Non-functional Requirements	10003256030/SPC	1.0	10-Nov-2017	–	✓	✓	✓
22	FRS Production Execution	10003526985/SPC	1.0	10-Nov-2017	–	✓	✓	✓
23	FRS Recipe and Workflow Management	10003236392/SPC	1.1	15-Dec-2017	–	✓	✓	✓
24	FRS Review and Approval	10003224065/SPC	1.0	10-Nov-2017	–	✓	✓	✓

#	Title	DIR	Revision ¹⁰	Approval ¹¹	S. Muench	M. Irmisch	S. Landes	M. Dittmer
25	FRS Runtime Data Management	10003253129/SPC	1.0	10-Nov-2017	–	✓	✓	✓
26	FRS DCS Phases	10003246501/SPC	1.1	18-Dec-2017	–	✓	✓	✓
27	FRS Dispense and Inline Weighing	10003224073/SPC	1.0	10-Nov-2017	–	✓	✓	✓
28	FRS EBR Phases	10003303200/SPC	1.0	10-Nov-2017	–	✓	✓	✓
29	FRS Equipment Automation Phases	10003303213/SPC	1.0	10-Nov-2017	–	✓	✓	✓
30	FRS Equipment Tracking Phases	10003303214/SPC	1.0	10-Nov-2017	–	✓	✓	✓
31	FRS IPC Phases	10003303216/SPC	1.1	15-Dec-2017	–	✓	✓	✓
32	FRS Material Tracking Phases	10003303218/SPC	1.0	10-Nov-2017	–	✓	✓	✓
33	FRS Output Weighing	10003224077/SPC	1.0	10-Nov-2017	–	✓	✓	✓
34	FRS Workflow Phases	10003447532/SPC	1.0	10-Nov-2017	–	✓	✓	✓
35	DCS 2.1 – FRS DCS Adapter	10003246500/SPC	1.1	15-Dec-2017	–	✓	✓	✓
36	Artifact Specific Risk Assessment	10003202617/REV	1.2	21-Dec-2017	–	✓	✓	✓
37	System Impact Analysis (SIA) & Risk Assessment (RA)	10003202619/REV	1.2	21-Dec-2017	–	✓	✓	✓
38	Risk Evaluation for MS SQL Server	10003626658/REV	1.0	19-Dec-2017	–	✓	✓	–
39	Assessment of Tools & Libraries	10003320986/REV	1.0	19-Dec-2017	–	✓	✓	–
40	Implementation of 21 CFR Part 11	10003305443/PUB	1.0	18-Dec-2017	–	✓	✓	✓
41	Implementation of Annex 11	10003305444/PUB	1.0	18-Dec-2017	–	✓	✓	✓
42	Implementation of Data Integrity	10003305445/PUB	1.0	18-Dec-2017	–	✓	✓	✓
43	TM Installation - Express Edition	10003224082/PUB	1.0	21-Dec-2017	–	✓	✓	✓

#	Title	DIR	Revision ¹⁰	Approval ¹¹	S. Muench	M. Irmisch	S. Landes	M. Dittmer
44	TM Installation - Enterprise Edition	10003224081/PUB	1.2	21-Dec-2017	–	✓	✓	✓
45	TM Installation - Building Blocks	10003224080/PUB	1.2	18-Dec-2017	–	✓	✓	✓
46	TM Installation – Upgrade	10003224083/PUB	1.0	21-Dec-2017	–	✓	✓	✓
47	TM Administration	10003236393/PUB	1.0	12-Dec-2017	–	✓	✓	✓
48	TM Configuration and Extension - Volume 1	10003224078/PUB	1.0	18-Dec-2017	–	✓	✓	✓
49	TM Configuration and Extension - Volume 2	10003224078/PUB	1.0	18-Dec-2017	–	✓	✓	✓
50	TM Configuration and Extension - Volume 3	10003224078/PUB	1.0	18-Dec-2017	–	✓	✓	✓
51	TM Configuration and Extension - Volume 4	10003224078/PUB	1.0	18-Dec-2017	–	✓	✓	✓
52	TM Configuration and Extension - Volume 5	10003224078/PUB	1.0	18-Dec-2017	–	✓	✓	✓
53	TM Developing System Building Blocks	10003323083/PUB	1.0	18-Dec-2017	–	✓	✓	✓
54	TM Phases of the Dispense Package	10003285729/PUB	1.0	12-Dec-2017	–	✓	✓	✓
55	TM Phases of the DCS Package	10003246504/PUB	1.0	12-Dec-2017	–	✓	✓	✓
56	TM Phases of the IPC Package	10003516655/PUB	1.0	12-Dec-2017	–	✓	✓	✓
57	DCS 2.1 - TM DCS Adapter	10003246498/PUB	1.0	12-Dec-2017	–	✓	✓	✓
58	UM Data Manager	10003224100/PUB	1.0	15-Dec-2017	–	✓	✓	✓
59	UM Production Management	10003547197/PUB	1.0	15-Dec-2017	–	✓	✓	✓
60	UM Production Execution	10003547199/PUB	1.0	15-Dec-2017	–	✓	✓	✓
61	UM Production Responses	10003547200/PUB	1.0	15-Dec-2017	–	✓	✓	✓
62	UM Recipe and Workflow Designer – Volume 1	10003547201/PUB	1.0	15-Dec-2017	–	✓	✓	✓

#	Title	DIR	Revision ¹⁰	Approval ¹¹	S. Muench	M. Irmisch	S. Landes	M. Dittmer
63	UM Recipe and Workflow Designer – Volume 2	10003547201/PUB	1.0	15-Dec-2017	–	✓	✓	✓
64	UM Recipe and Workflow Designer – Volume 3	10003547201/PUB	1.0	15-Dec-2017	–	✓	✓	✓
65	UM Recipe and Workflow Designer – Volume 4	10003547201/PUB	1.0	15-Dec-2017	–	✓	✓	✓
66	UM DCS Phases	10003246499/PUB	1.0	18-Dec-2017	–	✓	✓	✓
67	UM Dispense and Weighing Phases	10003224102/PUB	1.0	18-Dec-2017	–	✓	✓	✓
68	UM EBR Phases	10003547148/PUB	1.0	18-Dec-2017	–	✓	✓	✓
69	UM Equipment Phases	10003547202/PUB	1.0	18-Dec-2017	–	✓	✓	✓
70	UM IPC Phases	10003547203/PUB	1.0	18-Dec-2017	–	✓	✓	✓
71	UM Material Tracking Phases	10003547205/PUB	1.0	18-Dec-2017	–	✓	✓	✓
72	Release Notes	10003224098/PUB	1.1	15-Dec-2017	–	✓	✓	✓
73	Traceability Matrix	10003202613/FRM	1.0	22-Dec-2017	–	✓	✓	–
74	Final Test Status Report	10003547212/REV	1.0	22-Dec-2017	–	✓	–	–
75	Test Summary Report	10003547210/REV	1.0	22-Dec-2017	–	✓	✓	✓
76	Building Blocks - Compatibility Matrix	10003224093/FRM	1.0	22-Dec-2017	–	✓	✓	✓
77	Project Closure Report	10003547206/REV	1.0	22-Dec-2017	–	✓	✓	✓
78	Quality Document (→ this document)	10003305482/PUB	1.0	22-Dec-2017	–	✓	✓	✓
79	Quality Certificate	10003305481/PUB	1.0	22-Dec-2017	–	✓	✓	✓
Total number of documents signed (per approver)					1	78	78	62

Table 1: List of all documents and final approvals for PharmaSuite 8.4

Handwritten Signatures for Release Documents

With my handwritten signature, I certify that

- the data shown in the above Table 1 is correct and complete with respect to the documents that I have approved.
- I have electronically approved the latest revision of the respective documents, as listed in Table 1, through the ‘click signature’ implemented in Rockwell Software’s Document Management System [6].
- I have held the role(s) as indicated for release PharmaSuite 8.4.

Approver	Role(s)	Date	Signature
Stefan Muench	Project Manager (until Sep 2017) Test Manager (until Sep 2017)		
Martin Irmisch	Project Manager (since Oct 2017) Test Manager (since Oct 2017)		
Steffen Landes	Development Manager Information Development Manager		
Martin Dittmer	Product Manager		

Table 2: Handwritten signatures of all signees of PharmaSuite 8.4 release documents

Handwritten Signatures for Executed Tests

The following Table 3 and Table 4 provide the handwritten signature per tester. The total number of documents signed in Table 3 and Table 4 reflects the data for release PharmaSuite 8.4.

With my handwritten signature, I certify that

- the data shown in the Final Test Status Report [22], section 2.2 Test Case Details, is correct and complete with respect to the tests that I have executed within the course of PharmaSuite 8.4.
- I have personally executed the respective test cases, as listed in [22], section 2.2 Test Case Details. These executions have been recorded by ApTest Manager [7], documented in column “Actual Tester” in [22], and approved by the Test Manager.

Tester	# of test cases executed	Date	Signature
Svilen Angelov	2		
Hubert Bach	3		
Annette Blob	20		
Andreas Grossmann	4		
Bernd Lorenz	39		
Alexander Mattes	35		

Table 3: Handwritten signatures of signees 1-7 of PharmaSuite 8.4 executed tests

With my handwritten signature, I certify that

- the data shown in the Final Test Status Report [22], section 2.2 Test Case Details, is correct and complete with respect to the tests that I have executed within the course of PharmaSuite 8.4.
- I have personally executed the respective test cases, as listed in [22], section 2.2 Test Case Details. These executions have been recorded by ApTest Manager [7], documented in column “Actual Tester” in [22], and approved by the Test Manager.

Tester	# of test cases executed	Date	Signature
Shiva Prasad Namani	43		
Basavalingamma Sahu	69		
Janmejaya Sahu	49		
Satya Sivanjali	46		
Ignaz Wangler	31		
Total (Table 3 and Table 4)	341		

Table 4: Handwritten signatures of all signees 8-13 of PharmaSuite 8.4 executed tests